UnityConnect

Reference Guide

October 2017 Version 2.0



Contents

Overview of UnityConnect	1
Introduction	1
Illustration of UnityConnect	2
About This Guide	3
How Will You Use UnityConnect in Your Laboratory?	3
Data File Transfer with Desktop Software	3
Data File Transfer with Web-based Software	4
Data Stream with Desktop Software	5
Data Stream with Web-based Software	
Organization of UnityConnect	6
How UnityConnect Works	7
Example of Matching Information	8
Why UnityConnect?	g
Working with the Customized Transformer	10
Copying a Transformer	10
Editing a Transformer	11
Renaming a Transformer	11
Editing the Input File Type	13
Deleting a Transformer	15
Configuring a Transformer	16
Migrating Settings from a Previous Version	17
UnityConnect Tabs	21
Data File Options Tab	22
Example Data File	22
Traditional Configuration	23
Example Data File Options Tab—Traditional Configuration	24
Universal Flex File Configuration	28
Example of Data File Columns	28
Example Data File Options Tab—Universal Flex File	29
Overview of Other Configuration Tabs	29
Test Code Definitions Tab	30
Example Test Code Definitions Tab	30
Lots Tab	31
Example Lots Tab	32
Special Cases Tab	38
Instruments Tab	38
Example Instruments Tab	38
Tests Tab	46
Example Tests Tab	46
Qualitative/Semi-Quantitative Tab	52
Config Gens Tab	52

Working with QC Data Files	. 53
Overview	. 53
QC Data File Content	. 54
QC Data File Essentials	. 55
Maintaining QC Data File Integrity	. 55
Multiple QC Data Files	. 55
Single Point Data versus Summary Data	. 56
Single Point Data	. 56
Summary Data	. 56
Unfiltered Data versus Filtered Data	. 56
Unfiltered Data	. 56
Filtered Data	. 57
Creating a QC Data File	. 58
How Often to Create a QC Data File	. 58
Best Practices—Transform Data Often	. 58
Copying a QC Data File	. 59
Configuring UnityConnect for Data File Transfer with Desktop Software	62
Before You Begin	. 62
Configuration Flow Chart	. 63
Output Path	. 64
Configuring Watch Folder and Automatic Import	. 64
Configuring a Watch Folder	. 64
Configuring Automatic Import	
Configuring Import Options	. 68
Import Options—Using QC Data Management Software for SPC Rule	
Evaluation	. 69
Import Options—Not Using QC Data Management Software for SPC Ru	ıle
Evaluation	. 71
Downloading a Transformer	. 73
Configuring E-Mail Notification	. 74
Configuring the Default Date Settings	. 75
Configuring UnityConnect for Data File Transfer with Web-based Softw 76	are
Before You Begin	76
Configuration Flow Chart	
Setting Up UnityConnect with Web-based Software	
Output Path	
Configuring a Watch Folder	
Configuring Internet Options	
Testing the Internet Connection	
If the Internet Connection Fails	
Configuring Auto Post	
	. От

Downloading the Database	86
Downloading Lab Setup and Code Lists	87
Downloading a Transformer	88
Configuring E-Mail Notification	89
Configuring the Default Date Settings	90
Configuring Import Settings	
Import Settings—Using QC Data Management Software for SPC Evalu 92	ation
Import Settings—Not Using QC Data Management Software for SPC F Evaluation	
Configuring UnityConnect for a Data Stream with Desktop Software	
Before You Begin	
Configuration Flow Chart	
Output Path	
Configuring a Watch Folder and Automatic Report	
Configuring Automatic Impact	
Configuring Automatic Import	
Import Options—Using QC Data Management Software for SPC Rule	102
Evaluation	103
Import Options—Not Using QC Data Management Software for SPC R	
Evaluation	
Downloading a Transformer	
Configuring E-Mail Notification	
Configuring UnityConnect for a Data Stream with Web-based Software	э
110 Refere Year Regin	110
Before You Begin	
Configuration Flow Chart	
Setting Up UnityConnect for Web-based Software	
Output Path Configuring a Watch Folder	
Configuring Internet Options	
Testing the Internet Connection	
If the Internet Connection Fails	
Configuring Auto Post	
Downloading the Database	
Downloading the Database	
Downloading a Transformer	
Configuring E-Mail Notification	
Configuring Import Settings	
Import Settings—Using QC Data Management Software for SPC Rule	
Evaluation	197



Import Settings—Not Using QC Data Management Software for S	SPC Rule
Evaluation	
Data Acquisition Manager	129
Before You Begin	129
Overview	130
Installing the Service	130
Working with Devices	130
Adding a Device	130
Deleting a Device	131
Modifying a Device	
Connecting a Single Device	
Connecting All Devices	
Disconnecting a Single Device	
Disconnecting All Devices	
Starting and Stopping the Service	
Stopping the Service	
Exception Files	
Viewing Exception Files	
Deleting Exception Files.	
Device Logs.	
Viewing Device Logs Deleting Device Logs	
Deleting Device Logs	109
To a facility Data for Dealth and fine	1.40
Transforming Data for Desktop Software	
Overview	
Best Practices—Transform Data Often	
Transforming Data — Data File Transfer Method	
Before You Begin	
Reprocessing QC Data Files	
Best Practicies—Review the Number of Data Records Transformed	
Working with Multiple QC Data Files	
Configuring UnityConnect to Transform a Directory of Files	
Transforming a Directory of Files	
Best Practices—Review the Number of Data Records Transformed	
Transforming Data for Web-based Software	151
Overview	
Before You Begin	
Before You Begin	
Best Practices—Transform Data Often	
Transforming Data — Data File Transfer Method	
Reprocessing QC Data Files	

Best Practices—Review the Number of Data Records Transformed	158
Working with Multiple QC Data Files	159
Configuring UnityConnect to Transform a Directory of Files	159
Transforming a Directory of Files	160
Best Practices-Review the Number of Data Records Transformed	161
New Configuration in UnityConnect	162
Overview	
Summary of New Configuration	
Configuring a New Lot	
Configuring a New Lot When Using QC Data Management Software	
SPC Rule Evaluation	
Duplicating a Lot in Your QC Data Management Software	
Duplicating a Lot in Unity Desktop	
Duplicating a Lot in Unity Real Time	
Duplicating a Lot in UnityWeb 2.0	
Dulicating a Lot in Unity Real Time online	
Completing New Configuration in UnityConnect	
Configuring a New Lot When Not Using QC Data Management Softw	
for SPC Rule Evaluation	
Configuring Import Settings	
Completing New Configuration in UnityConnect	
Other Types of Configuration	
Configuring New Data File Options	
Configuring a New Instrument	
Configuring a New Test	
	100
VITROS Slide Generation Numbers and UnityConnect	201
Overview	
Slide Generation Number Decision Chart	
Example Config Gens Tab.	
Best Practices—Manage Manual VITROS Slide Generation Changes	
Managing Slide Generation Numbers in UnityConnect	
Slide Generation Worksheet	
Changing SPC Rules for a New VITROS Test	
For Unity Desktop and Unity Real Time	
For UnityWeb 2.0 and Unity Real Time online	
1 of Officyweb 2.0 and Officy Near Fifthe Offithe	214
Viewing the Rejection Log for Desktop Software	215
Overview	
Before You Begin	
-	
Using the Rejection Log Viewing the Rejection Log	
Best Practices — Review the Rejection Log by Error Type After Each	∠ 1 /



Transformation	217
Deleting the Rejection Log	220
Best Practices—Note the Issues Needing Adjustment and Delete the Reju	ection
Log	
Rejection Log Messages	
Error: Create new tests if necessary is disabled	221
Error: Data entry locked for this test	
Error: Date out of sequence	
Solution—Unity Desktop and Unity Real Time	226
Error: Lab closed	227
Error: Lot closed	227
Error: Lot expired	228
Error: Lot number undefined	228
Import Options—Using QC Data Management Software for SPC Rule	Э
Evaluation	229
Import Options—Not Using QC Data Management Software for SPC	Rule
Evaluation	233
Error: Result invalid	236
Error: Time out of sequence	237
Viewing the Rejection Log for Web-based Software	240
Overview	
Before You Begin	
Viewing the Rejection Log	
Best Practices—Review the Rejection Log by Error Type After Each	243
Transformation	040
Deleting the Rejection Log	
Best Practices—Note the Issues Needing Adjustment and Delete the Rej	
Log	
Rejection Log Messages	
Error: Create new tests if necessary is disabled	
Error: Data entry locked for this test	
Error: Date out of sequence	
Error: Lab closed	
Error: Lot closed	
Error: Lot expired	
Error: Lot number undefined	
Import Settings—Using QC Data Management Software for SPC Rule	
Evaluation	
Import Settings—Not Using QC Data Management Software for SPC	
Evaluation	
Error: Result invalid	
Error: Time out of sequence	257

Viewing Watch Folder Logs	259
Overview	259
Viewing Log Files	259
Viewing Processed Files	
Viewing Uploaded Files	265
Selecting Log Options	268
Cleaning Files Automatically	268
Limit Files to Options	268
Cleaning Log Files Manually	269
Limit Files to Options	
Backup and Restore	271
Overview	271
Backing Up Data	271
Restoring Data from a Backup	273
License Agreement, Warranty Information, and Trademark Notices	276
License Agreement	276
Warranty Information	277
Trademark Notices	277
Haderian Notices	211

Overview of UnityConnect

In This Chapter

Introduction	1
Illustration of UnityConnect	2
About This Guide	3
Organization of UnityConnect	6
How UnityConnect Works	7
Why UnityConnect?	9
Working with the Customized Transformer	10
UnityConnect Tabs	21

Introduction

UnityConnect is Bio-Rad's newest and most complete connectivity solution that allows "any laboratory" to get connected to the right Unity QC data management solution.


UnityConnect consists of Bio-Rad software and optional hardware to interface QC data from instruments, middleware, and Laboratory Information Systems (LIS) directly into your Unity QC data management software.



Note: Unity QC data management software includes desktop software (Unity Desktop and Unity Real Time) and Web-based software (UnityWeb 2.0 and Unity Real Time online).

UnityConnect provides effective QC data processing and automates the import of QC data, thereby eliminating manual data entry and errors.

Illustration of UnityConnect

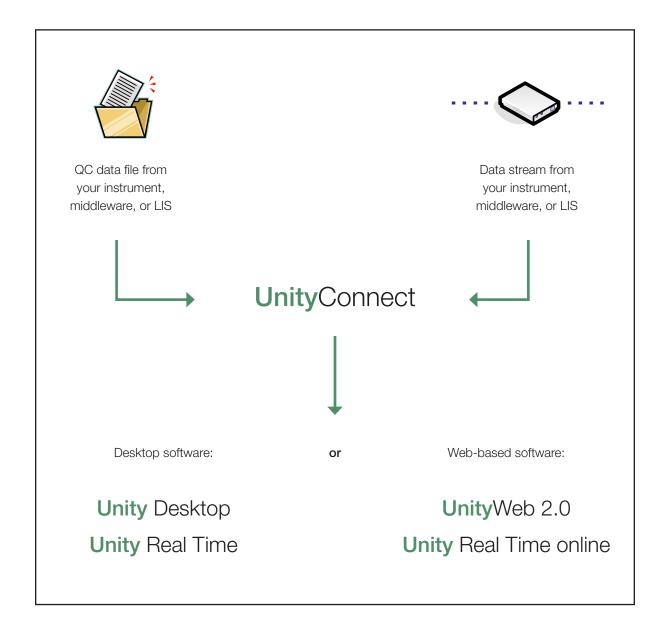

About This Guide


This guide contains information and instructions for all UnityConnect configurations. However, you only need to refer to chapters specific to the configuration used in your laboratory.

How Will You Use UnityConnect in Your Laboratory?

See the following sections for references to your specific configuration:

- Data file transfer with desktop software (page 3)
- Data file transfer with Web-based software (page 4)
- Data stream with desktop software (page 5)
- Data stream with Web-based software (page 5)


Data File Transfer with Desktop Software

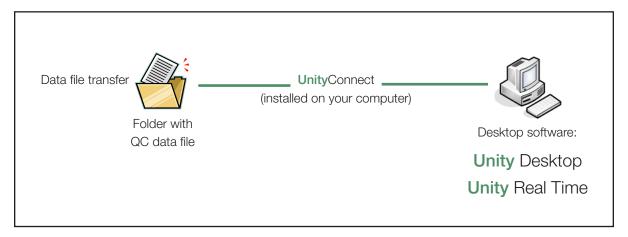

See the following chapters for more information about this configuration:

- Chapter 1, "Overview of UnityConnect" (this chapter)
- Chapter 2, "Working with QC Data Files"
- Chapter 3, "Configuring UnityConnect for Data File Transfer with Desktop Software"
- Chapter 8, "Transforming Data for Desktop Software"
- Chapter 10, "New Configuration in UnityConnect"
- Chapter 11, "VITROS Slide Generation Numbers and UnityConnect" (if you use slide generation numbers)
- Chapter 12, "Viewing the Rejection Log for Desktop Software"
- Chapter 14, "Viewing Watch Folder Logs"
- Chapter 15, "Backup and Restore"

Data File Transfer with Web-based Software

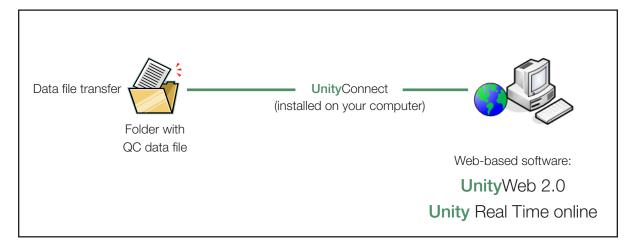

See the following chapters for more information about this configuration:

- Chapter 1, "Overview of UnityConnect" (this chapter)
- Chapter 2, "Working with QC Data Files"
- Chapter 4, "Configuring UnityConnect for Data File Transfer with Web-based Software"
- Chapter 9, "Transforming Data for Web-based Software"
- Chapter 10, "New Configuration in UnityConnect"
- Chapter 11, "VITROS Slide Generation Numbers and UnityConnect" (if you use slide generation numbers)
- Chapter 13, "Viewing the Rejection Log for Web-based Software"
- Chapter 14, "Viewing Watch Folder Logs"
- Chapter 15, "Backup and Restore"


Data Stream with Desktop Software



See the following chapters for more information about this configuration:

- Chapter 1, "Overview of UnityConnect" (this chapter)
- Chapter 5, "Configuring UnityConnect for a Data Stream with Desktop Software"
- Chapter 7, "Data Acquisition Manager"
- Chapter 10, "New Configuration in UnityConnect"
- Chapter 11, "VITROS Slide Generation Numbers and UnityConnect" (if you use slide generation numbers)
- Chapter 12, "Viewing the Rejection Log for Desktop Software"
- Chapter 14, "Viewing Watch Folder Logs"
- Chapter 15, "Backup and Restore"


Data Stream with Web-based Software

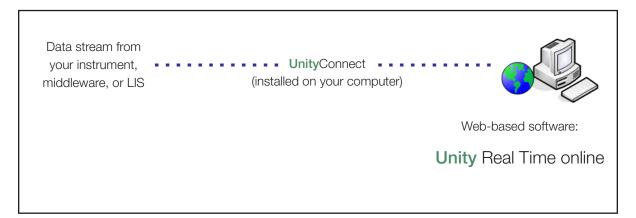


See the following chapters for more information about this configuration:

- Chapter 1, "Overview of UnityConnect" (this chapter)
- Chapter 6, "Configuring UnityConnect for a Data Stream with Web-based Software"
- Chapter 7, "Data Acquisition Manager"
- Chapter 10, "New Configuration in UnityConnect"
- Chapter 11, "VITROS Slide Generation Numbers and UnityConnect" (if you use slide generation numbers)
- Chapter 12, "Viewing the Rejection Log for Desktop Software"
- Chapter 14, "Viewing Watch Folder Logs"
- Chapter 15, "Backup and Restore"

Organization of UnityConnect

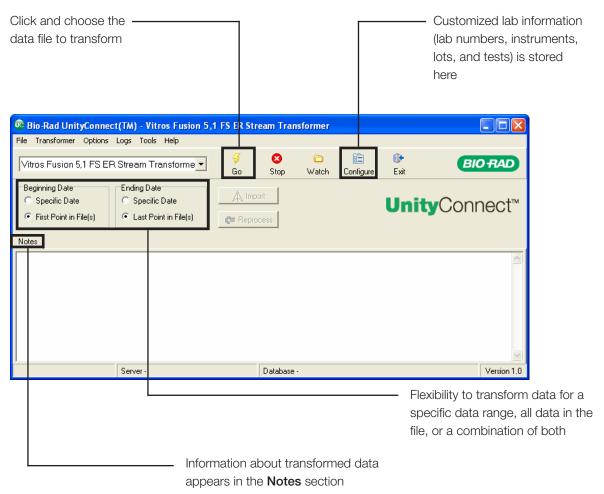

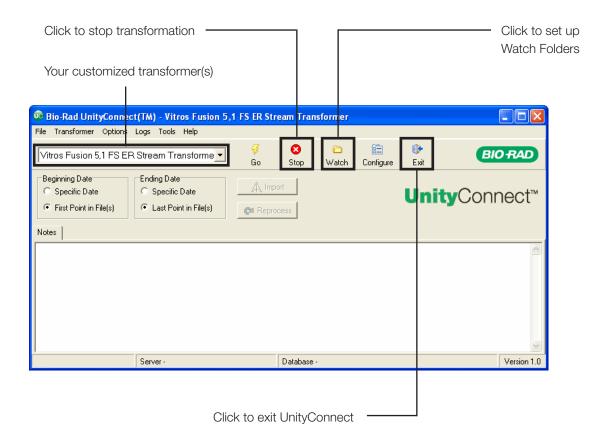

How UnityConnect Works

Laboratories around the world use unique names to identify information in their instruments, middleware, and LIS. For example, a laboratory may refer to the Siemens Diagnostics ADVIA Centaur as "Centaur 1." Everyone in the laboratory recognizes the instrument as Centaur 1 so the LIS is programmed with this information.

However, in order to compare QC data with other laboratories running tests on a Siemens Diagnostics ADVIA Centaur, standard identifiers are needed.

UnityConnect provides the ability to match the descriptions used for QC items (instruments, lots, test) in the LIS to the Bio-Rad description in your Unity software.

Example of Matching Information




Note: The descriptions in the following illustration are examples only. The descriptions in your QC data file will be different.

Bio-Rad descriptions in the QC data Descriptions in the QC data file management software Instrument names Instrument names AXSYM Abbot AxSYM D10 Bio-Rad D-10 RXL1 Siemens Dimension RxL Control names Control names Urine Chemistry UCHEM **CSF** Spinal Fluid IMMUN ◀ Immunoassy Plus Analyte names Analyte names **AMLY** Amylase **GLUC** Glucose VANC Vancomycin


Why UnityConnect?

The sole purpose of UnityConnect is to transform your QC data from the format of your instrument, middleware, or LIS into the Bio-Rad format. UnityConnect reads the unique information in the QC data file and transforms it so it can be added to the Unity peer group in a standardized manner.

Example of QC data file from LIS

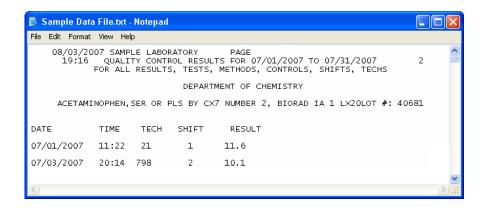

Example of Bio-Rad format

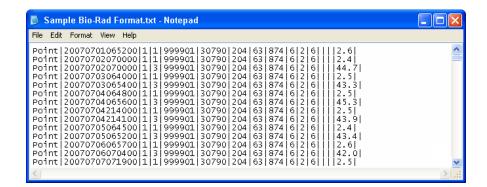

Data is automatically imported into your QC data management software after transformation

QC data management software

Working with the Customized Transformer


UnityConnect is a Bio-Rad connectivity solution consisting of one or more customized transformers created especially for your instrument, middleware, or LIS. UnityConnect provides the following capabilities for the transformer:

- Copy a transformer (page 10)
- Editing a transformer (page 11)
- Deleting a transformer (page 15)
- Configuring a transformer (page 16)
- Migrate settings from a previous version (page 17)

Copying a Transformer

Copying a transformer is helpful when you are using an instrument transformer and have multiple instruments of the same type.

1 Select the transformer you want to copy from the list of transformers.

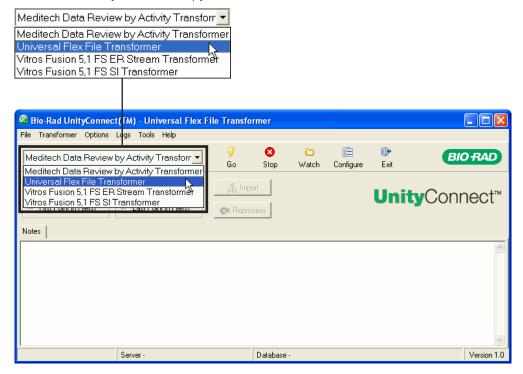

2 Click the **Transformer** menu and then click **Copy this Transformer**.



The Copy Transformer dialog box appears.



3 Click ✓ OK.



Tip: You can rename the copy of the transformer for easy identification. See "Renaming a Transformer" on page 11 for more information.

Editing a Transformer


You can edit a transformer by:

- Renaming a transformer (page 11)
- Editing the input file type (page 13)

Renaming a Transformer

Renaming a transformer is useful if you have more than one transformer and want the display name for the transformer the same as the instrument, middleware, or LIS name used in your laboratory.

1 Select the transformer you want to edit from the list of transformers.

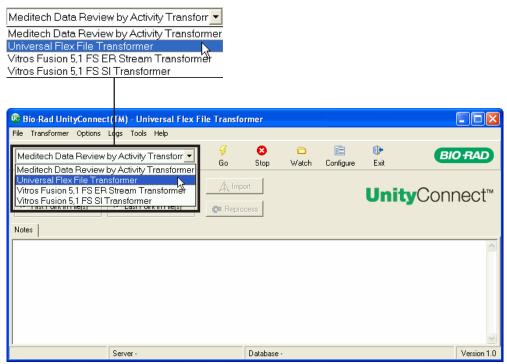

2 Click the **Transformer** menu and then click **Edit this Transformer**.




3 Click in the **Name** field and type the new name for the transformer.



4 Click VOK.

Editing the Input File Type

1 Select the transformer you want to edit from the list of transformers.

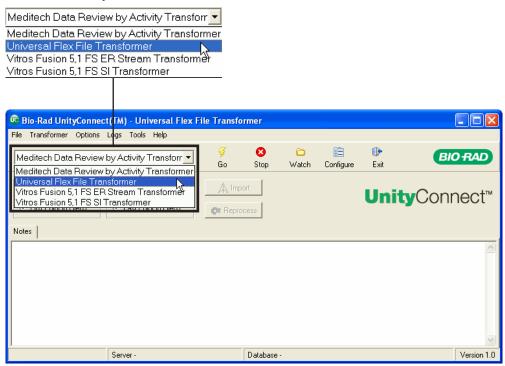

2 Click the Transformer menu and then click Edit this Transformer.



3 Select an option for the input file.





Note: Some Laboratory Information Systems (LIS) output multiple files when creating QC data files. One example is the Orchard LIS. Processing a directory is useful in this instance.


4 Click VOK.

Deleting a Transformer



Important: Deleting a transformer permanently removes the transformer from the UnityConnect software. You must contact Bio-Rad to use the transformer in the future. Do not delete a transformer unless requested to do so by Bio-Rad.

1 Select the transformer you want to delete from the list of transformers.

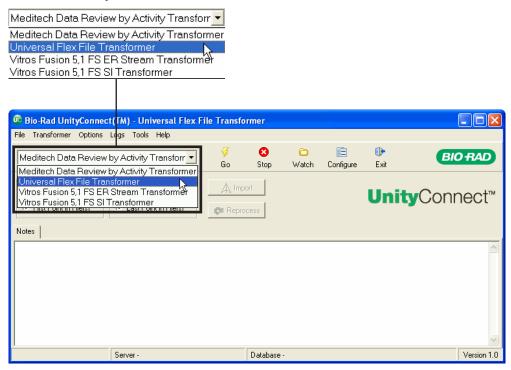

2 Click the Transformer menu and then click Delete this Transformer.



A message appears asking for confirmation.




3 Click <u>o</u>k .

Configuring a Transformer




Important: Only users with a valid database login can change the configuration for a transformer. See Chapter 10, "New Configuration in UnityConnect" for more information.

1 Select the transformer you want to configure from the list of transformers.

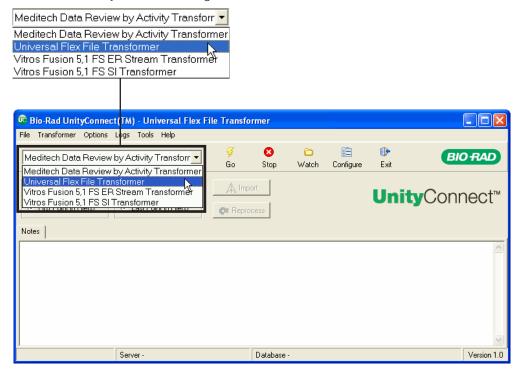

2 Click the **Transformer** menu and then click **Configure this Transformer**.

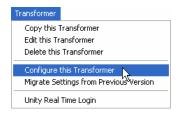

3 Log in to your QC data management software database if prompted.

The **Data File Options** tab appears.




Note: See "UnityConnect Tabs" on page 21 for more information about configuration.

Migrating Settings from a Previous Version

The migrate settings feature provides the ability to upgrade to newer versions of transformers when available.

1 Select the transformer you want to migrate from the list of transformers.

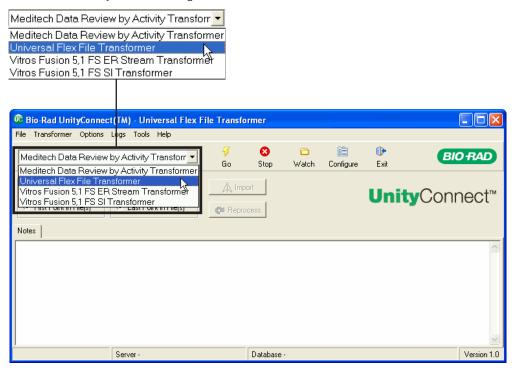

2 Click the **Transformer** menu and then click **Migrate Settings from Previous Transformer**.



The Migrate Settings from Previous Version dialog box appears.




3 Click (ellipsis button).


Migrate Settings from Previous Version [Universal Flex File Transformer]

Previous Version Setting File:

Migrate X Cancel

4 Select the configuration file.

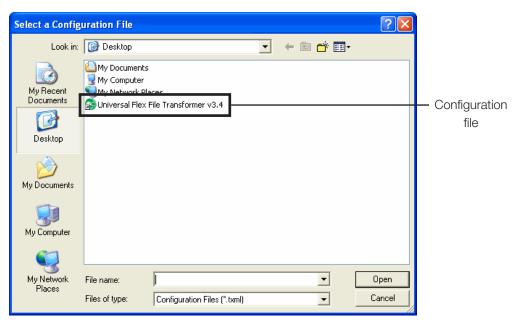

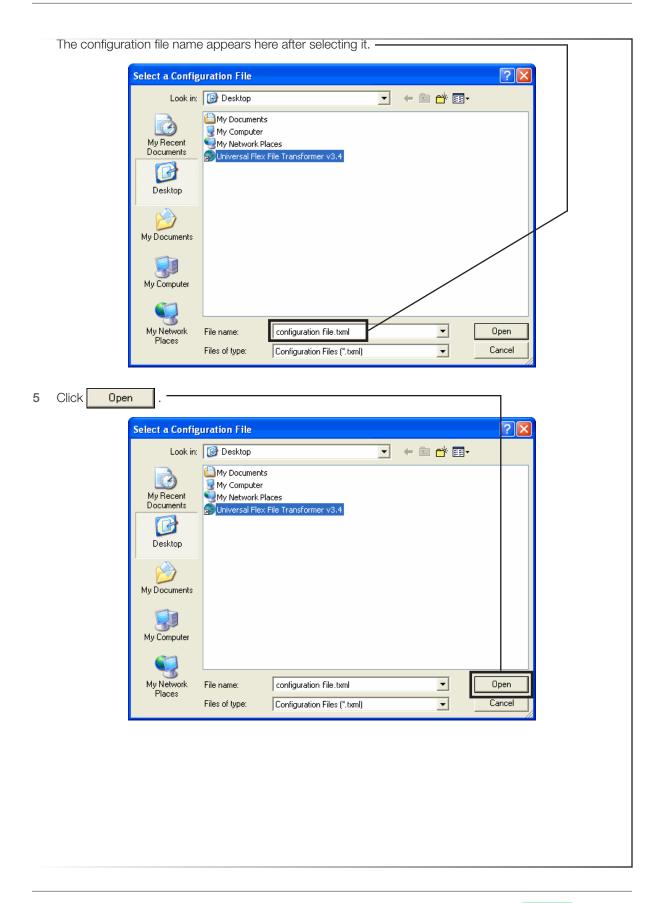

6 The Migrate Settings dialog box shows the selected file and path.



- 7 Click <u>Migrate</u>.
- 8 A message appears asking for confirmation.



9 Click <u>O</u>K .


UnityConnect Tabs

Bio-Rad created your customized UnityConnect software based on all the information in the original QC data file provided by your laboratory. After the initial configuration is complete, new configuration is required only if a new QC item (such as an instrument, lot, or test) is added to your QC data file.

UnityConnect processes the QC data file and, if necessary, displays messages prompting you to match or configure the new information in the QC data file with the information in your QC data management software. UnityConnect automatically determines the information needing configuration and presents the appropriate tab where the configuration can be completed.



Important: UnityConnect is a customized product. Therefore, your UnityConnect software may not present all of the tabs shown in this section. Click Plelp located in the lower right corner of the window to view customized online help for your transformer.

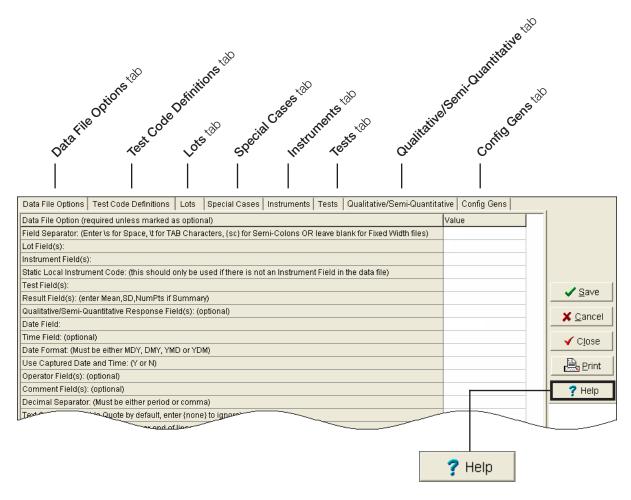

The following sections provide more information about the different tabs in UnityConnect:

- Data File Options tab (page 22)
- Test Code Definitions tab (page 30)

- Lots tab (page 31)
- Special Cases tab (page 38)
- Instruments tab (page 38)
- Tests tab (page 46)
- Qualitative/Semi-Quantitative tab (page 52)
- Config Gens tab (page 52)



Note: New configuration is discussed in more detail in Chapter 10, "New Configuration in UnityConnect".


Data File Options Tab

The **Data File Options** tab is the most important tab in UnityConnect. The settings on this tab tell UnityConnect where to find important information in your QC data file critical for successful data transformation.

Bio-Rad created your customized UnityConnect software based on the QC data file your laboratory provided. Your laboratory's QC data file provided the "blueprint" for Bio-Rad to determine the unique data identifiers used in your laboratory and their specific location within the QC data file. Your unique settings drove the UnityConnect software design process and the options and settings available on the **Data File Options** tab.

Data file formats vary from system to system. However, the essential components are the same as shown in the example below.

Example Data File

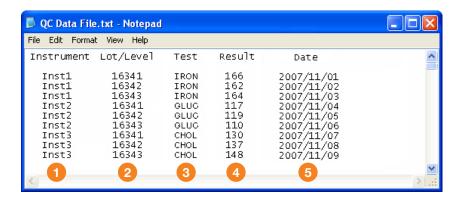

- 1 Unique designation for each instrument
 - Each instrument has a unique designation in the example: Inst1, Inst2, and Inst3.
- 2 Unique designation for the control and level

 Each control and level has a unique designation in the example: 16341 (lot 16340-level 1), 16342 (lot 16340-level 2), and 16343 (lot 16340-level 3).

3 Unique test designation

Each test has a unique designation in the example: IRON, GLUC, and CHOL.

4 Test result

Each result has a valid numeric value in the example.



Note: UnityConnect cannot transform data with a value of 0 (zero) or a negative numeric value.

Date of each test result

Each result has an associated date. In the example above the format is year/month/day.

The settings on the **Data File Options** tab reflect your unique format. Once configured, these settings remain unchanged as long as the QC data files used to transform QC data are created in the identical manner as the QC data file used to create your UnityConnect software.

In addition, adding new instruments, lots, or tests does not disrupt the functionality of your UnityConnect software as long as you set up the new QC item (such as an instrument, lot, or test) with a unique identifier and in the same format and location in the QC data file as existing QC items.

The Data File Options tab can be organized into one of two general types within UnityConnect:

- Traditional configuration (this page)
- Universal Flex File configuration (page 28)




Important: UnityConnect is a customized product. Therefore, the example information shown in this section may not be identical to your UnityConnect. Click Help located in the lower right corner of the window to view customized online help for your transformer.


Traditional Configuration

Traditional configuration uses fields and check boxes to configure the Data File Options tab.

Example Data File Options Tab—Traditional Configuration

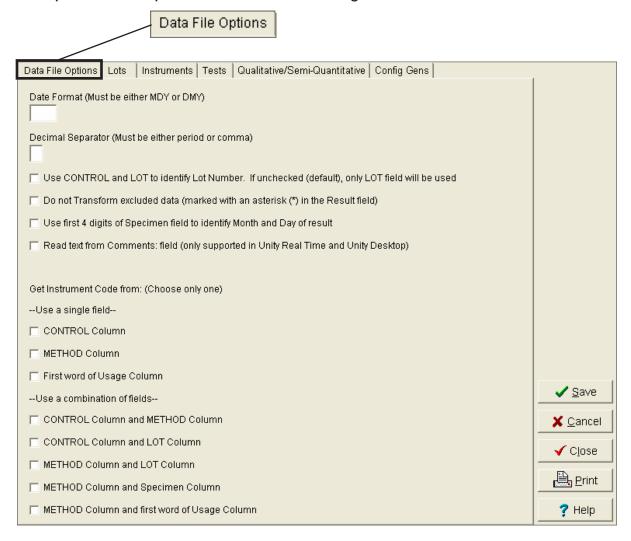

Example Data File Options Tab - Traditional Configuration (continued)

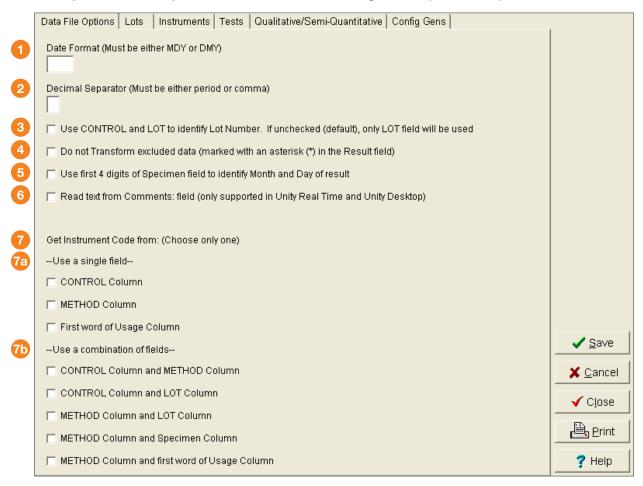

1 Date Format

This field identifies how the date is specified in the QC data file. Date formats can use a two digit or four digit year. Date format options are:

- MDY (month-day-year)
- **DMY** (day-month-year)

Important: The date format must be in uppercase letters. For example, MDY or DMY.

2 Decimal Separator

The decimal separator must be designated as a period (,) or a comma (,). It is the character used as the decimal separator on the local computer.

Example: 1.25 (decimal) or 1,25 (comma).



Note: The standard separator for United States customers is a period.

- Use CONTROL and LOT to identify Lot Number. If unchecked (default), only LOT field will be used UnityConnect reads the lot number from the Lot field if this check box is not selected. UnityConnect reads the lot number from the Control and Lot fields if this check box is selected.
- Do not Transform excluded data (marked with an asterisk (*) in the Result field)

 UnityConnect excludes outliers in the Results field that are marked with an asterisk (*) if this check box is selected.
- Use first 4 digits of Specimen field to identify Month and Day of result

 UnityConnect determines the month and day of the result using the first four digits of the Specimen field if this check box is selected.
- 6 Read text from Comments: field (only supported in Unity Desktop and Unity Real Time)

 UnityConnect transforms and imports text from the Comments field into the Comments field of your QC data management software if this check box is selected.



Note: This feature is not compatible with all software.




Get Instrument Code from (Choose only one)

These settings are used to identify which field(s) in the QC data file contain(s) the instrument designation.



Important: One check box must be selected so UnityConnect knows which field(s) to read from the original QC data file.

The instrument designation is in the Method field in the example below.

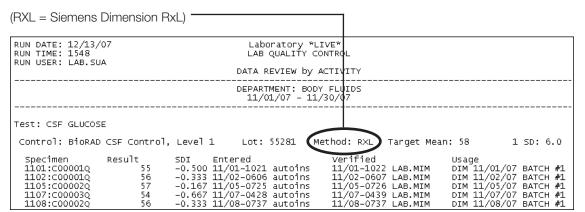



--Use a single field--

CONTROL Column

UnityConnect reads the instrument designation from the Control field in the QC data file if this check box is selected.

METHOD Column

UnityConnect reads the instrument designation from the Method field in the QC data file if this check box is selected.

First word of Usage Column

UnityConnect reads the instrument designation from the first word of the Usage column field in the QC data file if this check box is selected.



--Use a combination of fields--

CONTROL Column and METHOD Column

UnityConnect reads the instrument designation from a combination of the Control and Method fields in the QC data file if this check box is selected.

CONTROL Column and LOT Column

UnityConnect reads the instrument designation from a combination of the Control and Lot fields in the QC data file if this check box is selected.

METHOD Column and LOT Column

UnityConnect reads the instrument designation from a combination of the Method and Lot fields in the QC data file if this check box is selected.

METHOD Column and Specimen Column

UnityConnect reads the instrument designation from a combination of the Method and Specimen fields in the QC data file if this check box is selected.


METHOD Column and first word of Usage Column

UnityConnect reads the instrument designation from a combination of the Method and the first word of the Usage fields in the QC data file if this check box is selected.


Universal Flex File Configuration

The Universal Flex File configuration uses column numbers or fixed width columns to configure the **Data File Options** tab. Although column numbers do not actually appear in the QC data file, the columns are discernible as vertically arranged information. Columns are numbered beginning with one; column one is the leftmost column in the data file.

Example of Data File Columns

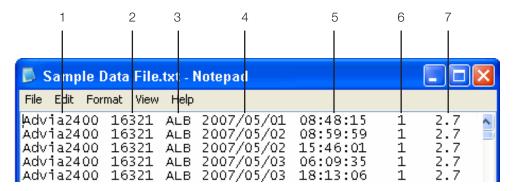

The column numbers from the QC data file are used in the appropriate fields on the **Data File Options** tab to designate corresponding information between the QC data file and UnityConnect. An example is shown in the illustration below.

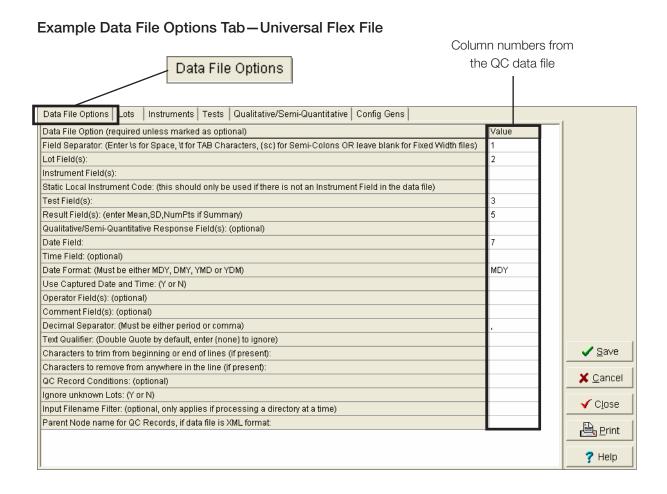



Note: Some fields on the Data File Options tab are optional.

If more than one field is used as a designation in the QC data file, both column numbers are used and separated by a comma (,).

Overview of Other Configuration Tabs

After the initial configuration is complete, new configuration is required only if a new QC item (such as an instrument, lot, or test) is added to the QC data file.

UnityConnect automatically determines the information needing configuration and presents the appropriate tab where the configuration can be completed. UnityConnect guides you to the appropriate tab where you complete configuration.



Important: UnityConnect is a customized product. Therefore, your UnityConnect software may not present all of the tabs shown in this section. Click located in the lower right corner of the window to view customized online help for your transformer.


Test Code Definitions Tab

The **Test Code Definitions** tab only appears for instrument transformers that define the analyte as a number. UnityConnect converts the number of the analyte to the analyte name and shows the name on the **Test** tab. For VITROS instrument transformers, the **Test Code Definitions** tab is also used to define wet tests versus dry tests.




Important: Test configuration in UnityConnect varies according to specific information in the QC data file. Therefore, the example configuration information shown in this section may not be identical to your software. Click | Pelp | located in the lower right corner of the window to view customized online help for your transformer.

Example Test Code Definitions Tab

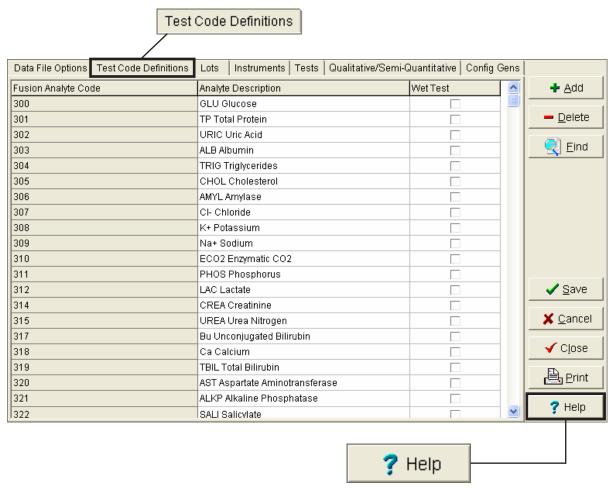

Example Test Code Definitions Tab (continued)

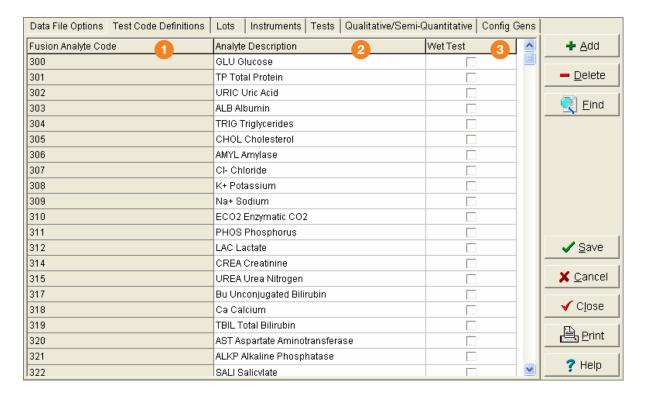

1 Analyte Code

The analyte designation in the QC data file.

2 Analyte Description

The instrument manufacturer's description of the analyte.

3 Wet Test


Select the check box in the Wet Test column to indicate a VITROS wet test.

Lots Tab


The Lots tab contains all of your specific lot information.



Important: Lot configuration in UnityConnect varies according to specific information in the QC data file. Therefore, the example configuration information shown in this section may not be identical to your software. Click located in the lower right corner of the window to view customized online help for your transformer.

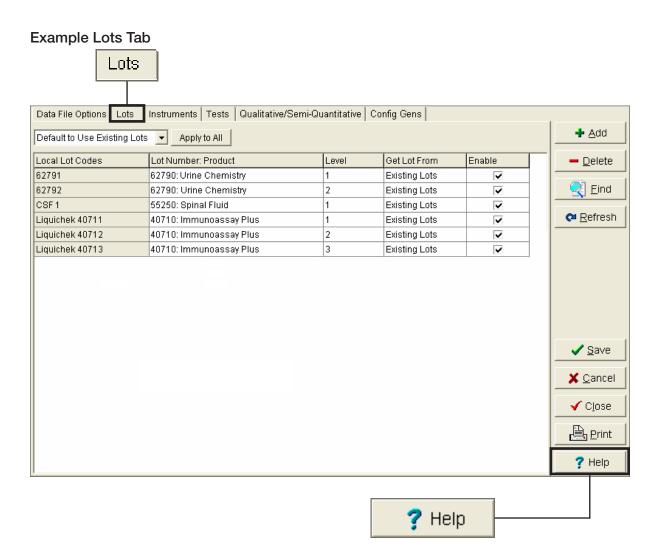


Example Lots Tab (continued)

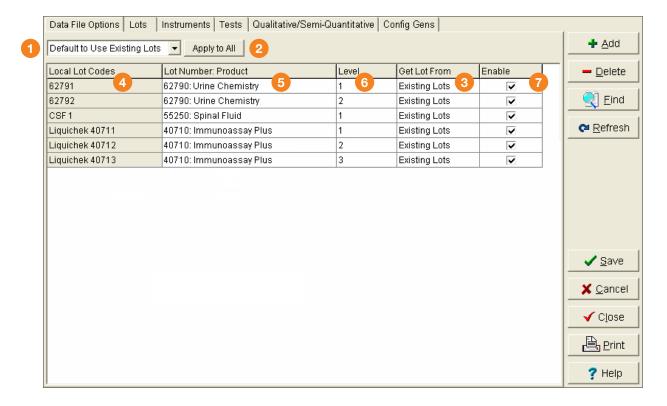

1 Database View Setting

Select an option from the Default to Use list to determine the default setting for the Get Lot From column . The setting in the Get Lot From column determines the information appearing in the Lot Number: Product list for an individual row.

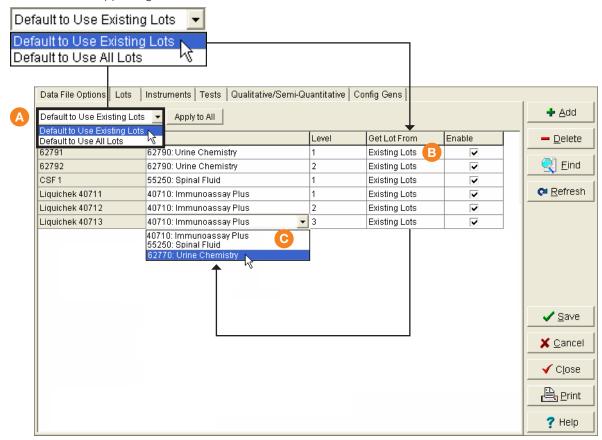

There are two options available:

Default to Use Existing Lots

Only lots set up in your QC data management software appear in the list in the **Lot Number: Product** column if this option is selected.

Default to Use All Lots

All lots in the Unity code list appear in the list in the **Lot Number: Product** column if this option is selected.




Tip: The process for configuring a new lot depends on your use of your QC data management software for SPC Rule Evaluation. See Chapter 10, "New Configuration in UnityConnect" for more information about new configuration of lots.


2 Apply to All

Click **Apply to All** to apply the **Default to Use** selection to all lots. Use the **Get Lot From** list to change the selection for individual rows if necessary (see 3).

Get Lot From


The default setting is determined by the **Database View Setting** selected (see 1 on page 34). Select another option from the list to change the information appearing in the **Lot Number: Product** list for an individual row.

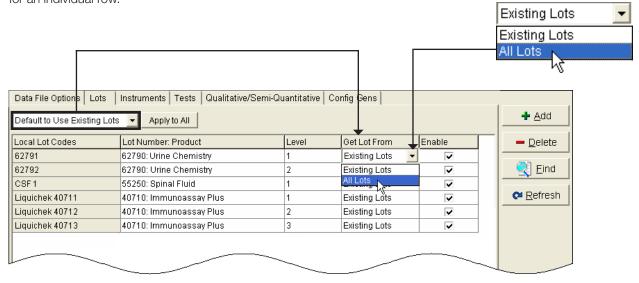


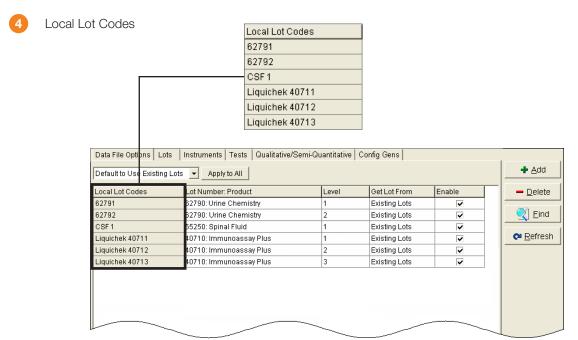

5 Lot Number: Product

Select the corresponding Bio-Rad master lot number and control product name from the list. The contents of the list are determined by the **Database View Setting** and/or the **Get Lot From** setting (see 3 on page 35)

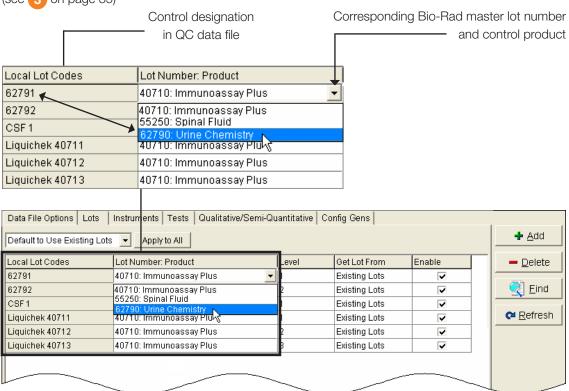

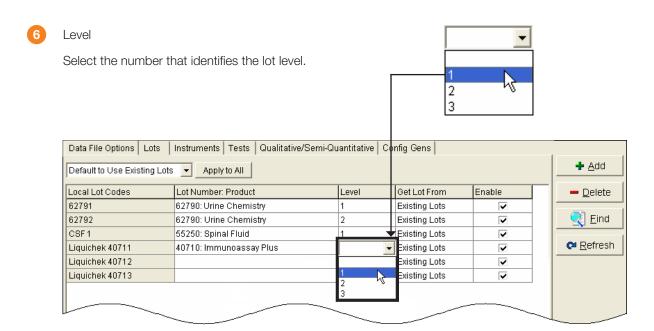


7 Enable

The **Enable** check box is selected by default so UnityConnect processes data for the lot/level.



Important: Click the **Enable** check box to clear it and UnityConnect will exclude the lot from current and future transforming. DO NOT delete the row. Deleting the row will cause the lot to appear again for future transforming.

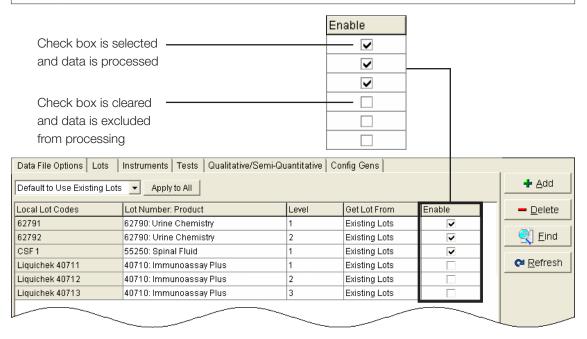

Special Cases Tab



Important: Due to the unique configuration of your customized UnityConnect software, talk to your Bio-Rad representative if you need to use the **Special Cases** tab.


Instruments Tab

The **Instruments** tab contains all of your specific instrument information.




Important: Instrument configuration in UnityConnect varies according to specific information in the QC data file. Therefore, the example configuration information shown in this section may not be identical to your software. Click Pelp located in the lower right corner of the window to view customized online help for your transformer.

Example Instruments Tab

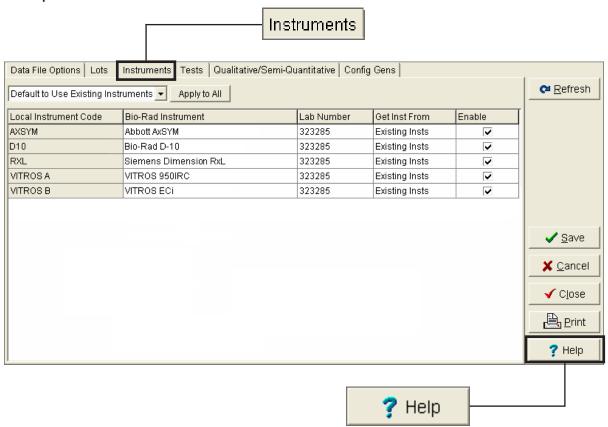


Example Instruments Tab (continued)

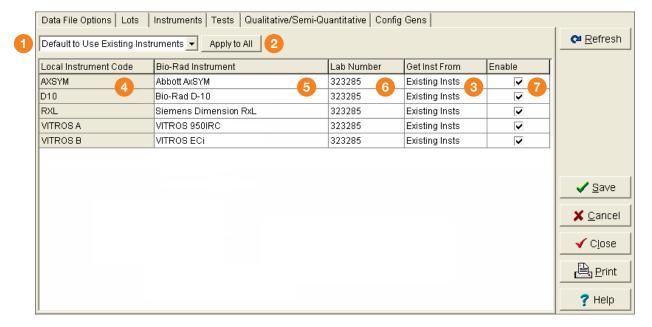

1 Database View Settings

Select an option from the **Default to Use** list odetermine the default setting for the **Get Inst From** column The setting in the **Get Inst From** column determines the information appearing in the **Bio-Rad Instruments** list for an individual row.

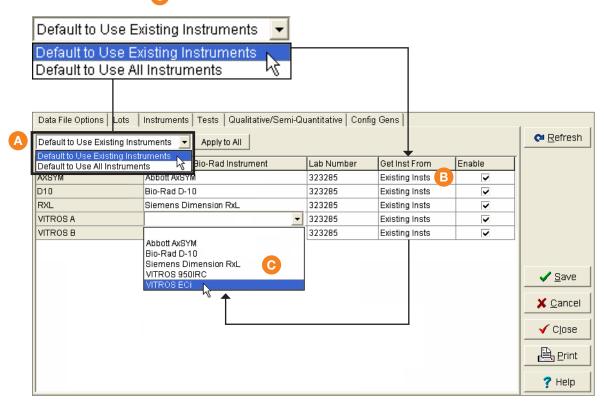

There are two options available:

Default to Use Existing Instruments


Only instruments set up in your QC data management software appear in the list in the **Bio-Rad Instrument** column if this option is selected.

Default to Use All Instruments

All instruments in the Unity code list appear in the list in the **Bio-Rad Instrument** column if this option is selected.


2 Apply to All


Click **Apply to All** to apply the **Default to Use** selection. Use the **Get Inst From** list to change the selection for individual rows if necessary (see 3 on page 41).

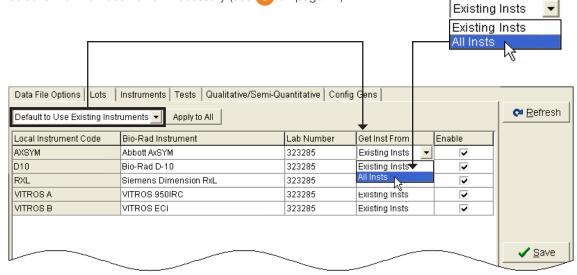

Get Inst From


The default setting in this list is determined by the **Database View Setting** selected (see 1 on page 40). Select another option from this list to change the information appearing in the **Bio-Rad Instrument** list for an individual row.

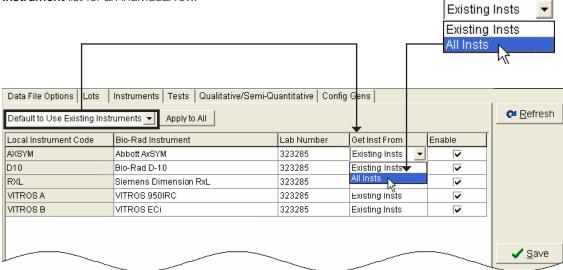

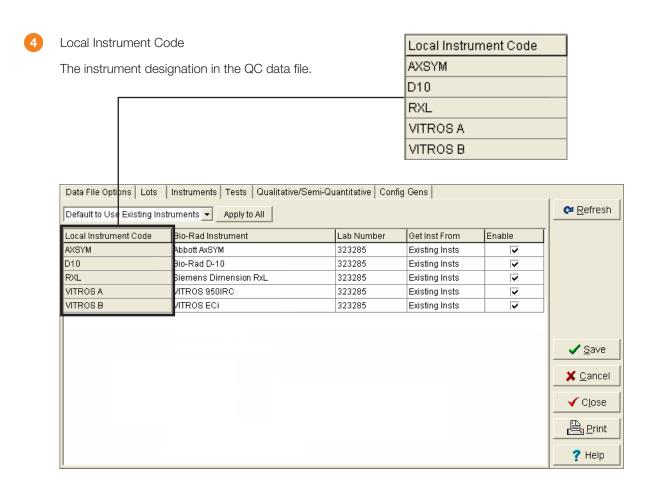

Bio-Rad Instrument


Select the corresponding Bio-Rad instrument name from the list. The contents of the list are determined by the Database View Settings and/or the Get Inst From setting (see 3 on page 41).

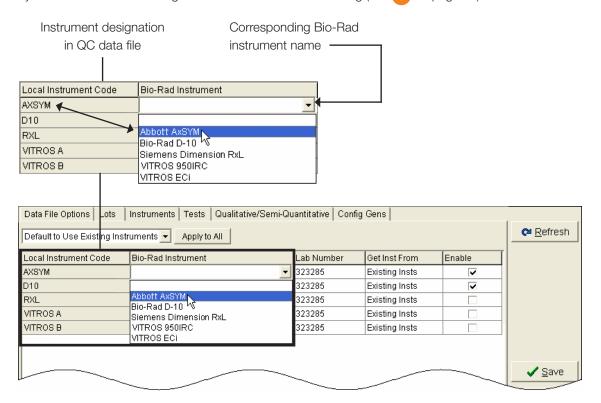

6

Lab Number

Select the Bio-Rad lab number assigned to the instrument. Only lab numbers set up in your QC data management software appear in the list.

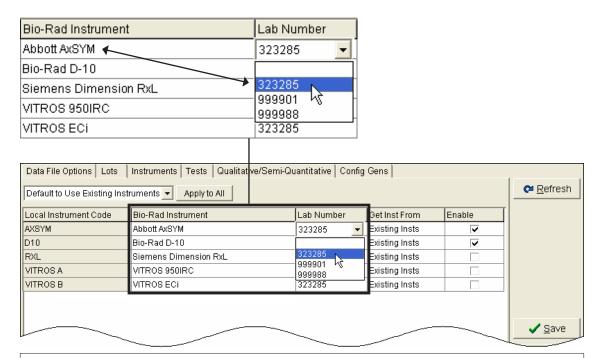



Important: Two instruments of the same type must be set up in separate lab numbers. A lab number can only be assigned by Bio-Rad. Contact your Bio-Rad QC Program Representative if you need additional lab numbers.




Enable

The **Enable** check box is selected by default so UnityConnect processes data for the instrument.



Important: Click the **Enable** check box to clear it and UnityConnect will exclude the instrument from current and future transforming. DO NOT delete the row. Deleting the row will cause the instrument to appear again for future transforming.

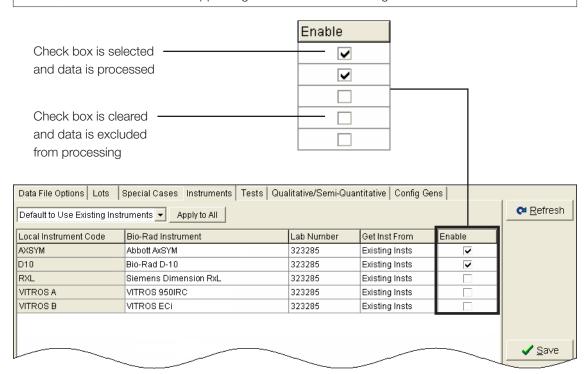


Tests Tab


The **Tests** tab contains all of your specific test information.



Important: Test configuration in UnityConnect varies according to specific information in the QC data file. Therefore, the example configuration information shown in this section may not be identical to your software. Click located in the lower right corner of the window to view customized online help for your transformer.

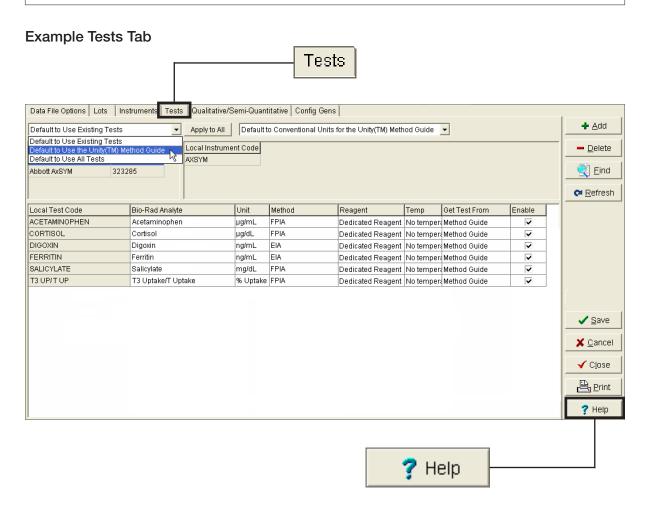


Example Tests Tab (continued)

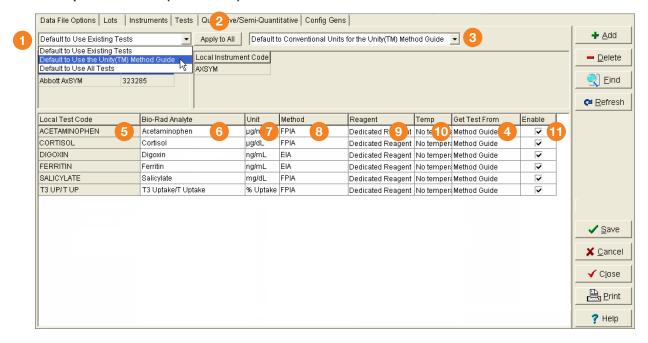

1 Database View Settings

Select an option from the **Default to Use** list A to determine the default setting for the **Get Test From** column B The setting in the **Get Test From** column determines the information appearing in the **Bio-Rad Analyte**, **Unit**, **Method**, **Reagent**, and **Temp** columns for an individual row.

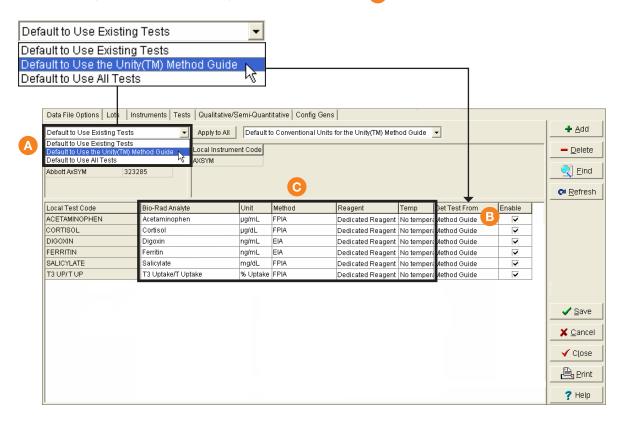

There are three options available:

Default to Use Existing Tests

Analyte information for tests set up in your QC data management software and run on the specified instrument appear in the lists in the **Bio-Rad Analyte**, **Unit**, **Method**, **Reagent**, and **Temp** columns if this option is selected.

Default to Use Method Guide

Only the analyte information for tests defined in the *Unity Method Guide for Selected Instruments* appear in the **Bio-Rad Analyte**, **Unit**, **Method**, **Reagent**, and **Temp** columns if this option is selected.

Default to Use All Tests


All analyte information in the Unity code list appear in the **Bio-Rad Analyte**, **Unit**, **Method**, **Reagent**, and **Temp** columns if this option is selected. This is the most useful method if tests are not defined in the *Unity Method Guide for Selected Instruments*.

2 Apply to All

Click Apply to All to apply the Database View Setting and/or the Default to Use Units selection(s). Use the Get Test From list to change the selection Database View Setting for individual rows if necessary (see 4 on page 49).


Operault to Use Units

Select the units (conventional or SI) to apply the units selection to all tests.

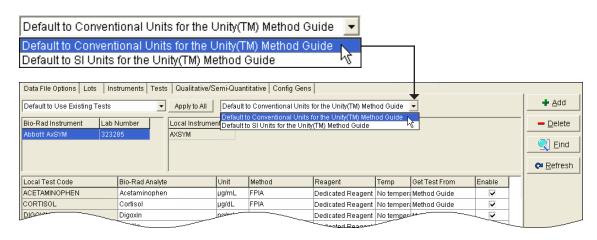


4 Get Test From

The default setting is determined by the **Database View Setting** selected (see 1 on page 48). Select another option from this list to change the information appearing in the **Bio-Rad Analyte**, **Unit**, **Method**, **Reagent**, and **Temp** lists for an individual row.

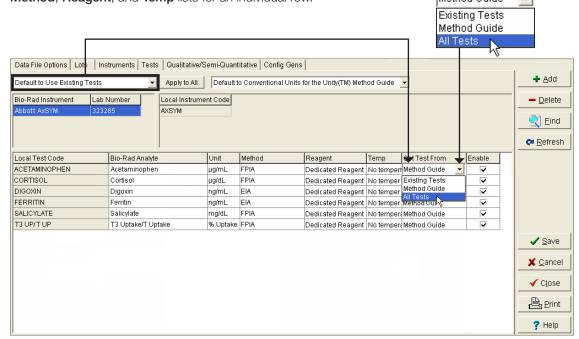


5 Local Test Code

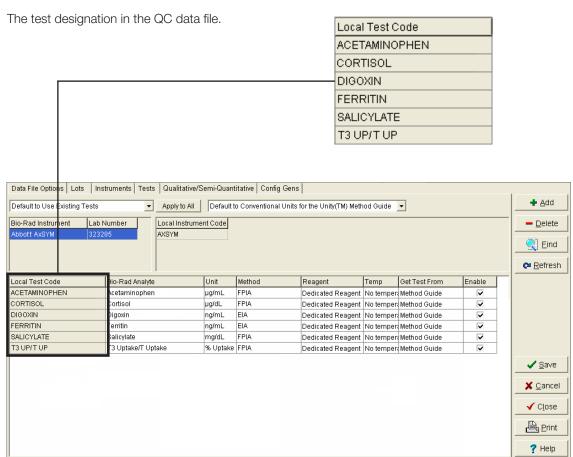

6 Bio-Rad Analyte

Select the corresponding Bio-Rad analyte name from the list. The contents of the list are determined by the **Database View Setting** and/or the **Get Test From** setting (see 4) on page 49).

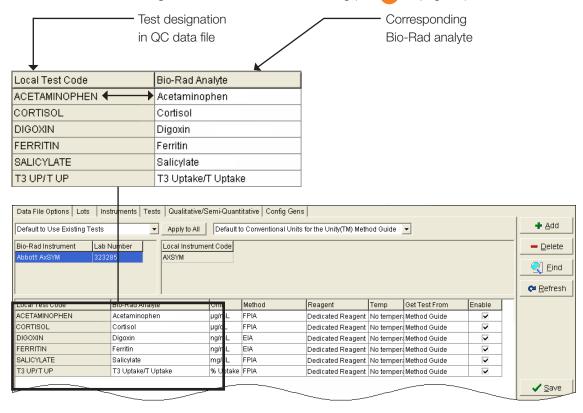

7 Unit

Select the corresponding unit from the list. The contents of the list are determined by the **Database**View Setting and/or the Get Test From setting (see 4) on page 49).

8 Method

Select the corresponding method from the list. The contents of the list are determined by the **Database View Setting** and/or the **Get Test From** setting (see 4 on page 49).

9 Reagent

Select the corresponding reagent from the list. The contents of the list are determined by the **Database**View Setting and/or the Get Test From setting (see 4) on page 49).

10 Temp

Select the corresponding temperature from the list. The contents of the list are determined by the **Database View Setting** and/or the **Get Test From** setting (see 4 on page 49).



Note: Temperature applies to enzymes only. For all other analytes, **No Temperature** is the only available option.




Enable

The Enable check box is selected by default so UnityConnect processes data for the test.



Important: Click the **Enable** check box to clear it and UnityConnect will exclude the lot from current and future transforming. DO NOT delete the row. Deleting the row will cause the lot to appear again for future transforming.

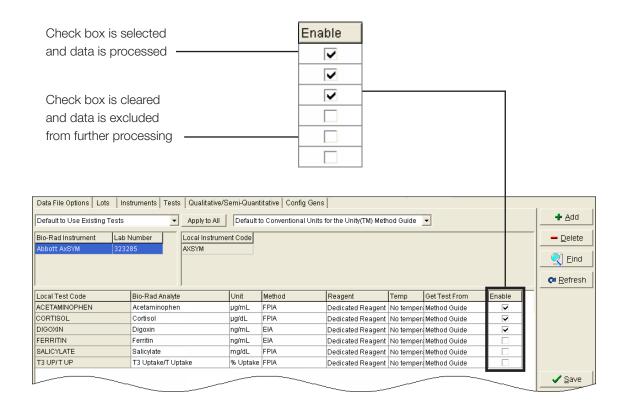

Qualitative/Semi-Quantitative Tab



Important: Due to the unique configuration of your customized UnityConnect software, talk to your Bio-Rad representative if you need to use the **Qualitative/Semi-Quantitative** tab.

Config Gens Tab



Important: The **Config Gens** tab only applies to customers with specific VITROS instruments. See Chapter 11, "VITROS Slide Generation Numbers and UnityConnect" on page 201 for more information about configuring slide generation numbers.

Working with QC Data Files

In This Chapter

Overview	53
QC Data File Content	54
QC Data File Essentials	
Creating a QC Data File	
Copying a QC Data File	
Copyling a QC Data File	59

Overview



Note: This chapter discusses QC data files used in the **data file transfer method** with UnityConnect.

A QC data file contains the quality control testing results performed on an instrument. The QC data file can be retrieved from an instrument, middleware, or a Laboratory Information System (LIS).


QC data management software cannot import QC data files in their original format. UnityConnect transforms the QC data file from the original format into the Bio-Rad format for automatic importing into your QC data management software.



Note: See Chapter 3, "Configuring UnityConnect for Data File Transfer with Desktop Software" for information about configuring UnityConnect for automatic importing. QC data management software includes desktop software (Unity Desktop and Unity Real Time), and Web-based software (UnityWeb 2.0 and Unity Real Time online).

QC Data File Content

The QC data file is the most critical component of the UnityConnect process. Although the look of QC data files varies according to the instrument, middleware, or LIS, the essential components are the same.

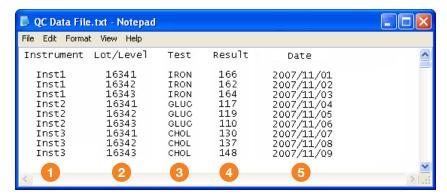

A proper QC data file contains unique designations for each component of the data file.

- Unique designation for each instrument

 Each instrument has a unique designation in the example above: Inst1, Inst2, and Inst3.
- Unique designation for the lot number and level

 Each lot number and level has a unique designation in the example above: 16341 (lot 16340-level 1), 16342 (lot 16340-level 2), 16343 (lot 16340-level 3).
- Unique test name
 Each test has a unique designation in the example above: IRON, GLUC, and CHOL.
- 4 Test result

Each result has a valid numeric value in the example above.



Note: UnityConnect cannot transform data with a 0 (zero) or negative numeric value.

Date of each test result

Each result has an associated date. The format used in the example above is year/month/day.

QC Data File Essentials

A properly formatted QC data file from your instrument, middleware, or LIS is essential to the UnityConnect process. A good data file allows UnityConnect to transform the QC data into a format that can be automatically imported into your QC data management software.



Note: See Chapter 3, "Configuring UnityConnect for Data File Transfer with Desktop Software" for information about configuring UnityConnect for automatic importing.

- The QC data file must be comprehensive and contain a unique designation for each instrument as well as all lots, levels, and tests run on Bio-Rad controls.
- The data file must include a value and a date for all QC data to be submitted to Bio-Rad for peer group comparison.
- The QC data file must be in a stable format that can be reproduced as often as necessary. You must create a new QC data file every time you use UnityConnect to transform and import the QC data into your QC data management software.



Important: Bio-Rad strongly recommends carefully documenting the steps taken to create the QC data file as this will become part of your regular QC process.

Maintaining QC Data File Integrity

Maintaining the integrity of the QC data file is very important. Never open the original QC data file to be used for transforming data. In some cases, opening a file alters the formatting and causes the QC data file to become useless.



Important: Do not open or manually alter the data in the original QC data file. For example, manually deleting outliers can corrupt the file.

If it is necessary to open a QC data file, make a copy of the file and open the copy. See "Copying a QC Data File" on page 59 for more information.

Multiple QC Data Files

For ease of use and simplicity, Bio-Rad recommends creating one QC data file containing all QC data. However, some instruments, middleware, and LIS separate QC data into multiple files. For example, customers using an Orchard LIS or a Beckman Coulter Access instrument are required to create multiple QC data files for different control products, levels, or different tests. See the following section according to the QC data management software you use:

- Desktop software
 - See "Working with Multiple QC Data Files" on page 148.
- Web-based software

See "Working with Multiple QC Data Files" on page 159.

Single Point Data versus Summary Data

UnityConnect can process single point data or summary data.

Single Point Data

Single point data files contain an individual value for each QC result. Single point data is typically preferred since it facilitates the use of advanced features of your QC data management software such as Levey-Jennings Charts, SPC (Westgard) rules, and the optional Westgard Advisor.

Summary Data


Summary data files contain the mean, standard deviation (SD), and number of points collected for each test over one month. Using summary data limits the availability of advanced features of your QC data management software.

Unfiltered Data versus Filtered Data


UnityConnect can process unfiltered QC data or filtered QC data. It is important to understand the difference between the two types of data.

Unfiltered Data

A QC data file is considered unfiltered when no outliers have been removed or excluded in the instrument, middleware, or LIS. The following illustration shows an example of an unfiltered QC data file.

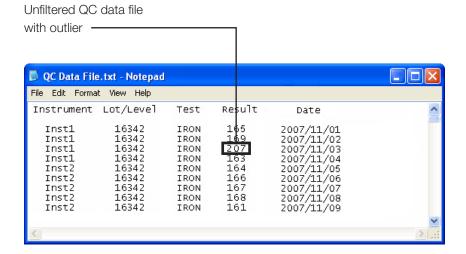

It is important to understand how your instrument, middleware, or LIS handles outliers in the QC data file. For example, the next illustration shows the same outlier, but in this case the LIS recognized the outlier and required a comment be added.

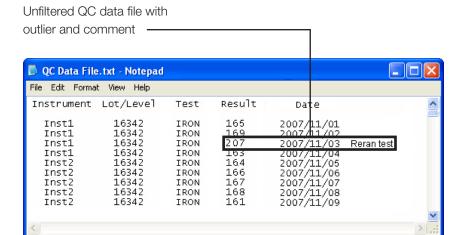

Although the LIS required a comment be added to the outlier, the outlier was not removed and therefore the QC data file is still considered unfiltered. To determine if a QC data file is unfiltered, make a copy of the data file, open the file, and search for a known outlier. (See "Copying a QC Data File" on page 59 for more information.)



Important: Do not open or manually alter the data in the QC data file. For example, manually deleting outliers can corrupt the file.


Filtered Data

A QC data file is considered filtered when the outliers have been removed or excluded in the instrument, middleware, or LIS. The following example shows a comparison between an unfiltered QC data file and a filtered QC data file.

QC data file—unfiltered showing outlier and comment QC Data File.txt - Notepad File Edit Format View Help Instrument Lot/Level Test Result Date Inst1 16342 IRON 165 2007/11/01 2007/11/02 16342 IRON Inst1 16342 IRON 207 2007/11/03 Reran test Inst1 Inst1 16342 IRON 16342 IRON 164 Inst2 2007/11/05 16342 Inst2 IRON 166 2007/11/06 2007/11/07 2007/11/08 Inst2 16342 IRON 167 Inst2 IRON 168 Inst2 16342 IRON 161 2007/11/09

Same QC data file-filtered

The outlier (207) and comment have been removed or excluded in the instrument, middleware, or LIS.

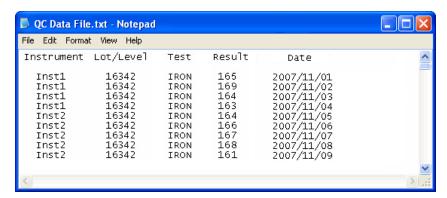

Creating a QC Data File

The steps for creating a data file vary according to the instrument, middleware, or LIS. Bio-Rad has directions for creating QC data files for most of the common instruments, middleware, and LIS.



Important: If Bio-Rad does not have directions for creating a QC data file on your instrument, middleware, or LIS, contact your laboratory's instrument, middleware, or LIS specialist or contact the instrument, middleware, or LIS manufacturer.

How Often to Create a QC Data File

You must create a new QC data file every time you use UnityConnect to transform and import QC data into your QC data management software. This must be done at least on a monthly basis in order to submit the QC data to Bio-Rad for peer group comparison and Unity Interlaboratory Reports.

Some customers find that 30 days of data is a very large file. In such cases, a QC data file can be created, transformed, and imported into your QC data management software more frequently as needed.




Best Practices - Transform Data Often

Bio-Rad recommends transforming data as often as possible throughout the month. This helps you become more familiar with the transformation process and alleviates the stress and pressure at the end of the month.


Copying a QC Data File

1 Right click on the original file and then click Copy. -

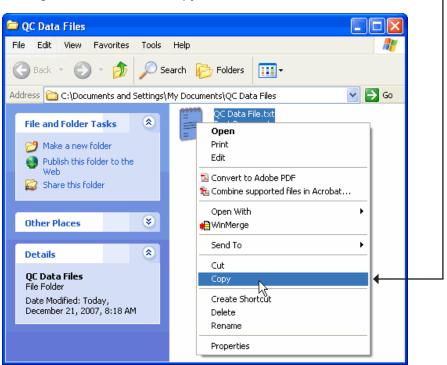


2 Move the mouse away from the original data file.

3 Right click in the window and then click Paste.

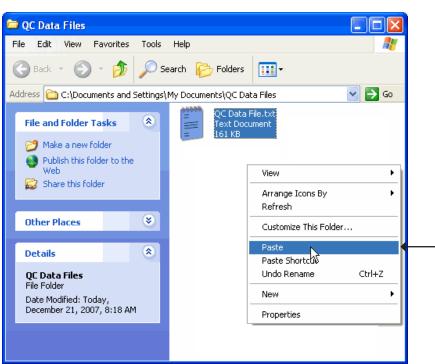

A copy of the data file appears in the window.

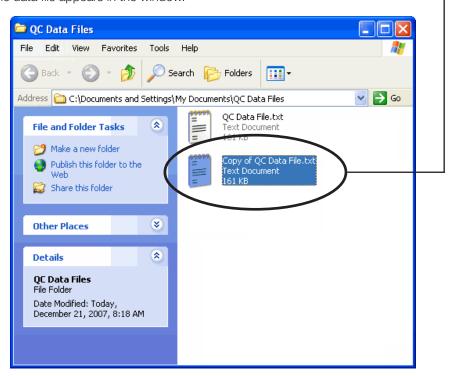

4 The copy of the QC data file can be opened and reviewed.



Important: Only use an original QC data file that has not been opened to transform your QC data.

Configuring UnityConnect for Data File Transfer with Desktop Software

In This Chapter

Before You Begin	62
Configuration Flow Chart	
Output Path	
Configuring Watch Folder and Automatic Import	
Configuring Import Options	
Downloading a Transformer	
Configuring E-Mail Notification	74
Configuring the Default Date Settings	75


Before You Begin



This chapter is only for customers using UnityConnect with the data file transfer method and any of the following desktop software:

Unity Desktop
Unity Real Time

Configuration Flow Chart

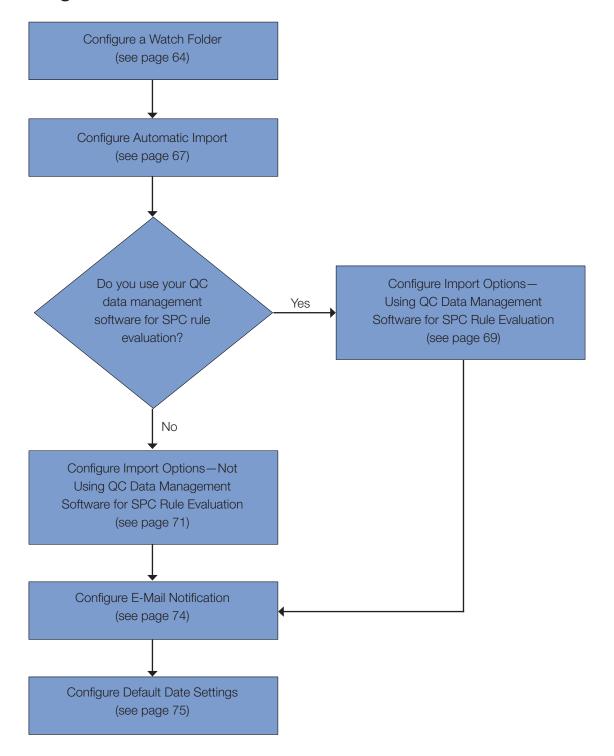

Output Path

Each installed and licensed transformer in UnityConnect requires an output path. The output path is the location where your transformed QC data file is stored.

In most instances, the output path is the same location for all transformers. However, depending on your laboratory set up, multiple output paths may be useful. For example, if your laboratory uses multiple transformers, it is useful to have a separate output path for each transformer. Consult Bio-Rad for more information.



Important: The output path is pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, do not change the output path. Talk to your Bio-Rad representative for more information about output paths.

Configuring Watch Folder and Automatic Import

UnityConnect can be configured to "watch" for data entering a "Watch Folder" and then automatically import the data into your QC data management software.


Configuring a Watch Folder

- 1 Create a new folder on the computer desktop where UnityConnect is installed.
- 2 Name the folder "Watch Folder."




Note: A Watch Folder must be created for each transformer. Give each Watch Folder a unique name if you create more than one Watch Folder. For example, "Watch Folder Meditech" and "Watch Folder RXL." This helps easily identify the Watch Folder for each transformer.


- 3 Start UnityConnect.
- 4 Click the **Options** menu and then click **Settings**.

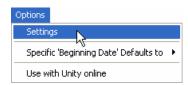

5 Click the Watch Folder tab.

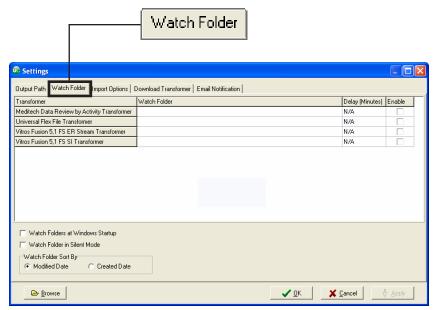


6 Click in the Watch Folder field for the transformer you want to configure.

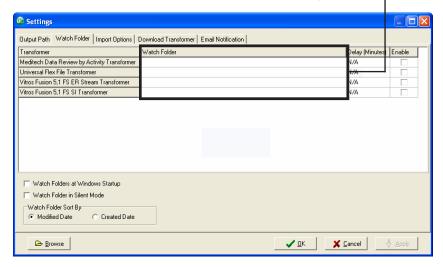

7 Click Browse

- 8 Select the Watch Folder you created on your computer desktop.
- 9 Click OK .
- 10 Select the **Enable** check box for the transformer.

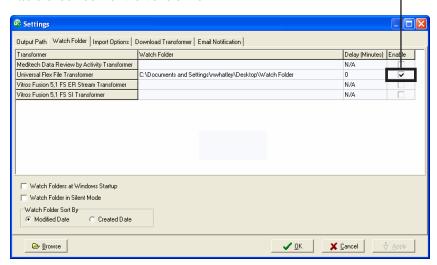

11 Select the **Watch Folders at Windows Startup** check box if you want the Watch Folder to be active each time the computer starts up.



Note: You must click on the UnityConnect main window to make the Watch Folder active if you do not select the **Watch Folders at Windows Startup** check box.

12 Select the Watch Folder in Silent Mode check box if you want the Watch Folder to run in silent mode.



Note: Silent mode processes all data in the QC data file that UnityConnect recognizes (that is, all data that has been configured in UnityConnect). If UnityConnect does not recognize any item(s) in the QC data file, a copy of the original QC data is saved for later reprocessing when the unrecognized item(s) can be configured. UnityConnect can also be configured to send an e-mail notification when there are QC data files requiring reprocessing. See "Configuring E-Mail Notification" on page 74 for more information.


- 13 Select a Watch Folder Sort By option:
 - Modified Date
 - Created Date
- 14 Click ♦ Apply
- 15 Continue with the next section, "Configuring Automatic Import."

Configuring Automatic Import




Important: Make sure you configure a Watch Folder before configuring automatic import. See "Configuring a Watch Folder" on page 64 for more information.

1 Click the **Import Options** tab.

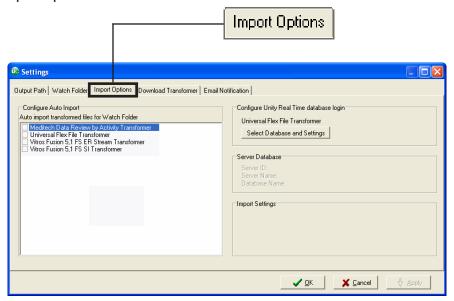

2 Select the check box for the transformer you want to set up for automatic importing.

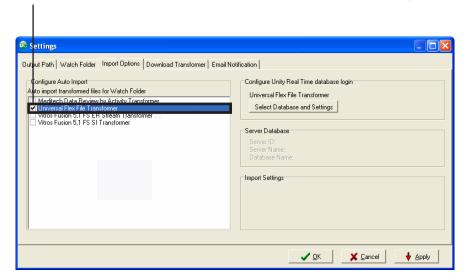

- **3** Click the appropriate button:
 - Click v <u>DK</u> to apply the selected settings and close the **Settings** dialog box.
 - Click to apply the selected settings without closing the Settings dialog box. This is convenient if you want to continue configuring settings on other tabs.


Configuring Import Options

Transformed QC data is imported into your QC data management software based on the import options defined in UnityConnect. Select import options based on your use of your QC data management software for SPC rule evaluation.



Note: The import options only need to be configured one time. You do not need to configure the settings each time you transform data.

Use the import options described in the appropriate section according to how you use your QC data management software for SPC rule evaluation.

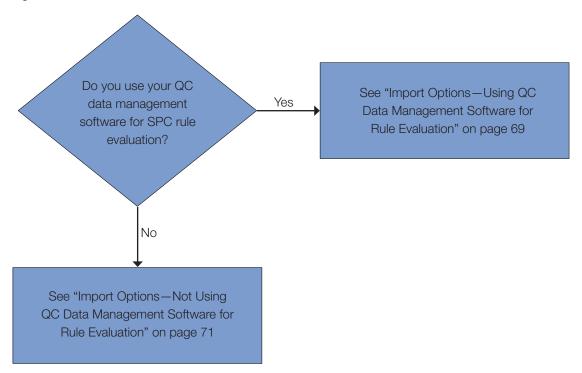


Import Options—Using QC Data Management Software for SPC Rule Evaluation


1 Click the **Options** menu and then click **Settings**.



2 Click the **Import Options** tab.

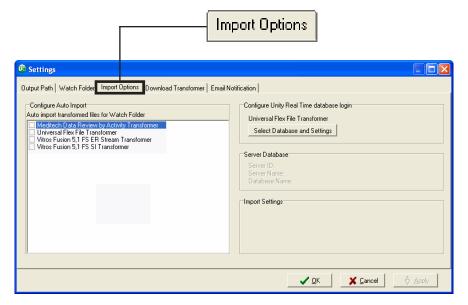

3 Select the transformer for which you want to configure import options.

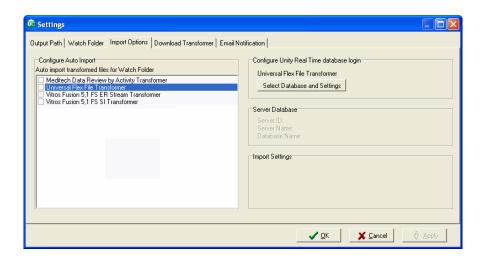



Note: Do not select the check box for the transformer unless you want to set up automatic import. See "Configuring Automatic Import" on page 67 for more information.

Output Path | Watch Folder | Import Options | Download Transformer | Email Notification |

Configure Auto Import | Configure Auto Import International Review by Activity Transformer | Universal Flex File Transformer | Universal Flex File Transformer | Universal Flex File Transformer | Vitros Fusion 5,1 FS SI Transformer | Vitros Fusion 5,1 FS SI Transformer | Output Path | Vitros Fusion 5,1 FS SI Transformer | Vitros Fusion 5,1 FS SI Transformer | Output Path | Vitros Fusion 5,1 FS SI Transformer | Vitros Fusion 5,1 FS SI Transformer | Output Path | Vitros Fusion 5,1 FS SI Transformer | Vitros Fusion 5,1 FS SI Transformer | Output Path | Vitros Fusion 5,1 FS SI Transformer | Vitros Fusion 5,1 FS SI Transfor

5 Make sure the Create New Lots if Necessary check box is not selected.

The following illustration shows recommended import settings for customers using their QC data management software for SPC rule evaluation.

✓ <u>0</u>K

X Cancel

The Create New Lots if Necessary check box is not selected. 👵 Unity Real Time Database Login Import Settings Transformer. Universal Flex File Transformer Decimal Point: Run Length: 15 \$ Minute(s) SQL Server Database Connect using Windows Authentication Create New Lots if Necessary Create New Tests if Necessary SQL Server ID: ✓ Sort Import Files SQL Server Password: Redundant Data Filter Disable All Rules Violation ▼ SQL Server Name: Delete Rejection Log Records after 60 🕏 Day(s) • Database Name: Operator ▼ User Name: User Password: **✓** <u>0</u>K 💢 <u>C</u>ancel


This ensures UnityConnect will not create news lots when transforming and overwrite your current SPC rule settings in your QC data management software.

6 Click ✓ <u>0</u>K .


To duplicate your current SPC rule settings from the current lot to the new lot, leave the **Create New Lots if Necessary** check box blank and duplicate the lot in your QC data management software before transforming data for the lot for the first time. See "Configuring a New Lot" on page 164 for more information.

Import Options—Not Using QC Data Management Software for SPC Rule Evaluation


1 Click the **Options** menu and then click **Settings**.

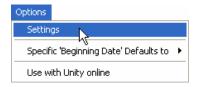

2 Click the **Import Options** tab.

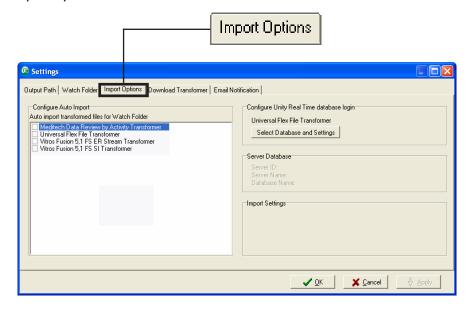

3 Select the transformer for which you want to configure import options.

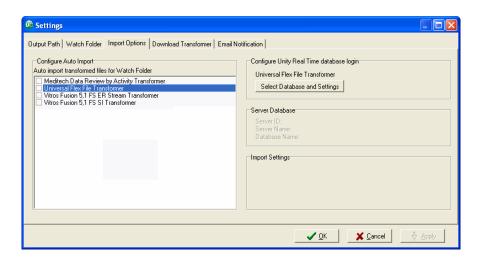

Important: Do not select the check box for the transformer unless you want to set up automatic import. See "Configuring Automatic Import" on page 67 for more information.

4 Click Select Database and Settings

Output Path Watch Folder Import Options Download Transformer Email Notification

Configure Auto Import
Auto import transformed files for Watch Folder

Meditech Data Review by Activity Transformer
Universal Flex File Transformer
Select Database and Settings

Server Database
Server D.
Server Database
Server D.
Server Name.
Database Name.


Import Settings

Import Settings

5 Make sure the Create New Lots if Necessary check box is selected.

The following illustration shows recommended import settings for customers not using their QC data management software for SPC rule evaluation.

Create New Lots if Necessary check box is selected.

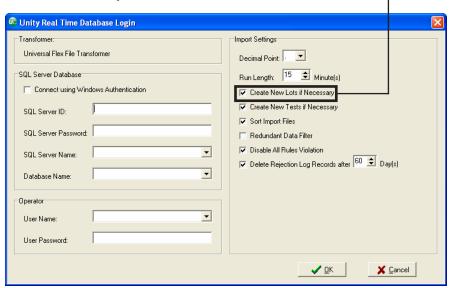

6 Click VOK

A new lot is automatically created when UnityConnect finds a new lot in your QC data file.

Downloading a Transformer

There are many instrument, middleware, and LIS transformers available for UnityConnect that you can download and install via the Internet.



Important: Talk to your Bio-Rad representative for more information about licensing and downloading a transformer.


Configuring E-Mail Notification

UnityConnect has the ability to send e-mail notifications to one or more e-mail addresses when:

- There are QC data files requiring reprocessing.
- There is an error during automatic import.
- If a Rejection Log has been created in UnityConnect.
- 1 Start UnityConnect.
- 2 Point to the Options menu and then click Settings.



3 Click the **Email Notification** tab.

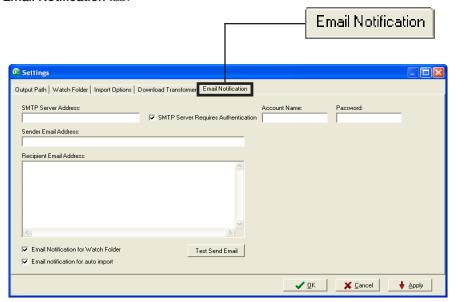

- 4 Type your SMTP server address in the **SMTP Server Address** field.
- 5 Select the **SMTP Server Requires Authentication** check box if your server requires authentication to send e-mail messages.
- Type the account name in the **Account Name** field if you selected the **SMTP Server Requires Authentication** check box.
- 7 Type the password in the Password field if you selected the **SMTP Server Requires Authentication** check hox
- 8 Type a sender e-mail address in the Sender Email Address field.



Tip: This helps prevent the sender e-mail address from being filtered as spam by your e-mail server.

- 9 Type each e-mail address you want to receive notification in the **Recipient Email Address** field. Type each e-mail address on a separate line.
- 10 Click Test Send Email to test the e-mail settings.

A test e-mail is sent to each recipient.

- 11 Select the **Email Notification for Watch Folder** check box if you want recipients to receive an e-mail when there is a file waiting for reprocessing.
- 12 Select the **Email notification for auto import** if you want recipients to receive an e-mail when an error or rejection occurs during import.
- **13** Click the appropriate button:
 - Click to apply the selected settings and close the Settings dialog box.
 - Click to apply the selected settings without closing the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.

Configuring the Default Date Settings

You can customize the default beginning date for the QC data file to transform. The beginning date you select is automatically set each time you open UnityConnect to transform your QC data file.



Note: You can always manually change the date range for the QC data file to transform. QC data can be transformed for a particular date range, from the first point in the QC data file to the last point in the QC data file, or a combination of both.

Click the Options menu, point to Specific Beginning Date Defaults to, and then click the option you want:

- First Day of Current Month
- Current Date



Configuring UnityConnect for Data File Transfer with Web-based Software

In This Chapter

Before You Begin	76
Configuration Flow Chart	77
Setting Up UnityConnect with Web-based Software	
Output Path	78
Configuring a Watch Folder	79
Configuring Internet Options	82
Configuring Auto Post	84
Downloading the Database	
Downloading Lab Setup and Code Lists	
Downloading a Transformer	
Configuring E-Mail Notification	
Configuring the Default Date Settings	90
Configuring Import Settings	91

Before You Begin




This chapter is only for customers using UnityConnect with the data file transfer method and any of the following desktop software:

UnityWeb

Unity Real Time online

Configuration Flow Chart

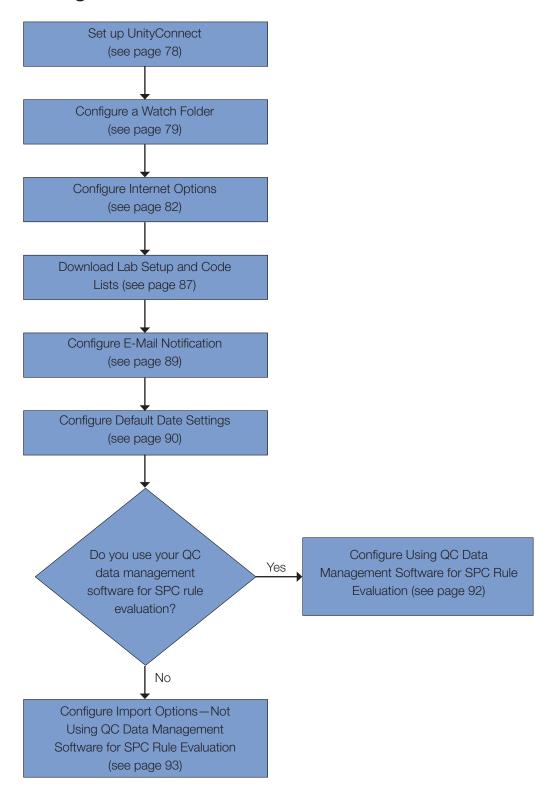

Setting Up UnityConnect with Web-based Software

- 1 Start UnityConnect.
- 2 Click the Options menu and then click Use with Unity online.



A check mark appears to the left indicating the **Use with Unity online** option is selected.



Output Path

Each installed and licensed transformer in UnityConnect requires an output path. The output path is the location where your transformed QC data file is stored.

In most instances, the output path is the same location for all transformers. However, depending on your laboratory set up, multiple output paths may be useful. For example, if your laboratory uses multiple transformers, it is useful to have a separate output path for each transformer. Consult Bio-Rad for more information.



Important: The output path is pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, do not change the output path. Talk to your Bio-Rad representative for more information about output paths.

Configuring a Watch Folder


UnityConnect can be configured to "watch" for data entering a "Watch Folder."

- 1 Create a new folder on the computer desktop where UnityConnect is installed.
- 2 Name the folder "Watch Folder."




Note: A Watch Folder must be created for each transformer. Give each Watch Folder a unique name if you create more than one Watch Folder. For example, "Watch Folder Meditech" and "Watch Folder RXL." This helps easily identify the Watch Folder for each transformer.


- 3 Start UnityConnect.
- 4 Click the **Options** menu and then click **Settings**.

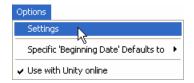

5 Click the Watch Folder tab.

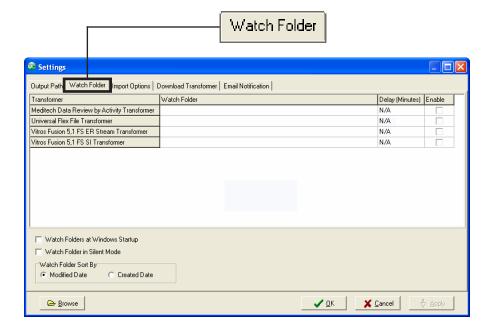

6 Click in the Watch Folder field for the transformer you want to configure.

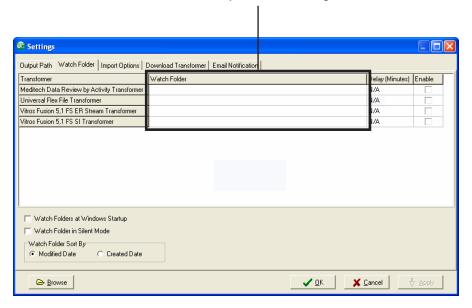

- 7 Click Browse .
- 8 Select the Watch Folder you created on your computer desktop.



9 Click OK

Settings Output Path Watch Folder Import Options Download Transformer Email Notification Transformer Watch Folder Delay (Minutes) Ena Meditech Data Review by Activity Transformer C:\Documents and Settings\nwhatley\Desktop\Watch Folder 0 Universal Flex File Transformer Vitros Fusion 5,1 FS ER Stream Transformer N/A N/A Vitros Fusion 5.1 FS SI Transformer Watch Folders at Windows Startup Watch Folder in Silent Mode Watch Folder Sort By Modified Date C Created Date <u>B</u>rowse **✓** <u>o</u>K X Cancel

10 Select the Enable check box for the transformer.

11 Select the **Watch Folders at Windows Startup** check box if you want the Watch Folder to be active each time the computer starts up.



Note: You must click on the UnityConnect main window to make the Watch Folder active if

you do not select the Watch Folders at Windows Startup check box.

12 Select the Watch Folder in Silent Mode check box if you want the Watch Folder to run in silent mode.



Note: Silent mode processes all data in the QC data file that UnityConnect recognizes (that is, all data that has been configured in UnityConnect). If UnityConnect does not recognize any item(s) in the QC data file, a copy of the original QC data is saved for later reprocessing when the unrecognized item(s) can be configured. UnityConnect can also be configured to send an e-mail notification when there are QC data files requiring reprocessing. See "Configuring E-Mail Notification" on page 89 for more information.

- 13 Select a Watch Folder Sort By option:
- Modified Date
- Created Date
- **14** Click the appropriate button:

 - Click to apply the selected settings and leave the **Settings** dialog box option. This is convenient if you want to continue configuring settings on other tabs.

Configuring Internet Options




Important: The Internet options are pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, talk to your Bio-Rad representative if you need further assistance configuring Internet options.

- 1 Start UnityConnect.
- 2 Click the **Options** menu and then click **Settings**.



3 Click the **Internet Options** tab.

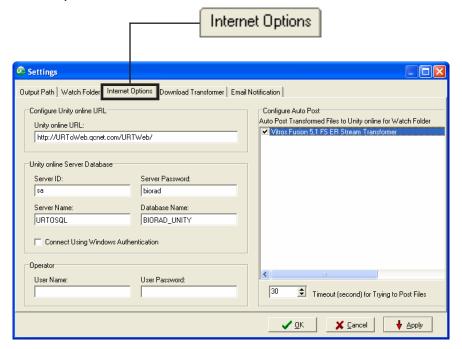

4 Verify the settings in the following fields:

Field	Setting
Unity online URL:	http://URToWeb.qcnet.com/URTWeb/
Server ID:	sa
Server Password:	biorad
Server Name:	URTOSQL
Database Name:	BIORAD_UNITY




Note: The settings shown above are the suggested settings at the time of this publication.

- 5 Type the user name you use when logging on to Unity Real Time online or UnityWeb 2.0 in the **User Name** field.
- Type the password you use when logging on to Unity Real Time online or UnityWeb 2.0 in the **User Password** field.
- 7 Click the appropriate button:
 - Click of to apply the selected settings and close the **Settings** dialog box.
 - Click to apply the selected settings without leaving the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.

Testing the Internet Connection


You can test the Internet connection after configuration to ensure the configuration is correct.

1 Click the File menu and then click Download Database.

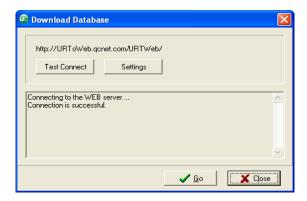

2 Click Test Connect

A message appears when the connection is successful.

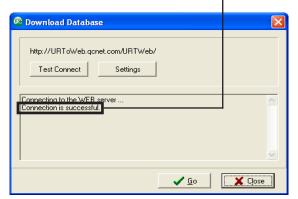

3 Click X Close .

If the Internet Connection Fails

- 1 Click Settings .
- 2 Review the configuration and make any changes as required.
- 3 Repeat the testing as described on the previous page.
- 4 Contact Bio-Rad Software Support for assistance if the testing continues to fail.

Configuring Auto Post




Note: Auto post automatically uploads your transformed QC data to Unity Real Time online or UnityWeb 2.0.


- 1 Start UnityConnect.
- 2 Click the Options menu and then click Settings.



3 Click the Internet Options tab.

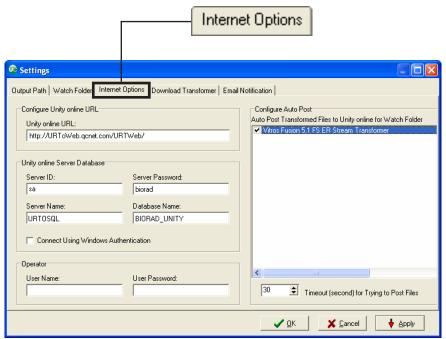

4 Select the check box for the transformer you want to configure for auto post. -

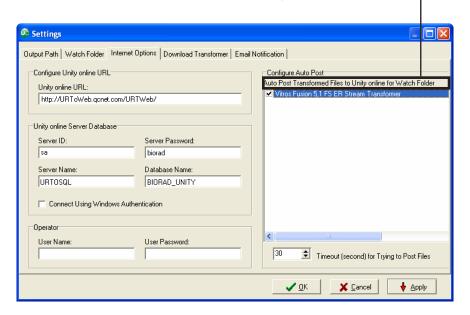



Note: The settings shown in the example above are the suggested settings at the time of this publication.

- 5 Type the user name you use when logging on to Unity Real Time online in the **User Name** field.
- 6 Type the password you use when logging on to Unity Real Time online in the **User Password** field.

- 7 Click the appropriate button:
 - Click V DK to apply the selected settings and close the **Settings** dialog box.
 - Click to apply the selected settings without leaving the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.


Downloading the Database



Note: The database only needs to be download one time. See "Downloading Lab Setup and Code Lists" on page 87 for information about keeping the database up-to-date.

- 1 Click the File menu and then click Download Database.
- 2 Click 🗸 🗓 .

A message appears when the download completes 3 successfully.

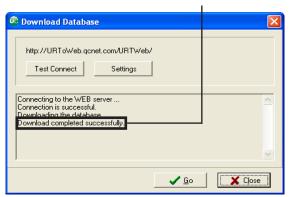

3 Click Close

Downloading Lab Setup and Code Lists




Important: Perform the steps in this section any time you make a configuration change to your UnityWeb 2.0 or Unity Real Time online software. In addition, Bio-Rad recommends performing these steps at least monthly to ensure your software always has the most up-to-date codes lists and other information.

1 Click the File menu and then click Download Lab Setup/Code Lists.




- 2 Select the check box for each item you want to download:
 - Download Lab Setup
 - Download Code Lists

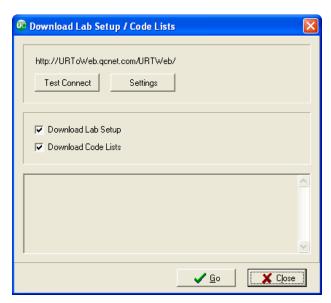

3 Click ✓ Go

A message appears when the download completes successfully.

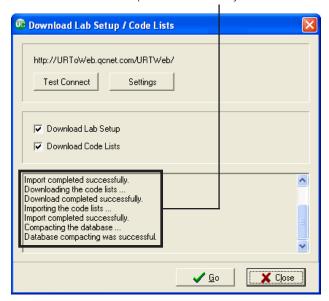

4 Click Close .

Downloading a Transformer

There are many instrument, middleware, and LIS transformers available for UnityConnect that you can download and install via the Internet.




Important: Talk to your Bio-Rad representative for more information about licensing and downloading a transformer.

Configuring E-Mail Notification

UnityConnect has the ability to send e-mail notifications to one or more e-mail addresses when:

- There are QC data files requiring reprocessing.
- There is an error during automatic upload.
- If a Rejection Log has been created in UnityConnect.
- Start UnityConnect.
- 2 Point to the **Options** menu and then click **Settings**.



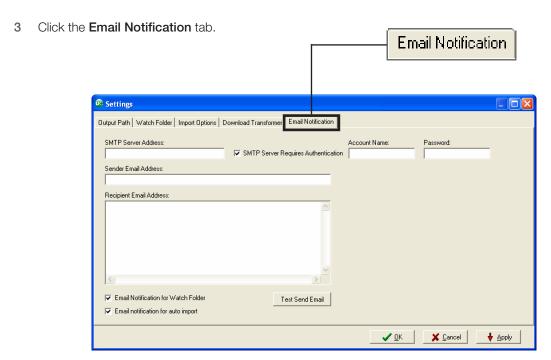

- 4 Type your SMTP server address in the **SMTP Server Address** field.
- 5 Select the **SMTP Server Requires Authentication** check box if your server requires authentication to send e-mail messages.
- Type the account name in the **Account Name** field if you selected the **SMTP Server Requires Authentication** check box.
- 7 Type the password in the **Password** field if you selected the **SMTP Server Requires Authentication** check box.

8 Type a sender e-mail address in the **Sender Email Address** field.



Tip: This helps prevent the sender e-mail address from being filtered as spam by your e-mail server.

- 9 Type each e-mail address you want to receive notification in the **Recipient Email Address** field. Type each e-mail address on a separate line.
- 10 Click Test Send Email to test the e-mail settings.

A test e-mail is sent to each recipient.

- 11 Select the **Email Notification for Watch Folder** check box if you want recipients to receive an e-mail when a file is waiting for reprocessing.
- 12 Select the **Email Notification for Auto Post** if you want recipients to receive an e-mail when an error has occurred during auto-post.
- 13 Click the appropriate button:
 - Click ✓ <u>□</u>K to apply the selected settings and close the **Settings** dialog box.
 - Click Apply to apply the selected settings and leave the **Settings** dialog box option. This is convenient if you want to continue configuring settings on other tabs.

Configuring the Default Date Settings

You can customize the default beginning date for the QC data file to transform. The beginning date you select is automatically set each time you open UnityConnect to transform your QC data file.



Note: You can always manually change the date range for the QC data file to transform. QC data can be transformed for a particular date range, from the first point in the QC data file to the last point in the QC data file, or a combination of both.

Click the Options menu, point to Specific Beginning Date Defaults to, and then click the option you want:

- First Day of Current Month
- Current Date




Configuring Import Settings

Transformed QC data is imported into your QC data management software based on the import settings defined in the software. Select import settings based on your use of your QC data management software for SPC rule evaluation.



Note: The import settings only need to be configured one time. You do not need to configure the settings each time you transform data.

Use the import settings described in the appropriate section according to how you use your QC data management software for SPC rule evaluation.

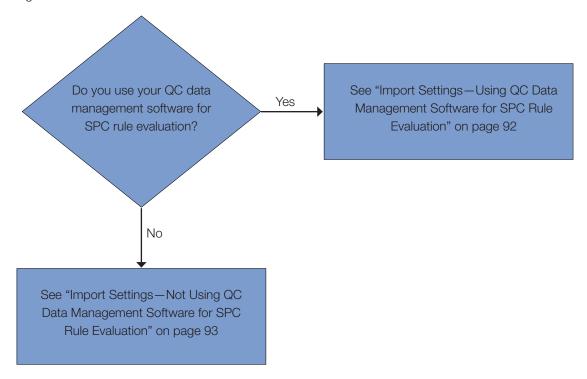


Import Settings—Using QC Data Management Software for SPC Evaluation

These settings ensure UnityConnect will not create new lots when transforming and overwrite your current SPC rule settings in your QC data management software.

To duplicate your current SPC rule settings from the current lot to the new lot, leave the **Create new Lots** check box blank and duplicate the lot in your QC data management software before transforming data for the lot for the first time. See "Configuring a New Lot" on page 164 for more information.

- 1 Log in to the UnityWeb 2.0 or Unity Real Time online software.
- 2 Click the Tools tab.
- 3 Click Utilities.
- 4 Click Import Settings.


The following illustration shows recommended import settings for customers using UnityWeb 2.0 or Unity Real Time online for SPC rule evaluation.

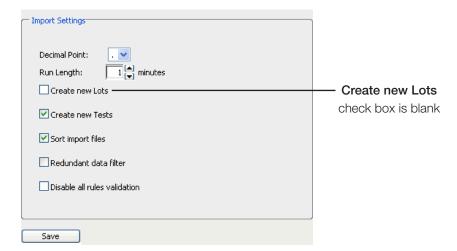

Import Settings—Not Using QC Data Management Software for SPC Rule Evaluation

- 1 Log in to the UnityWeb 2.0 or Unity Real Time online software.
- 2 Click the **Tools** tab.
- 3 Click Utilities.
- 4 Click Import Settings.

The following illustration shows recommended import settings for customers who do not use UnityWeb 2.0 or Unity Real Time online for SPC rule evaluation.

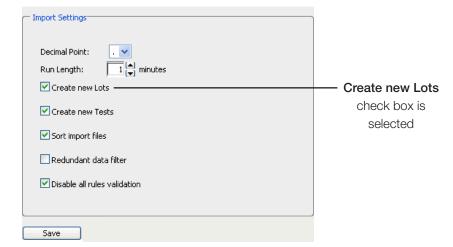

Configuring UnityConnect for a Data Stream with Desktop Software

In This Chapter

Before You Begin	94
Configuration Flow Chart	95
Output Path	96
Configuring a Watch Folder and Automatic Report	96
Configuring Import Options.	
Downloading a Transformer	
Configuring E-Mail Notification	

Before You Begin




This chapter is only for customers using a **data stream** from an instrument, middleware, or LIS with UnityConnect and any of the following **desktop** software:

Unity DesktopUnity Real Time

Configuration Flow Chart



Note: Installation, licensing, and downloading of one or more transformers is already completed by Bio-Rad.

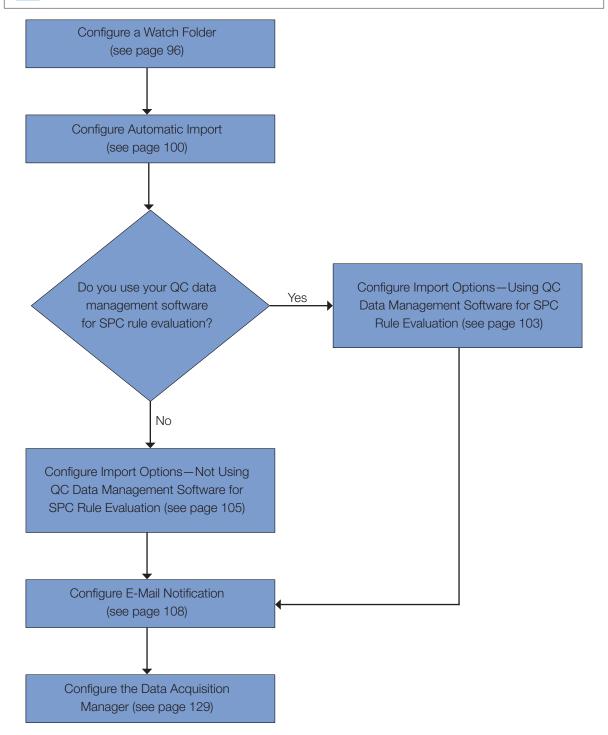

Output Path

Each installed and licensed transformer in UnityConnect requires an output path. The output path is the location where your transformed QC data is stored.

In most instances, the output path is the same location for all transformers. However, depending on your laboratory set up, multiple output paths may be useful. For example, if your laboratory uses multiple transformers, it is useful to have a separate output path for each transformer.



Important: The output path is pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, do not change the output path. Talk to your Bio-Rad representative for more information about output paths.

Configuring a Watch Folder and Automatic Report

UnityConnect can be configured to "watch" for data entering a "Watch Folder" and automatically transform and import the data into your QC data management software.

Configuring a Watch Folder

- 1 Create a new folder on the computer desktop where UnityConnect is installed.
- 2 Name the folder "Watch Folder."



Note: A Watch Folder must be created for each transformer. Give each Watch Folder a unique name if you create more than one. For example, "Watch Folder Meditech" and "Watch Folder RXL." This helps easily identify the Watch Folder for each transformer.


- 3 Start UnityConnect.
- 4 Click the **Options** menu and then click **Settings**.



5 Click the Watch Folder tab.



6 Click in the Watch Folder field for the transformer you want to configure.

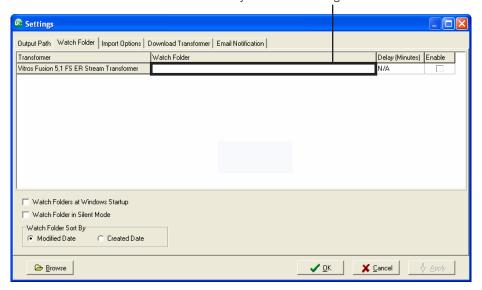


7 Click Browse


8 Select the Watch Folder you created on your computer desktop.




- 9 Click OK
- 10 Select the **Enable** check box for the transformer.

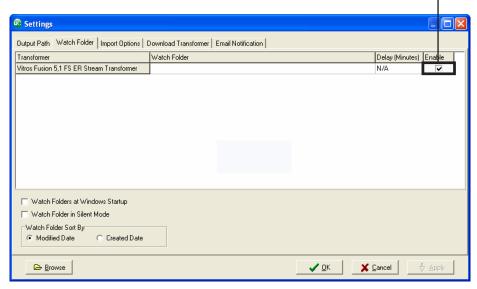


11 Select the **Watch Folders at Windows Startup** check box so the Watch Folder is active each time the computer starts.

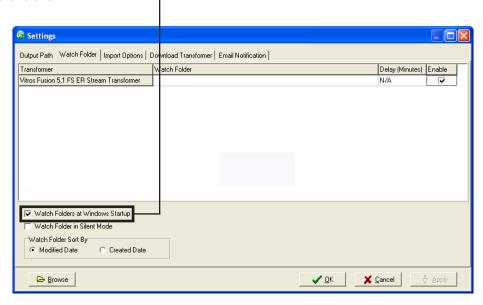

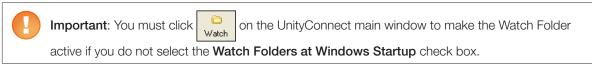

12 Select the Watch Folder in Silent Mode check box if you want the Watch Folder to run in silent mode.

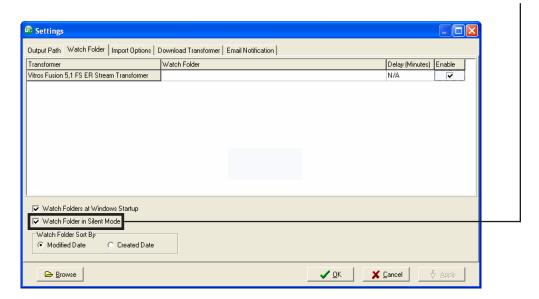



Note: Silent mode processes all data in the QC data stream that UnityConnect recognizes (that is, all data that has been configured in UnityConnect). If UnityConnect does not recognize any item(s) in the QC data stream, a copy of the original QC data is saved for later reprocessing when the unrecognized item(s) can be configured. UnityConnect can also be configured to send an e-mail notification when there is QC data requiring reprocessing. See "Configuring E-Mail Notification" on page 108 for more information.


- 13 Select a Watch Folder Sort By option:
 - Modified Date
 - Created Date
- 15 Continue with the next section, "Configuring Automatic Import."


Configuring Automatic Import



Important: Make sure you configure a Watch Folder before configuring automatic import. See "Configuring a Watch Folder" on page 96 for more information.

1 Click the **Import Options** tab.

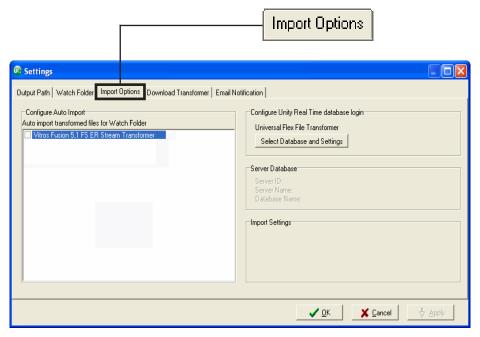

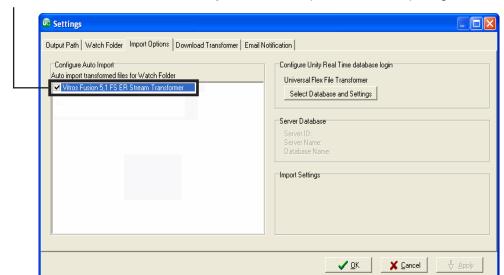

2 Select the check box for the transformer you want to set up for automatic importing.

- 3 Click the appropriate button:
 - Click v <u>DK</u> to apply the selected settings and close the **Settings** dialog box.
 - Click to apply the selected settings without closing the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.


Configuring Import Options

Transformed QC data is imported into your QC data management software based on the import options defined in UnityConnect. Select import options based on your use of your QC data management software for SPC rule evaluation.



Note: The import options only need to be configured one time. You do not need to configure the settings each time you transform data.

Use the import options described in the appropriate section according to how you use your your QC data management software for SPC rule evaluation.

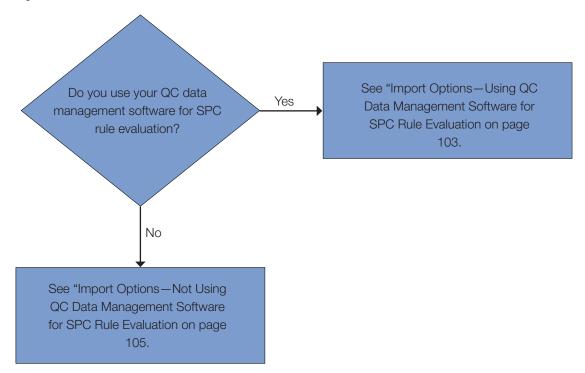


Import Options—Using QC Data Management Software for SPC Rule Evaluation


1 Click the **Options** menu and then click **Settings**.




2 Click the **Import Options** tab.

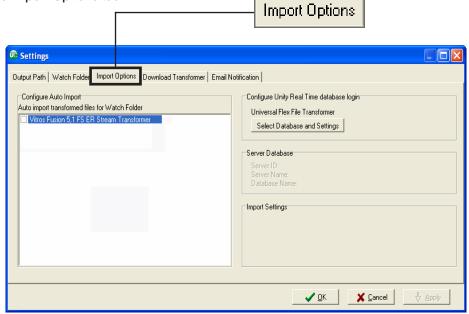

3 Select the transformer for which you want to configure import options.

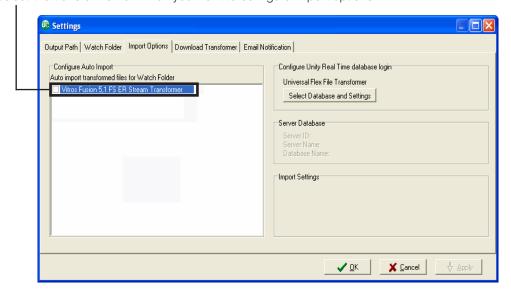

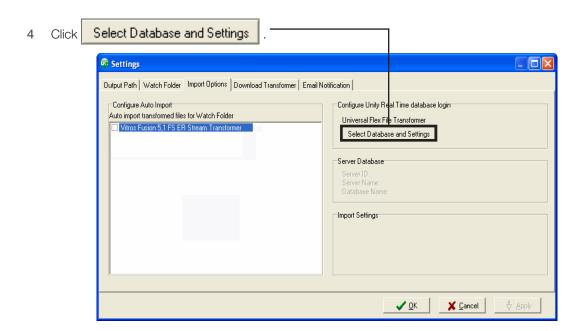

5 Make sure the **Create New Lots if Necessary** check box is not selected.

The following illustration shows recommended import settings for customers using their QC data management software for SPC rule evaluation.

The Create New Lots if Necessary check box is not selected.-👵 Unity Real Time Database Login Transformer: Import Settings Universal Flex File Transformer Decimal Point: SQL Server Database Run Length: 15 🕏 Minute(s) Connect using Windows Authentication Create New Lots if Necessary ▼ Create New Tests if Necessary SQL Server ID: Sort Import Files SQL Server Password: Redundant Data Filter Disable All Rules Violation SQL Server Name: ▼ Delete Rejection Log Records after 60 ♣ Day(s) Database Name: Operator ▾ User Name: User Password: X Cancel ✓ <u>o</u>K

This ensures UnityConnect will not create news lots when transforming and overwrite your current SPC rule settings in your QC data management software.

6 Click VOK.


To duplicate your current SPC rule settings from the current lot to a new lot, duplicate the lot in your QC data management software before transforming data for the lot for the first time. See "Configuring a New Lot" on page 164 for more information.

Import Options—Not Using QC Data Management Software for SPC Rule Evaluation


1 Click the **Options** menu and then click **Settings**.




2 Click the **Import Options** tab.

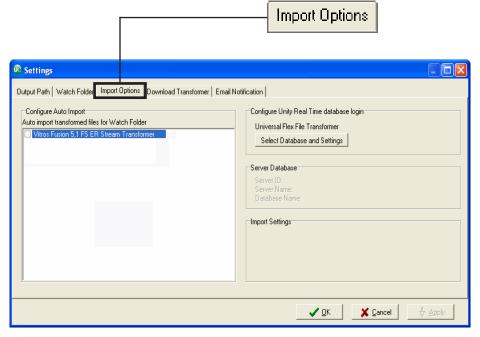

3 Select the transformer for which you want to configure import options. -

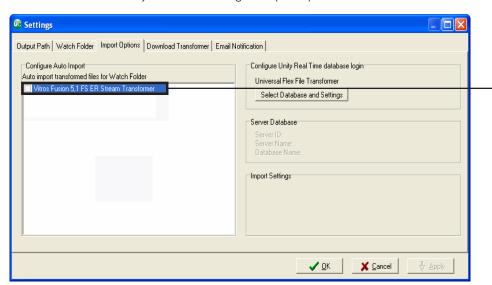


4 Click Select Database and Settings .

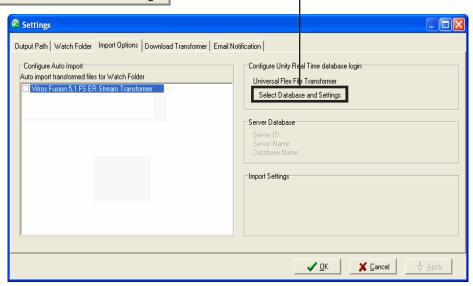

5 Make sure the **Create New Lots if Necessary** check box is selected.

The following illustration shows recommended import settings for customers not using their QC data management software for SPC rule evaluation.

Create new Lots if Necessary and Disable All Rules Violation check boxes are selected. -

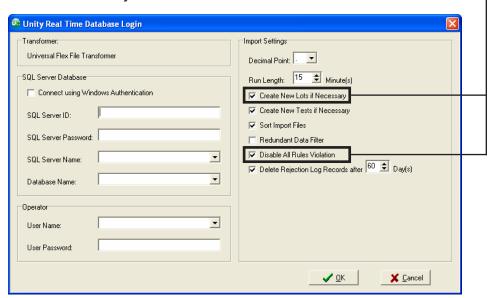

6 Click VOK.

A new lot is automatically created when UnityConnect finds a new lot in your QC data.

Downloading a Transformer

There are many instrument, middleware, and LIS transformers available for UnityConnect. Downloading a transformer is completed by Bio-Rad prior to configuring your UnityConnect software.





Important: Talk to your Bio-Rad representative for more information about licensing and downloading additional transformers.

Configuring E-Mail Notification

UnityConnect has the ability to send e-mail notifications to one or more e-mail addresses in the following circumstances:

- There is QC data requiring reprocessing.
- There is an error during automatic import.
- A Rejection Log has been created in UnityConnect.
- 1 Start UnityConnect.
- 2 Point to Options and then click Settings.

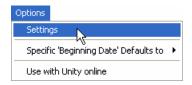

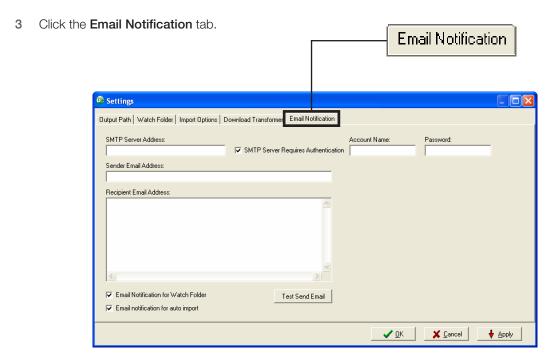

- 4 Type your SMTP server address in the **SMTP Server Address** field.
- 5 Select the **SMTP Server Requires Authentication** check box if your server requires authentication to send e-mail messages.
- Type the account name in the **Account Name** field if you selected the **SMTP Server Requires Authentication** check box.
- 7 Type the password in the **Password** field if you selected the **SMTP Server Requires Authentication** check box.

8 Type a sender e-mail address in the Sender Email Address field.



Tip: This helps prevent the sender e-mail address from being filtered as spam by your e-mail server.

- Type each e-mail address you want to receive notification in the **Recipient Email Address** field. Type each e-mail address on a separate line.
- 10 Click Test Send Email to test the e-mail settings.

A test e-mail is sent to each recipient.

- 11 Select the **Email Notification for Watch Folder** check box if you want recipients to receive an e-mail when there is QC data waiting for reprocessing.
- 12 Select the **Email notification for auto import** if you want recipients to receive an e-mail when an error or rejection occurs during import.
- 13 Click the appropriate button:
 - Click v <u>DK</u> to apply the selected settings and close the **Settings** dialog box.
 - Click to apply the selected settings without closing the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.

Configuring UnityConnect for a Data Stream with Web-based Software

In This Chapter

Before You Begin	110
Configuration Flow Chart	111
Setting Up UnityConnect for Web-based Software	112
Output Path	112
Configuring a Watch Folder	113
Configuring Internet Options	117
Configuring Auto Post	119
Downloading Lab Setup and Code Lists	122
Downloading a Transformer	123
Configuring E-Mail Notification	124
Configuring Import Settings	126

Before You Begin




This chapter is only for customers using a **data stream** from an instrument, middleware, or LIS with UnityConnect and the following **desktop** software:

Unity Real Time online

Configuration Flow Chart




Note: Installation, licensing, and downloading of one or more transformers is completed by Bio-Rad.

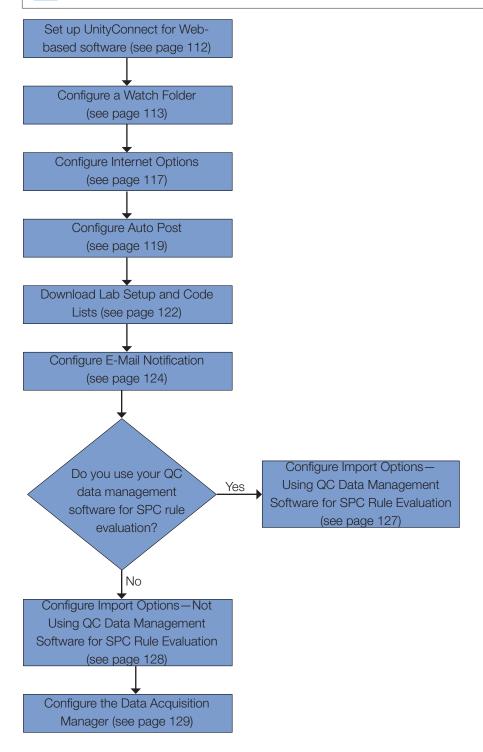

Setting Up UnityConnect for Web-based Software

- 1 Start UnityConnect.
- 2 Click the Options menu and then click Use with Unity online.

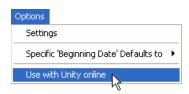

A check mark appears to the left indicating the **Use with Unity online** option is selected.



Output Path

Each installed and licensed transformer in UnityConnect requires an output path. The output path is the location where your transformed QC data is stored.

In most instances, the output path is the same location for all transformers. However, depending on your laboratory set up, multiple output paths may be useful. For example, if your laboratory uses multiple transformers, it is useful to have a separate output path for each transformer.



Important: The output path is pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, do not change the output path. Talk to your Bio-Rad representative for more information about output paths.

Configuring a Watch Folder


UnityConnect can be configured to "watch" for data entering a "Watch Folder."

- 1 Create a new folder on the computer desktop where UnityConnect is installed.
- 2 Name the folder "Watch Folder."




Note: A Watch Folder must be created for each transformer. Give each Watch Folder a unique name if you create more than one Watch Folder. For example, "Watch Folder Meditech" and "Watch Folder RXL." This helps easily identify the Watch Folder for each transformer.


- 3 Start UnityConnect.
- 4 Click the **Options** menu and then click **Settings**.

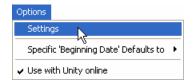

5 Click the Watch Folder tab.

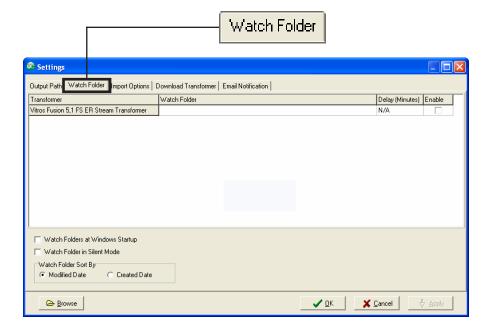

6 Click in the Watch Folder field for the transformer you want to configure.

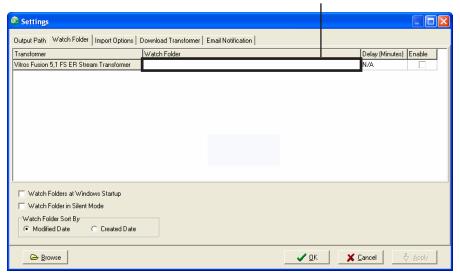


- 7 Click Browse .
- 8 Select the Watch Folder you created on your computer desktop.




9 Click OK .

10 Select the Enable check box for the transformer. -

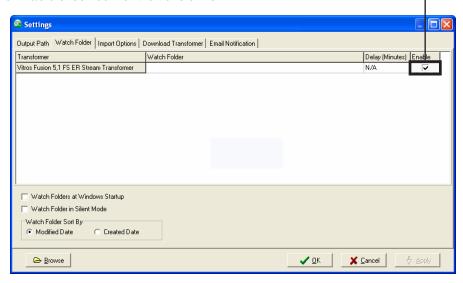


11 Select the Watch Folders at Windows Startup check box so the Watch Folder is active each time the computer starts. —

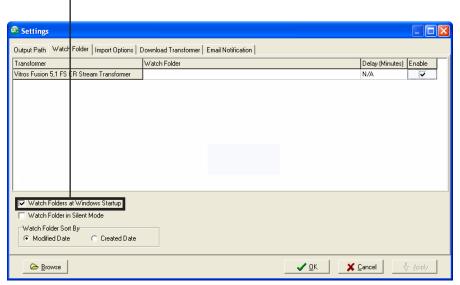



Note: You must click on the UnityConnect main window to make the Watch Folder active if you do not select the **Watch Folders at Windows Startup** check box.

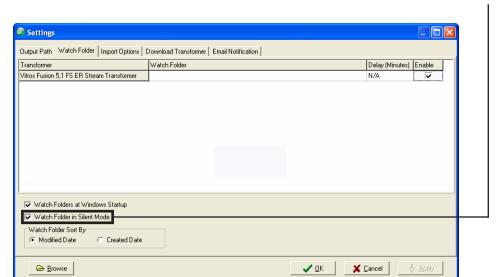

12 Select the Watch Folder in Silent Mode check box if you want the Watch Folder to run in silent mode.



Note: Silent mode processes all data in the QC data file that UnityConnect recognizes (that is, all data that has been configured in UnityConnect). If UnityConnect does not recognize any item(s) in the QC data file, a copy of the original QC data is saved for later reprocessing when the unrecognized item(s) can be configured. UnityConnect can also be configured to send an e-mail notification when there are QC data files requiring reprocessing. See "Configuring E-Mail Notification" on page 124 for more information.

- 13 Select a Watch Folder Sort By option:
 - Modified Date
 - Created Date
- **14** Click the appropriate button:
 - Click v <u>QK</u> to apply the selected settings and close the **Settings** dialog box.
 - Click to apply the selected settings without leaving the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.

Configuring Internet Options




Important: The Internet options are pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, talk to your Bio-Rad representative if you need further assistance with configuring Internet options.

- 1 Start UnityConnect.
- 2 Click the Options menu and then click Settings.



3 Click the **Internet Options** tab.

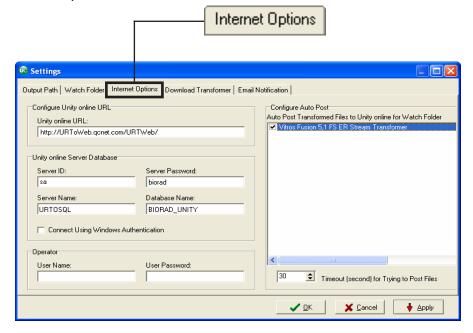

4 Verify the settings in the following fields:

Field	Setting
Unity online URL:	http://URToWeb.qcnet.com/URTWeb/
Server ID:	sa
Server Password:	biorad
Server Name:	URTOSQL
Database Name:	BIORAD_UNITY



Note: The settings shown above are the suggested settings at the time of this publication.

- 5 Type the user name you use when logging on to Unity Real Time online in the **User Name** field.
- 6 Type the password you use when logging on to Unity Real Time online in the User Password field.
- 7 Click the appropriate button:

 - Click to apply the selected settings without leaving the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.

Testing the Internet Connection

You can test the Internet connection after configuration to ensure the configuration is correct.


1 Click the File menu and then click Download Database.

The **Download Database** dialog box appears.



2 Click Test Connect

A message appears when when the connection is successful. -

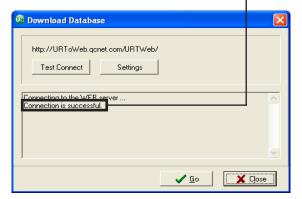

3 Click Close .

If the Internet Connection Fails

- 1 Click Settings .
- 2 Review the configuration and make any changes as required.
- 3 Repeat the testing as described on the previous page.
- 4 Contact Bio-Rad Software Support for assistance if the testing continues to fail.

Configuring Auto Post




Note: Auto post automatically uploads your transformed QC data to Unity Real Time online.


- 1 Start UnityConnect.
- 2 Click the Options menu and then click Settings.



3 Click the Internet Options tab.

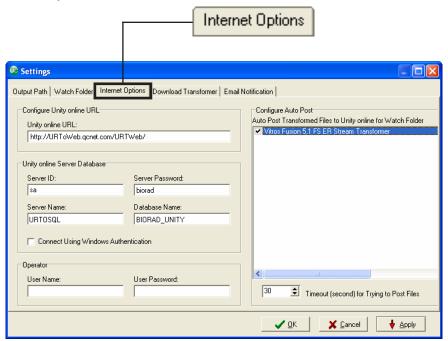

4 Select the check box for the transformer you want to configure for auto post.

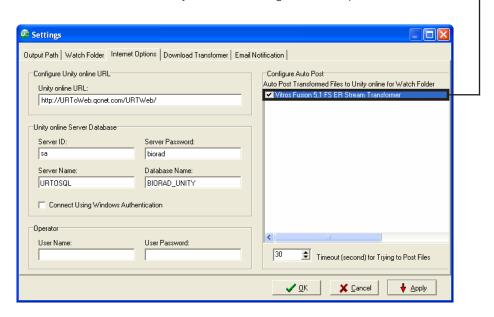



Note: The settings shown in the example above are the suggested settings at the time of this publication.

- 5 Type the user name you use when logging on to Unity Real Time online in the **User Name** field.
- 6 Type the password you use when logging on to Unity Real Time online in the User Password field.
- 7 Click the appropriate button:
 - Click to apply the selected settings and close the Settings dialog box.
 - Click to apply the selected settings without leaving the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.


Downloading the Database



Note: The database only needs to be download one time. See "Downloading Lab Setup and Code Lists" on page 122 for information about keeping the database up-to-date.

- 1 Click the **File** menu and then click **Download Database**.
- 2 Click 🗸 🗓 .

A message appears when when the download completes successfully.

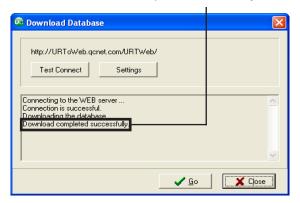


3 Click X Close

Downloading Lab Setup and Code Lists




Important: Perform the steps in this section any time you make a configuration change to your QC data management software. In addition, Bio-Rad recommends performing these steps at least monthly to ensure your software always has the most up-to-date codes lists and other information.

1 Click the File menu and then click Download Lab Setup/Code Lists.

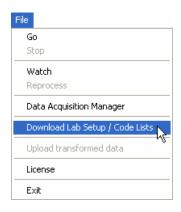


- 2 Select the check box for each item you want to download:
 - Download Lab Setup
 - Download Code Lists

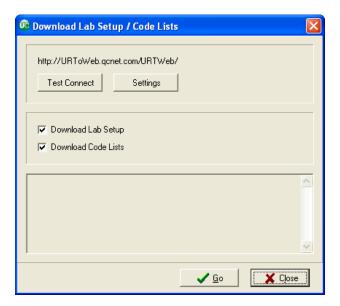

3 Click Go

A message appears when when the download completes successfully.

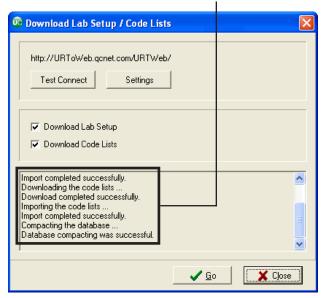

4 Click Close .

Downloading a Transformer

There are many instrument, middleware, and LIS transformers available for UnityConnect. Downloading a transformer is completed by Bio-Rad prior to configuring your UnityConnect software.





Important: Talk to your Bio-Rad representative for more information about licensing and downloading additional transformers.

Configuring E-Mail Notification

UnityConnect has the ability to send e-mail notifications to one or more e-mail addresses when:

- There are QC data files requiring reprocessing.
- There is an error during automatic posting.
- 1 Start UnityConnect.
- 2 Point to the **Options** menu and then click **Settings**.

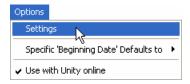

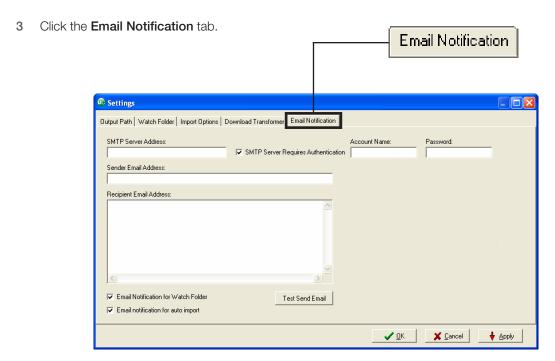

- 4 Type your SMTP server address in the **SMTP Server Address** field.
- 5 Select the **SMTP Server Requires Authentication** check box if your server requires authentication to send e-mail messages.
- Type the account name in the **Account Name** field if you selected the **SMTP Server Requires Authentication** check box.
- 7 Type the password in the **Password** field if you selected the **SMTP Server Requires Authentication** check box.

8 Type a sender e-mail address in the **Sender Email Address** field.



Tip: This helps prevent the sender e-mail address from being filtered as spam by your e-mail server.

- 9 Type each e-mail address you want to receive notification in the **Recipient Email Address** field. Type each e-mail address on a separate line.
- 10 Click Test Send Email to test the e-mail settings.

A test e-mail is sent to each recipient.

- 11 Select the **Email Notification for Watch Folder** check box if you want recipients to receive an e-mail when a file is waiting for reprocessing.
- 12 Select the Email Notification for Auto Post if you want recipients to receive an e-mail when an error has occurred during auto-post.
- 13 Click the appropriate button:
 - Click to apply the selected settings and close the **Settings** dialog box.
 - Click to apply the selected settings without leaving the **Settings** dialog box. This is convenient if you want to continue configuring settings on other tabs.


Configuring Import Settings

Transformed QC data is posted to your QC data management software based on the import settings defined in the software. Select import settings based on your use of your QC data management software for SPC rule evaluation.



Note: The import settings only need to be configured one time. You do not need to configure the settings each time you transform data.

Use the import settings described in the appropriate section according to how you use your QC data management software for SPC rule evaluation.

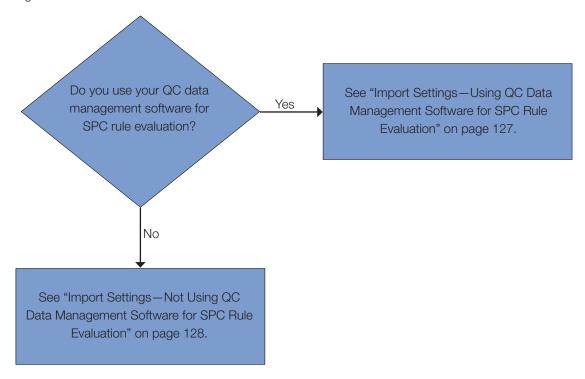


Import Settings—Using QC Data Management Software for SPC Rule Evaluation

These settings ensure UnityConnect will not create news lots when transforming and overwrite your current SPC rule settings in your QC data management software.

To duplicate your current SPC rule settings from the current lot to the new lot, leave the **Create new Lots** check box blank and duplicate the lot in your QC data management software before transforming data for the lot for the first time. See "Configuring a New Lot" on page 164 for more information.

- 1 Log in to Unity Real Time online.
- 2 Click the Tools tab.
- 3 Click Utilities.
- 4 Click Import Settings.


The following illustration shows recommended import settings for customers using Unity Real Time online for SPC rule evaluation.

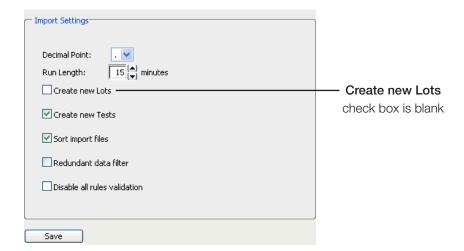

Import Settings—Not Using QC Data Management Software for SPC Rule Evaluation

- 1 Log in to Unity Real Time online.
- 2 Click the **Tools** tab.
- 3 Click Utilities.
- 4 Click Import Settings.

The following illustration shows recommended import settings for customers who do not use Unity Real Time online for SPC rule evaluation.

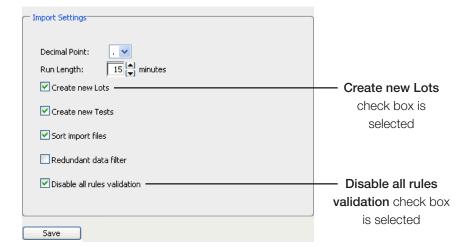

Data Acquisition Manager

In This Chapter

Before You Begin	129
Overview	130
Installing the Service	130
Working with Devices	
Starting and Stopping the Service	132
Exception Files	134
Device Logs	137

Before You Begin



This chapter is only for customers using a **data stream** from an instrument, middleware, or LIS with UnityConnect and any of the following QC data management software:

Unity Desktop
Unity Real Time
Unity Real Time online

Overview

The Data Acquisition Manager in UnityConnect provides the ability to collect data from an instrument, middleware, or LIS on a real-time basis when installed with the proper network or wireless hardware. The existing flow of data (QC, patient, and support) between the instrument, middleware, or LIS is not affected.

Installing the Service

The Data Acquisition Manager Service (Service) captures data from your instrument, middleware, or LIS (devices). Installing the Service activates the connections between the Data Acquisition Manager and devices.



Note: This is a one-time setup. The message "Service is running" appears if the Data Acquisition Manager Service is already installed.

- 1 Start UnityConnect.
- 2 Click the File menu and then click Data Acquisition Manager.
 - The Data Acquisition Manager dialog box appears.
- 3 Click the Service menu and then click Install Service.
 - The message "Service is installed but is not started" appears.
- 4 See "Starting and Stopping the Service" on page 132 for more information.

Working with Devices

The Data Acquisition Manager in UnityConnect supports instrument, middleware, and LIS devices.



Important: The device settings are pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, do not change the device settings. Talk to your Bio-Rad representative for more information if you need to add another device or change the device settings.

Adding a Device

- Start UnityConnect.
- 2 Click the File menu and then click Data Acquisition Manager.
 - The Data Acquisition Manager dialog box appears.
- 3 Click the **Device** menu and then click **Add**.



Important: The device settings are pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, do not change the device settings. Talk to your Bio-Rad representative for more information if you need to add another device or change the device settings.

Deleting a Device



Important: Deleting a device permanently removes the configuration from the Data Acquisition Manager. Talk to your Bio-Rad representative for more information before you delete a device.

- Start UnityConnect.
- 2 Click the File menu and then click Data Acquisition Manager.

The Data Acquisition Manager dialog box appears.

- 3 Select the device you want to delete.
- 4 Click the **Device** menu and then click **Delete**.

Modifying a Device



Important: The device settings are pre-configured by Bio-Rad. Due to the unique configuration of your customized UnityConnect software, do not modify the device settings. Talk to your Bio-Rad representative for more information if you need to modify the device settings.

Connecting a Single Device

- 1 Start UnityConnect.
- 2 Click the File menu and then click **Data Acquisition Manager**.

The **Data Acquisition Manager** dialog box appears.

- 3 Select the device you want to connect.
- 4 Click the **Device** menu and then click **Connect**.

Connecting All Devices

- 1 Start UnityConnect.
- 2 Click the File menu and then click Data Acquisition Manager.

The Data Acquisition Manager dialog box appears.

3 Click the **Device** menu and then click **Connect All Devices**.

Disconnecting a Single Device

- 1 Start UnityConnect.
- 2 Click the File menu and then click Data Acquisition Manager.
 - The Data Acquisition Manager dialog box appears.
- 3 Select the device you want to disconnect.
- 4 Click the **Device** menu and then click **Disconnect**.

Disconnecting All Devices

- 1 Start UnityConnect.
- 2 Click the File menu and then click Data Acquisition Manager.
 - The Data Acquisition Manager dialog box appears.
- 3 Click the **Device** menu and then click **Disconnect All Devices**.

Starting and Stopping the Service

Starting the Service initiates communication between the Data Acquisition Manager and your configured instrument, middleware, or LIS.

- 1 Start UnityConnect.
- 2 Click the File menu and then click Data Acquisition Manager.
 - The **Data Acquisition Manager** dialog box appears.
- 3 Click the Service menu and then click Start Service.

Symbol	Indicates
Green circle	Successful connection of the Data Acquisition Manager Service and data flow from the instrument, middleware, or LIS.
Red stop symbol	Successful connection of the Data Acquisition Manager Service but no data flow from the instrument, middleware, or LIS.
No symbol	No connection of the Data Acquisition Manager Service. The Service is not started or the device is disconnected.

Stopping the Service

- 1 Start UnityConnect.
- 2 Click the File menu and then click Data Acquisition Manager.
 - The **Data Acquisition Manager** dialog box appears.
- 3 Click the **Service** menu and then click **Stop Service**.


Symbol	Indicates
Green circle	Successful connection of the Data Acquisition Manager Service and data flow from the instrument, middleware, or LIS.
Red stop symbol	Successful connection of the Data Acquisition Manager Service but no data flow from the instrument, middleware, or LIS.
No symbol	No connection of the Data Acquisition Manager Service. The Service is not started or the device is disconnected.

Exception Files


An exception file is created each time there is communication between the Data Acquisition Manager and a device. Viewing exception files can be helpful if you are troubleshooting communication problems.

Viewing Exception Files


- 1 Click the File menu and then click Data Acquisition Manager.
- 2 Select the device you want to view.
- 3 Click the Exception tab.

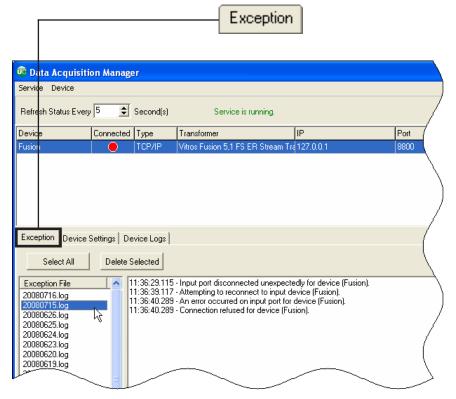

4 Select the date of the exception files you want to view. The most recent day is at the top of the list.

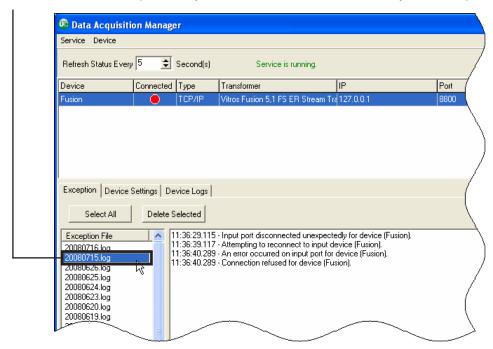


The details appear here.

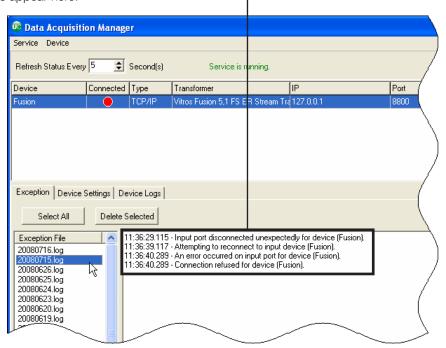

Deleting Exception Files

- 1 Click the File menu and then click Data Acquisition Manager.
- 2 Select the device with the exception files you want to delete.
- 3 Click the Exception tab.

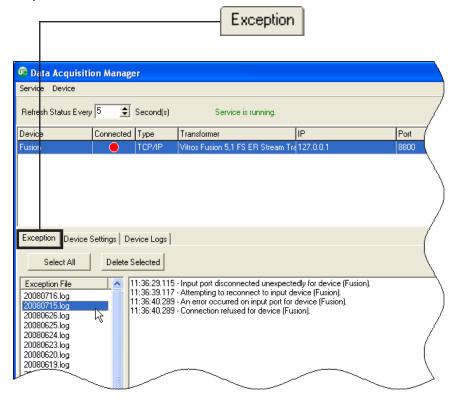


- 4 Select an individual exception file to delete or click Select All
- 5 Click Delete Selected .


Device Logs

Viewing Device Logs

- 1 Click the File menu and then click Data Acquisition Manager.
- 2 Click the **Device Logs** tab.

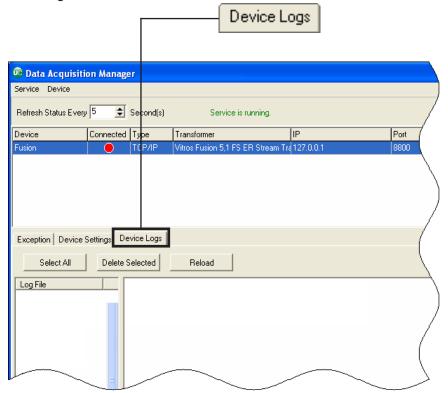


3 Select the device log you want to view.

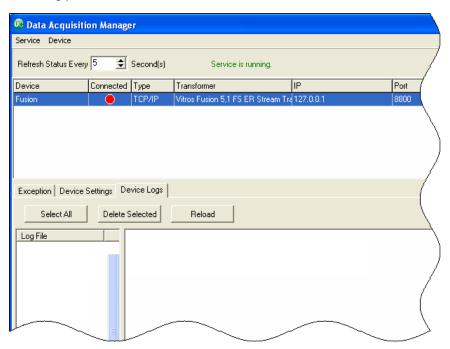

The details appear here. Data Acquisition Manager Service Device Refresh Status Every 5 Second(s) Service is running. Connected Type Device Transformer ĺΡ Port Vitros Fusion 5,1 FS ER Exception Device Settings Device Logs Select All Delete Selected Reload Log File

Deleting Device Logs

- 1 Click the File menu and then click Data Acquisition Manager.
- 2 Click the Device Logs tab.

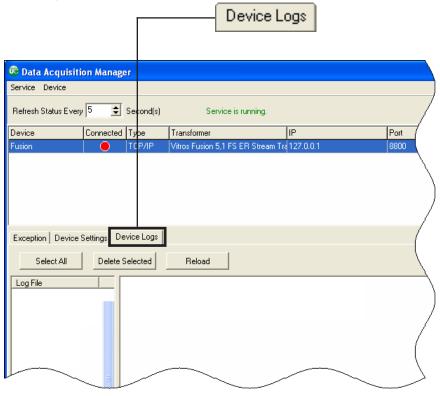

- 3 Select an individual device log to delete or click Select All
- 4 Click Delete Selected

Transforming Data for Desktop Software

In This Chapter

Overview	140
Transforming Data — Data File Transfer Method	141
Reprocessing QC Data Files	145
Working with Multiple QC Data Files	148

Overview



Note: This chapter describes the transforming process for the data file transfer method only.

Transforming data is the process of translating your QC data from the format of your instrument, middleware, or LIS into the Bio-Rad format. After transforming, data is automatically imported into your QC data management software.

Before You Begin



This chapter is only for customers using UnityConnect with the data file transfer method and any of the following desktop software:

Unity Desktop
Unity Real Time



Best Practices-Transform Data Often

Bio-Rad recommends transforming data as often as possible throughout the month. This helps you become more familiar with the transformation process and alleviates the stress and pressure at the end of the month.

Transforming Data – Data File Transfer Method



Before You Begin

- Place the QC data file in the QC Data Folder on your computer desktop or insert the CD or disk into the disk drive.
- 2 Give the QC data file a unique name so it can easily be identified.
- 3 Verify the import options as described in "Configuring Import Options" on page 68.




Note: The import options only need to be configured one time. You do not need to configure the settings each time you transform data.

- 1 Open UnityConnect.
- 2 Set the date range for the data to transform. Data can be transformed for a particular date range, from the first point in the QC data file to the last point in the QC data file, or a combination of both.







Note: See "Configuring the Default Date Settings" on page 75 for information about customizing the date range.

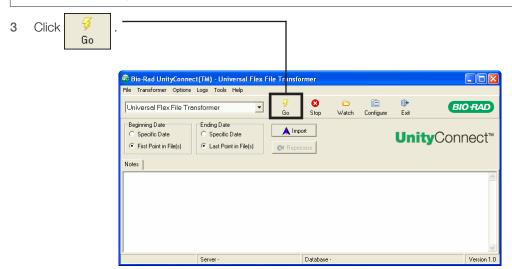


4 Select the QC data file to transform.

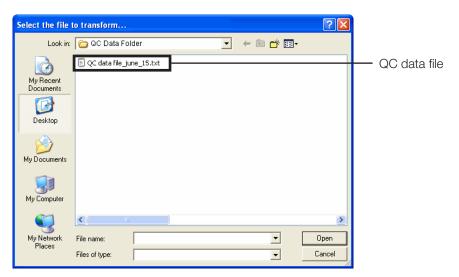


The QC data file name appears here after selecting it.

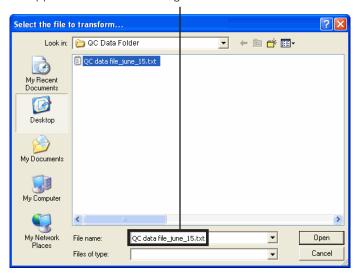

5 Click Open

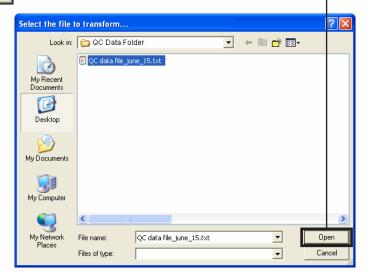


A message states the transformation is in progress.

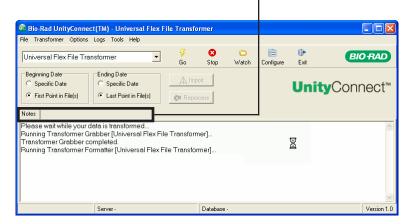



Note: A message appears if UnityConnect finds a new QC item requiring configuration. See Chapter 10, "New Configuration in UnityConnect" for more information.

A second message appears when transformation is complete.

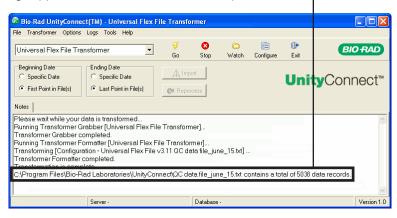

6 Read the information in the Notes section carefully and note the number of data records transformed.



Important: Transformed data is automatically imported into your QC data management software if configured to do so. See "Configuring Watch Folder and Automatic Import" on page 64 for more information.


Reprocessing QC Data Files

When UnityConnect is used in silent watch mode, files requiring configuration are stored for reprocessing at a later date. Data already configured is transformed and imported. See "Configuring Watch Folder and Automatic Import" on page 64 for more information.




Note: An e-mail is sent to the recipients listed on the **Email Notification** tab stating a file is waiting to be reprocessed if configured to do so. See "Configuring E-Mail Notification" on page 89 for more information. If an e-mail notification is not configured, the Reprocess button in UnityConnect indicates a file is ready to be reprocessed.

1 Start UnityConnect.

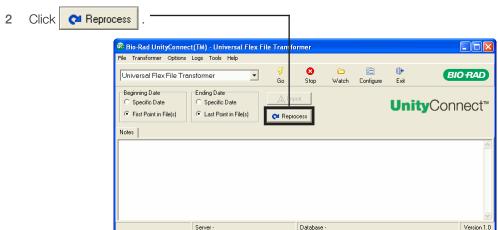

A message appears when UnityConnect finds the item(s) requiring configuration.

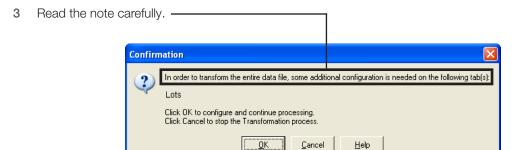



Important: This is not an error message. The message is simply stating there are new codes in the QC data file that have not been matched to Bio-Rad codes in your QC data management software.


- 4 Click <u>OK</u>
- 5 Log in to your QC data management software database if prompted.

UnityConnect presents the tab requiring configuration.

6 Configure all tabs UnityConnect presents.




Tip: See Chapter 10, "New Configuration in UnityConnect" for more information.

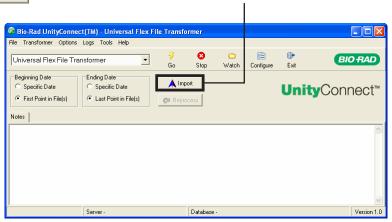




Note: Automatic import does not occur when reprocessing QC data files even if UnityConnect is configured for automatic import.

The **Import transformed data** dialog box appears.

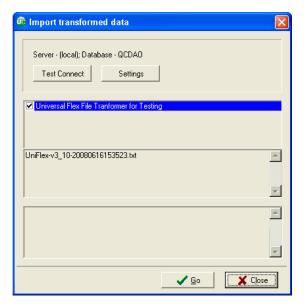

8 Click **4** <u>G</u>o .

A message appears when importing is complete.

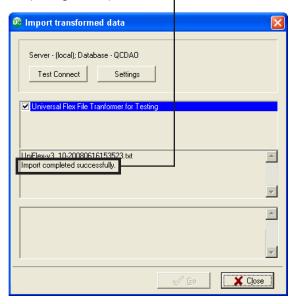

9 Click X Close

Best Practicies—Review the Number of Data Records Transformed

Always review and note the number of data records transformed. Customers who transform data on a regular basis at the same time each month will become familiar with the typical number of data points in their transformed QC data file. For example, if your transformed file normally contains 10,000 data points (records) for the month but the message indicates 3,000 records, further investigation may be required to ensure the QC data file you transformed contained all of your QC data for the month.

Working with Multiple QC Data Files



Important: Due to the unique configuration of your customized UnityConnect software, talk to your Bio-Rad representative about the best process for transforming multiple QC data files for your laboratory.


Some Laboratory Information Systems (LIS) output multiple files when creating QC data files. One example is the Orchard LIS. UnityConnect can be configured to automatically process one directory containing multiple files which simplifies the transformation process.

Configuring UnityConnect to Transform a Directory of Files



Note: You must perform a one-time configuration for each transformer.

1 Select the transformer you want to set up to process a directory of files.

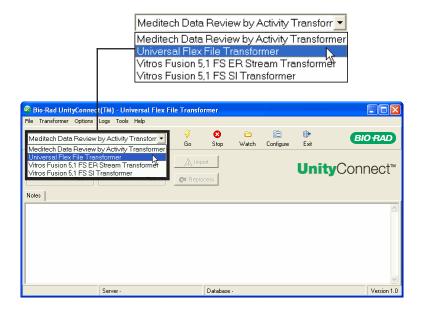

2 Click the **Transformer** menu and then click **Edit this Transformer**.




3 Select the Process One Directory at a Time option. -




4 Click VOK.


Transforming a Directory of Files

- 1 Open UnityConnect.
- 2 Set the date range for the data to transform. Data can be transformed for a particular date range, from the first point in the QC data file to the last point in the QC data file, or a combination of both.

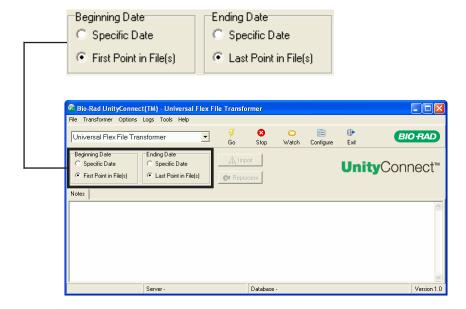

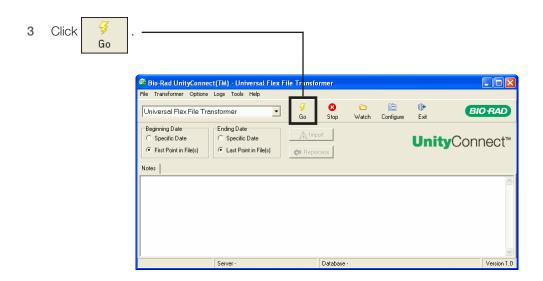

4 Select the directory (or folder) containing the multiple data files to transform.

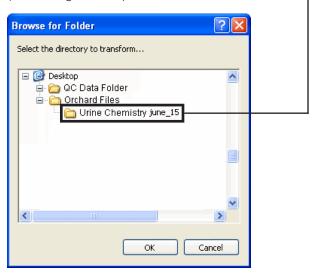

5 Click OK

Best Practices—Review the Number of Data Records Transformed

Always review and note the number of data records transformed. Customers who transform data on a regular basis at the same time each month will become familiar with the typical number of data points in their transformed QC data file. For example, if your transformed file normally contains 10,000 data points (records) for the month but the message indicates 3,000 records, further investigation may be required to ensure the QC data file you transformed contained all of your QC data for the month.



Transforming Data for Web-based Software

In This Chapter

Overview	151
Transforming Data—Data File Transfer Method	152
Reprocessing QC Data Files	156
Working with Multiple QC Data Files	159

Overview



Note: This chapter describes the transforming process for the data file transfer method only.

Transforming data is the process of translating your QC data from the format of your instrument, middleware, or LIS into the Bio-Rad format. After transforming, data is automatically imported into your QC data management software.



Before You Begin

This chapter is only for customers using UnityConnect with any of the following **Web-based** software:

UnityWeb 2.0

Unity Real Time online



Best Practices-Transform Data Often

Bio-Rad recommends transforming data as often as possible throughout the month. This helps you become more familiar with the transformation process and alleviates the stress and pressure at the end of the month.

Transforming Data – Data File Transfer Method




Before You Begin

- Place the QC data file in the QC Data Folder on your computer desktop or insert the CD or disk into the disk drive.
- 2 Give the QC data file a unique name so it can easily be identified.
- 3 Verify the import options as described in "Configuring Import Options" on page 68.




Note: The import options only need to be configured one time. You do not need to configure the settings each time you transform data.

- 1 Open UnityConnect.
- 2 Set the date range for the data to transform. Data can be transformed for a particular date range, from the first point in the QC data file to the last point in the QC data file, or a combination of both.

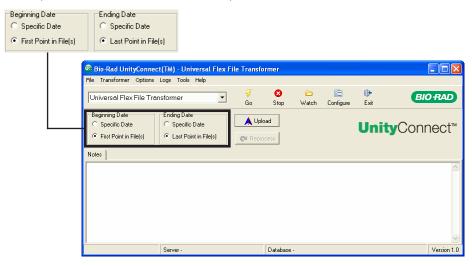





Note: See "Configuring the Default Date Settings" on page 75 for information about customizing the date range.

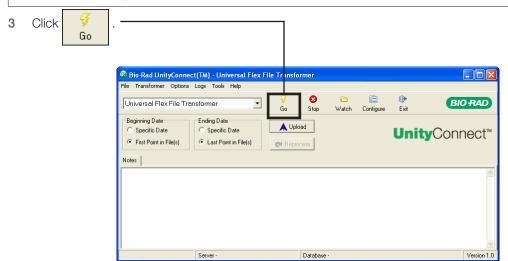


4 Select the QC data file to transform.

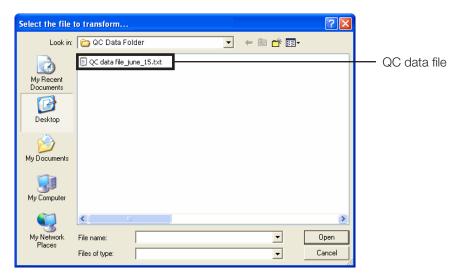


The QC data file name appears here after selecting it.

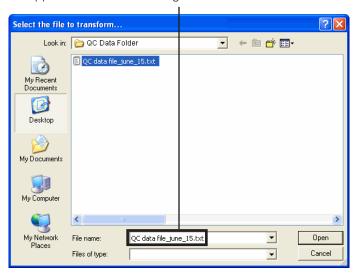

5 Click Open

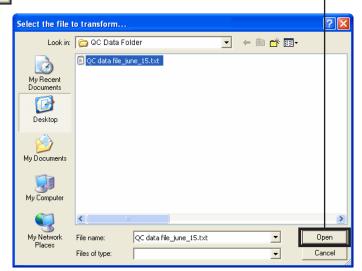


A message states the transformation is in progress. -

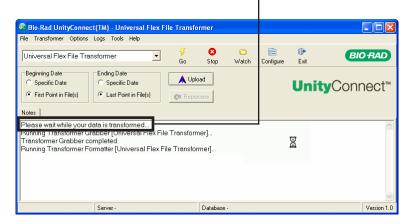



Note: A message appears if UnityConnect finds a new QC item requiring configuration. See Chapter 10, "New Configuration in UnityConnect" for more information.

A second message appears when transformation is complete.

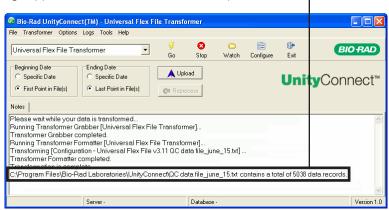

6 Read the information in the **Notes** section carefully and note the number of data records transformed.



Important: Transformed data is automatically imported into your QC data management software if configured to do so. See "Configuring E-Mail Notification" on page 74 for more information.

Reprocessing QC Data Files

When UnityConnect is used in silent watch mode, files requiring configuration are stored for reprocessing at a later date. Data already configured are transformed and imported.




Note: An e-mail is sent to the recipients listed on the **Email Notification** tab stating a file is waiting to be reprocessed if configured to do so. See "Configuring E-Mail Notification" on page 74 for more information. If an e-mail notification is not configured, the Reprocess button in UnityConnect indicates a file is ready to be reprocessed.

1 Start UnityConnect.



A message appears when UnityConnect finds the item(s) requiring configuration.

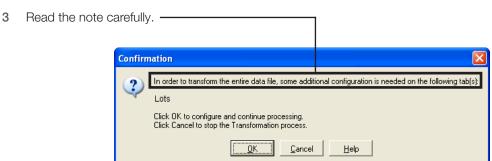

Important: This is not an error message. The message is simply stating there are new codes in the QC data file that have not been matched to Bio-Rad codes in your QC data management software.

- 4 Click OK
- 5 Log in to your QC data management software database if prompted.


UnityConnect presents the tab requiring configuration.

6 Configure all tabs UnityConnect presents.




Tip: See Chapter 10, "New Configuration in UnityConnect" for more information.

7 Click A Upload after the file has completed transforming.

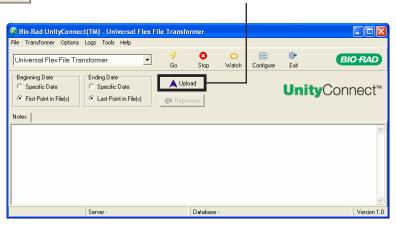


Note: Automatic upload does not occur when reprocessing QC data files even if UnityConnect is configured for automatic import.

The **Upload transformed data** dialog box appears.

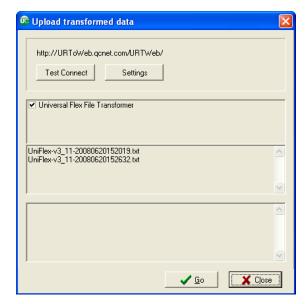

8 Click Go.

A message appears when the upload completes successfully.

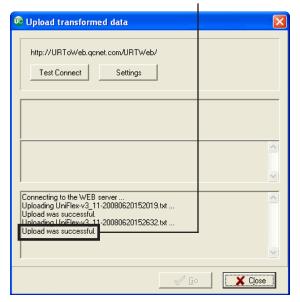

9 Click Close

Best Practices—Review the Number of Data Records Transformed

Always review and note the number of data records transformed. Customers who transform data on a regular basis at the same time each month will become familiar with the typical number of data points in their transformed QC data file. For example, if your transformed file normally contains 10,000 data points (records) for the month but the message indicates 3,000 records, further investigation may be required to ensure the QC data file you transformed contained all of your QC data for the month.

Working with Multiple QC Data Files



Important: Due to the unique configuration of your customized UnityConnect software, talk to your Bio-Rad representative about the best process for transforming multiple QC data files for your laboratory.


Some Laboratory Information Systems (LIS) output multiple files when creating QC data files. One example is the Orchard LIS. UnityConnect can be configured to automatically process one directory containing multiple files which simplifies the transformation process.

Configuring UnityConnect to Transform a Directory of Files



Note: You must perform a one-time configuration for each transformer.

1 Select the transformer you want to set up to process a directory of files.

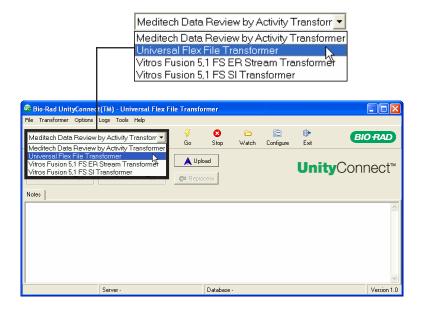

2 Click the Transformer menu and then click Edit this Transformer.




3 Select the Process One Directory at a Time option.




4 Click VOK.


Transforming a Directory of Files

- 1 Open UnityConnect.
- 2 Set the date range for the data to transform. Data can be transformed for a particular date range, from the first point in the QC data file to the last point in the QC data file, or a combination of both.

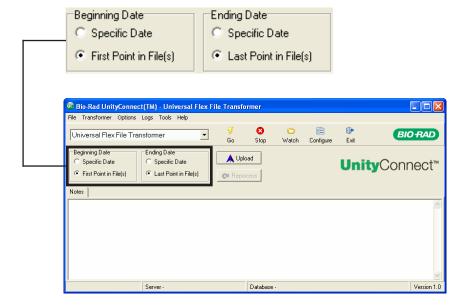

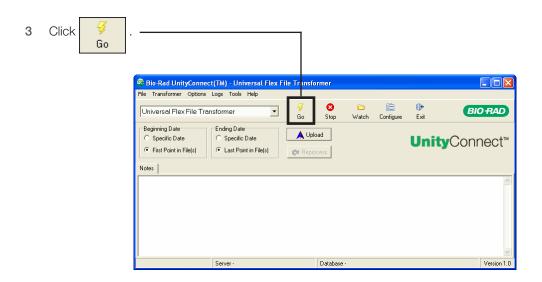

4 Select the directory (or folder) containing the multiple data files to transform.

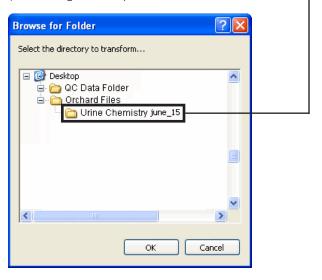

5 Click OK

Best Practices–Review the Number of Data Records Transformed

Always review and note the number of data records transformed. Customers who transform data on a regular basis at the same time each month will become familiar with the typical number of data points in their transformed QC data file. For example, if your transformed file normally contains 10,000 data points (records) for the month but the message indicates 3,000 records, further investigation may be required to ensure the QC data file you transformed contained all of your QC data for the month.

New Configuration in UnityConnect

In This Chapter

Overview	162
Configuring a New Lot	164
Other Types of Configuration	
Configuring New Data File Options	190
Configuring a New Instrument	190
Configuring a New Test	195

Overview

Bio-Rad customized the UnityConnect software based on all the information in the original QC data file provided by your laboratory. After the initial configuration is complete, new configuration is required only if a new QC item (such as an instrument, lot, or test) is added to the QC data file.



Important: Configuring information in UnityConnect varies according to the information in your QC data file and the specifics of your customized UnityConnect software. Therefore, the example configuration information shown in this chapter may not be identical.

UnityConnect processes the QC data file and, if necessary, displays messages prompting you to match or configure the new information in the QC data file with the information in your Unity software. UnityConnect automatically determines the information needing configuration and presents the appropriate tab where the configuration can be completed.

Summary of New Configuration




UnityConnect prompts
you to match the new
information in the QC
data file to the Bio-Rad
information.

Configuring the new QC item is a one-time process.

After the new QC item is configured, UnityConnect recognizes the item for all future transformations.

Configuring a New Lot

The process for configuring a new lot depends on your use of your QC data management software for SPC rule evaluation.

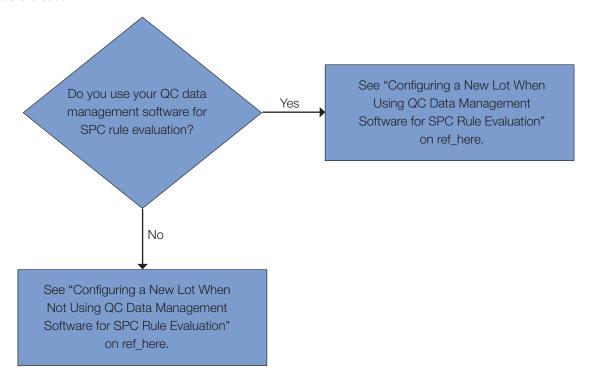

Configuring a New Lot When Using QC Data Management Software for SPC Rule Evaluation

Switching to a new lot of control material is one of the most common reasons to complete new configuration in UnityConnect. Configuring UnityConnect for a new lot is a simple two-step process.

- 1 Duplicate the current lot in your QC data management software.
- 2 Complete new configuration in UnityConnect.

Duplicating a Lot in Your QC Data Management Software



Tip: Duplicating a lot is helpful when you switch to a new control material lot number. All tests under the old lot are duplicated.

See the appropriate section according to your software for information on duplicating a lot:

- Unity Desktop (page 165)
- Unity Real Time (page 166)
- UnityWeb 2.0 (page 166)
- Unity Real Time online (page 167)

Duplicating a Lot in Unity Desktop

- 1 Use one of the following methods to open the **Lot** dialog box:
 - On the toolbar, click Lot
 - Click the **Select** menu and then click **Lot**.
 - Press CTRL+F3.

The **Lot** dialog box appears.

- 2 Select the lot you want to duplicate in the **Open lots** list.
- 3 Click Duplicate.

The **Duplicate Lot** dialog box appears.

4 Choose the **new lot number**.



Note: The **Fixed means**, **Fixed SDs**, and **Target values for Analytical Goals** options apply only to customers using Westgard (SPC) rules in Unity Desktop.

5 Click OK.

The new lot appears in the **Open lots** list.



Tip: The **Fixed means**, **Fixed SDs**, and **Target values for Analytical Goals** options apply only to customers using Westgard (SPC) rules in Unity Desktop.

6 Click Close.

7 Continue with the section, "Completing New Configuration in UnityConnect" on page 167.

Duplicating a Lot in Unity Real Time

- 1 Use one of the following methods to open the **Lot** dialog box:
 - On the toolbar, click
 - Click the **Select** menu and then click **Lot**.
 - Press CTRL+F3.

The Lot dialog box appears.

- 2 Select the lot you want to duplicate in the **Open lots** list.
- 3 Click Duplicate.

The **Duplicate Lot** dialog box appears.

- 4 Choose the new lot number.
- 5 Click OK.

The new lot appears in the **Open lots** list.



Tip: The **Fixed means**, **Fixed SDs**, and **Target values for Analytical Goals** options apply only to customers using Westgard (SPC) rules in Unity Real Time.

- 6 Click Close.
- 7 Continue with the section, "Completing New Configuration in UnityConnect" on page 167.

Duplicating a Lot in UnityWeb 2.0

- 1 Click the **Configure** tab.
- 2 Click Lot.
- 3 If you have more than one lab number, make sure the correct number appears in the Lab number list.
- 4 Select the lot number you want to duplicate in the **Open lots** list.
- 5 Click Duplicate.
- 6 Select the new lot number from the **New lot number** list.



Tip: Only unexpired lot numbers in the same product group appear in the list.

- 7 Select the appropriate check box for each parameter to duplicate:
 - Fixed means
 - Fixed SDs
- 8 Click OK.

The new lot number appears at the bottom of the Open lots list.

9 Continue with the section, "Completing New Configuration in UnityConnect" on ref. here.

Dulicating a Lot in Unity Real Time online

- Click the Configure tab.
- 2 Click Lot.
- 3 Make sure the correct number appears in the Lab number list if you have more than one lab number.
- 4 Select the lot number to duplicate in the **Open lots** list.
- 5 Click Duplicate.
- 6 Select the new lot number from the **New lot number** list.



Note: Only unexpired lot numbers in the same product group appear in the Lot number list.


7 Click OK.


The new lot number appears at the bottom of the **Open lots** list.

8 Continue with the section, "Completing New Configuration in UnityConnect" on ref_here.


Completing New Configuration in UnityConnect

- 1 Start UnityConnect.
- 2 Set the date range for the data to transform. Data can be transformed for a particular date range, from the first point in the QC data file to the last point in the QC data file, or a combination of both.

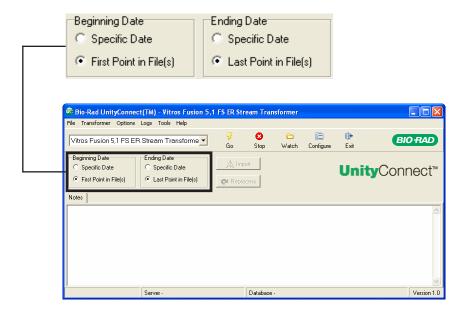

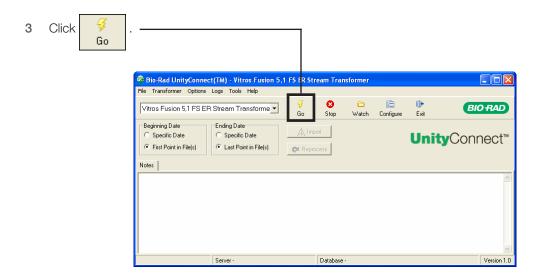


4 Select the QC data file to transform.

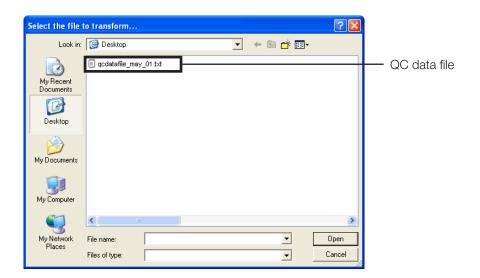




Tip: Bio-Rad recommends giving the QC data file a unique name so it can easily be identified.


The QC data file name appears here after selecting it.

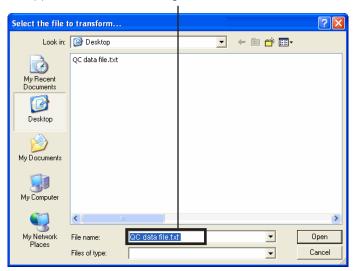


5 Click Open

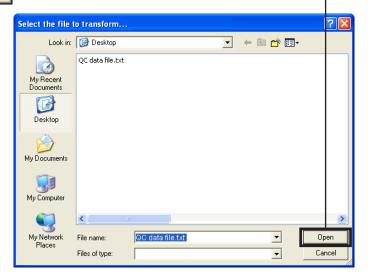


A message states the transformation is in progress. —

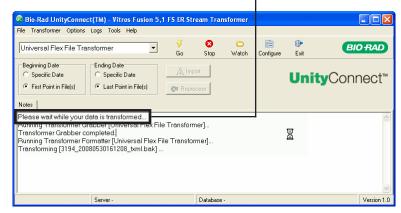


6 A message appears when UnityConnect finds the new lot. Read the note carefully.

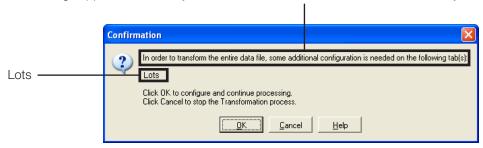

- Important: This is not an error message. The message is simply stating there are new codes in the QC data file that have not been matched to Bio-Rad codes in your QC data management software.
- 7 Click OK .
- 8 Log in to your QC data management software database if prompted.

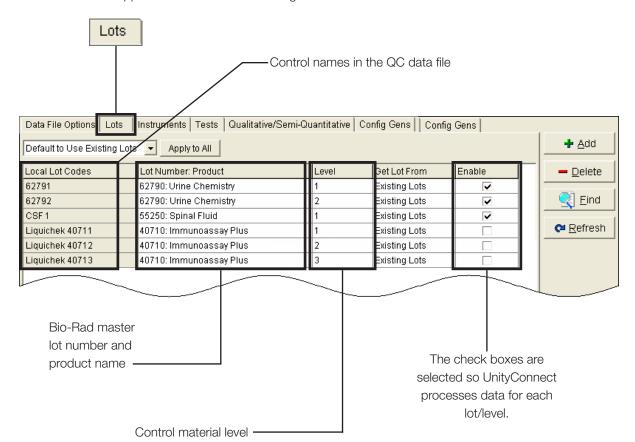


The **Lots** tab appears and shows the following information.

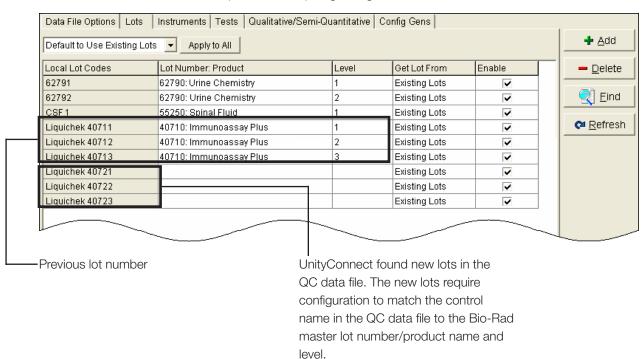


The **Lots** tab also shows the specific lots requiring configuration.

Notice the format of the lot number is consistent from the previous lot number to the new lot number.

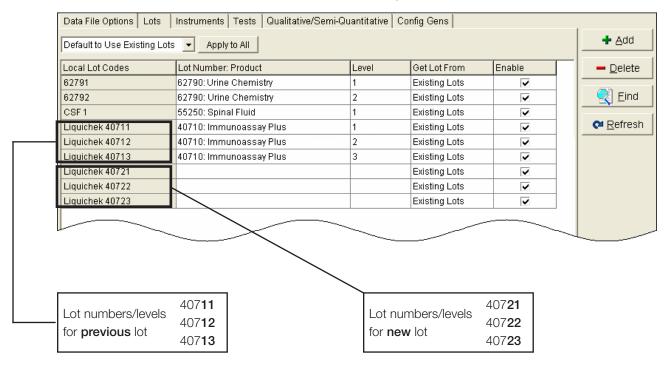

9 Select Default to Use Existing Lots.

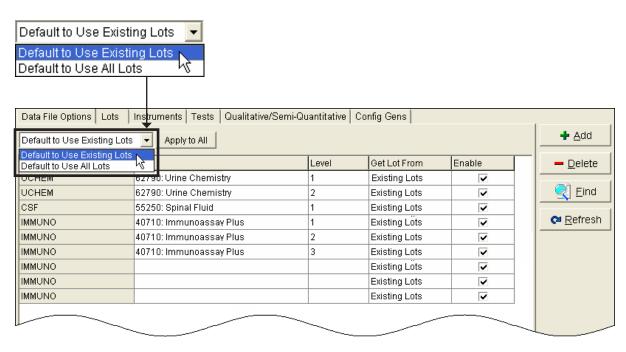


10 Click Apply to All

A message appears asking for confirmation.




11 Click <u>O</u>K

12 Click in the first field of the Lot Number: Product column requiring configuration and select the Bio-Rad master lot number and product name corresponding to the control in the Control Name column and the lot in the Local Lot Code column.

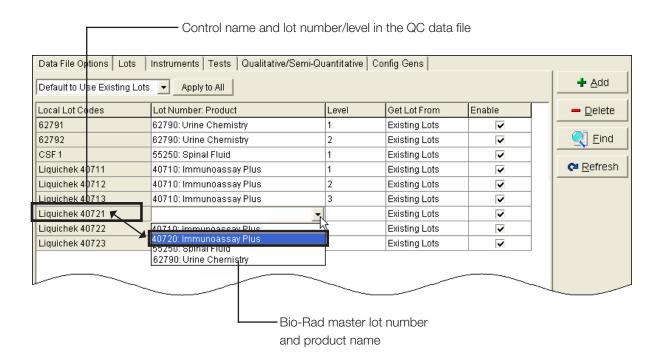




Note: The new lot number will not appear in the list of available lots if the old lot was not duplicated in your QC data management software.

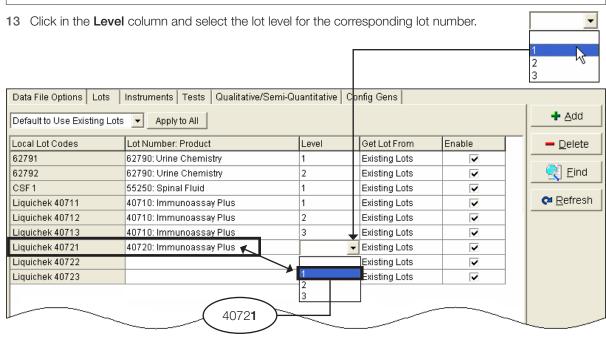

14 Make sure the Enable check box is selected for UnityConnect to process data for the lot.




Important: Click the **Enable** check box to clear it and UnityConnect will exclude the lot from current and future transforming. DO NOT delete the row. Deleting the row will cause the lot to appear again for future transforming.

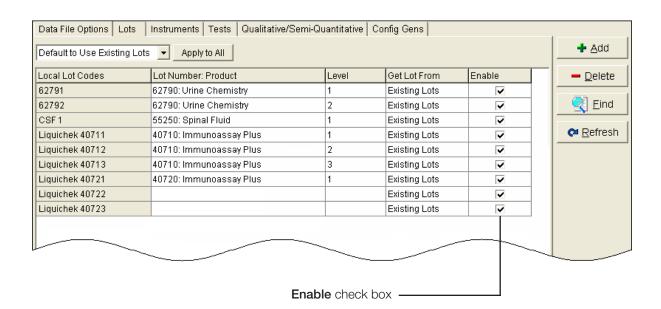

15 Configure any additional lots as necessary.

The following examples show configuring level 2 and level 3 for the new lot.

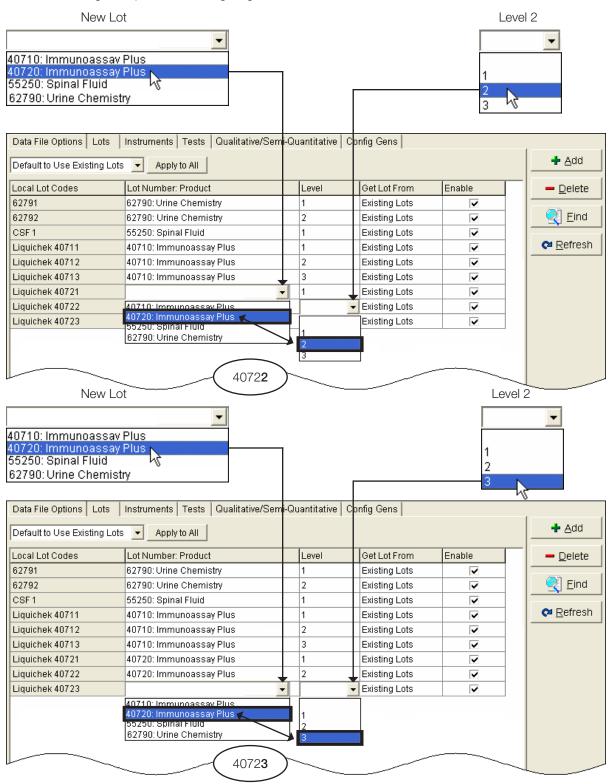



UnityConnect now has the information needed to transform data.

Server

A second message appears when transformation is complete.

Bio-Rad UnityConnect(TM) - Universal Flex File Transformer

File Transformer Options Logs Tools Help

Universal Flex File Transformer

Beginning Date

Specific Date

Specific Date

Last Point in File(s)

Notes

Please wait while your data is transformed...
Running Transformer Grabber completed.
Running Transformer Formatter Universal Flex File Transformer]...
Transformer Grabber completed.
Running Transformer Formatter Universal Flex File Transformer]...
Transformer Formatter Consoleted

Transformer Formatter Consoleted

Transformer Formatter Consoleted

Transformer Formatter Consoleted

CXProgram Files/Bio-RAD Al Laboratories\UnivyConnect\QC data file_june_15.txt contains a total of 5038 data records.

18 Read the information in the Notes section carefully and note the number of data records transformed.



Important: Configuring the new QC item is a one-time process. After the new item is configured, UnityConnect recognizes the item for all future transformations.

Database

Configuring a New Lot When Not Using QC Data Management Software for SPC Rule Evaluation

Switching to a new lot of control material is one of the most common reasons to complete new configuration in UnityConnect. Configuring UnityConnect for a new lot is a simple two-step process.

- 1 Configure the Import Settings to automatically create new lots in your QC data management software (this is a one-time setup).
- Complete new configuration in UnityConnect.

Configuring Import Settings


Configuring the import options/settings varies according to the type of QC data management software you use. See the appropriate section for more information:

- For Unity Desktop and Unity Real Time: See "Import Options-Not Using QC Data Management Software for SPC Rule Evaluation" on ref_here.
- For UnityWeb 2.0 and Unity Real Time online: See "Import Settings-Not Using QC Data Management Software for SPC Rule Evaluation" on ref. here.


Continue with the following section, "Completing New Configuration in UnityConnect" after configuring the import options/settings.

Completing New Configuration in UnityConnect

- 1 Start UnityConnect.
- 2 Set the date range for the data to transform. Data can be transformed for a particular date range, from the first point in the QC data file to the last point in the QC data file, or a combination of both.

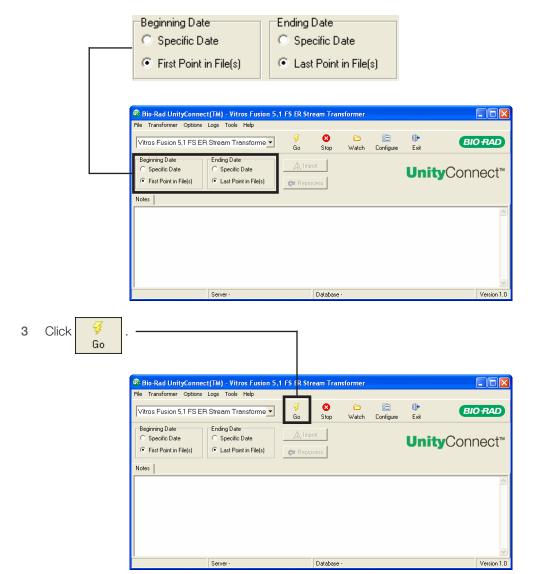

4 Select the QC data file to transform.

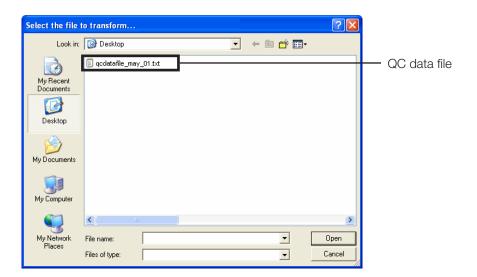


Tip: Bio-Rad recommends giving the QC data file a unique name so it can easily be identified.

Select the file to transform...


Look in: Desktop

QC data file.txt


My Recent Documents

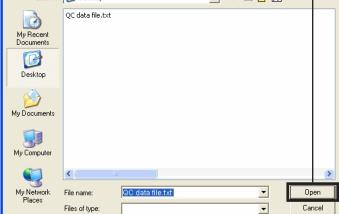

The QC data file name appears here after selecting it.

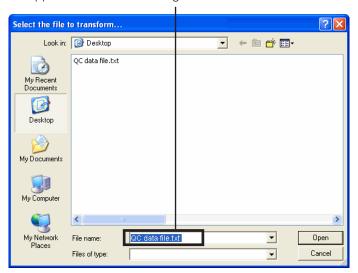


6 Click Open

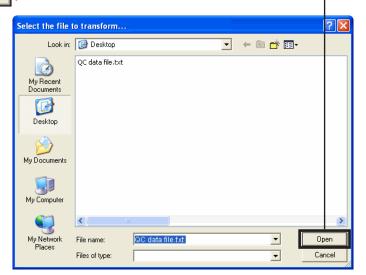


A message states the transformation is in progress.


Bio-Rad UnityConnect(TM) - Vitros Fusion 5,1 FS ER Stream Transformer
File Transformer Options Logs Tools Help

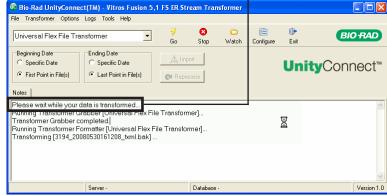

7 A message appears when UnityConnect finds the new lot. Read the note carefully.

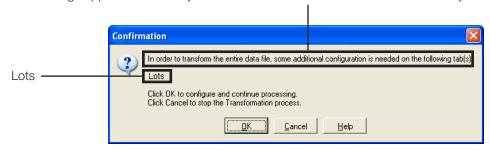


- **Important**: This is not an error message. The message is simply stating there are new codes in the QC data file that have not been matched to Bio-Rad codes in your QC data management software.
- 8 Click OK .
- 9 Log in to your QC data management software database if prompted.

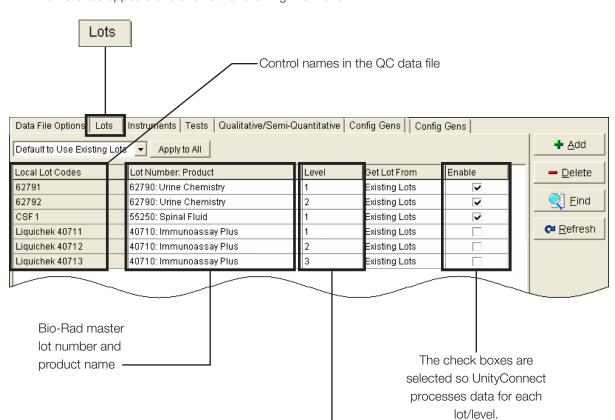


The **Lots** tab appears and shows the following information.

Control material level

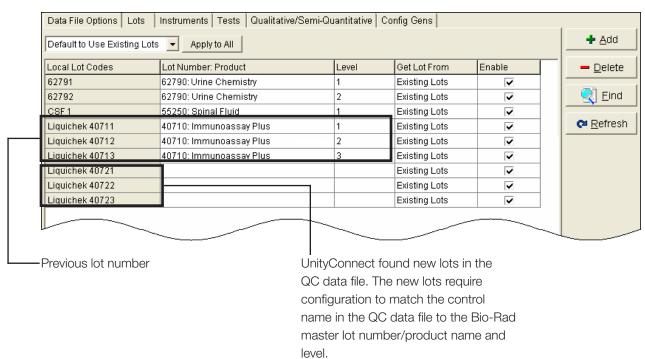


The **Lots** tab also shows the specific lots requiring configuration.

Notice the format of the lot number is consistent from the previous lot number to the new lot number.

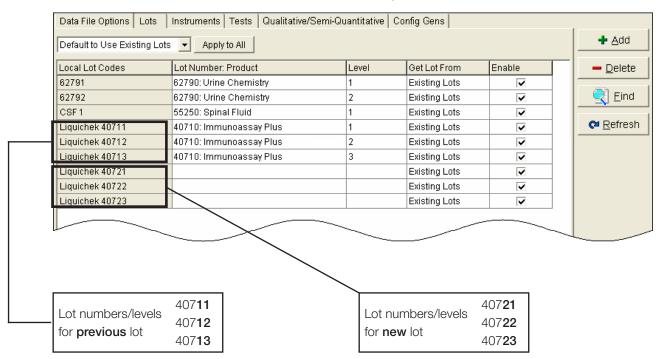


10 Select Default to Use All Lots.

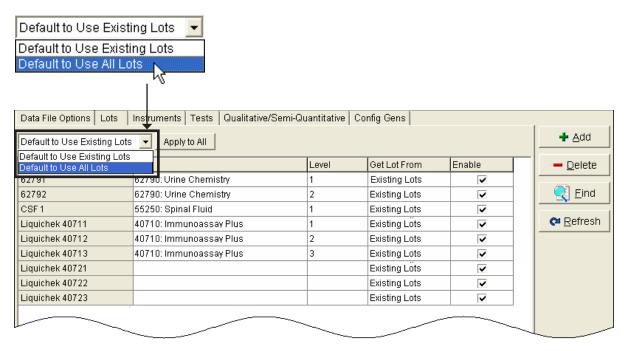


11 Click Apply to All


A message appears asking for confirmation.

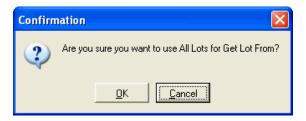

12 Click <u>Q</u>K

13 Click in the first field of the **Lot Number: Product** column requiring configuration and select the Bio-Rad master lot number and product name corresponding to the control in the **Control Name** column and the lot in the **Lot Number** column.

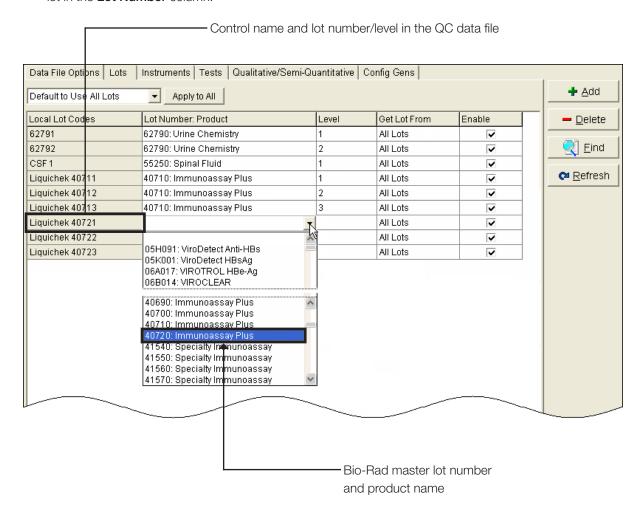

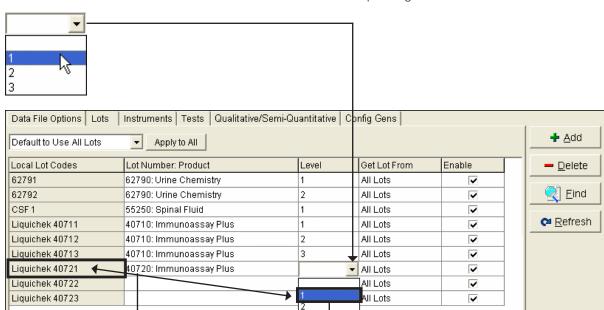

3


14 Click in the Level column and select the lot level for the corresponding lot number.

15 Make sure the **Enable** check box is selected for UnityConnect to process data for the lot.

4072**1**





Important: Click the **Enable** check box to clear it and UnityConnect will exclude the lot from current and future transforming. DO NOT delete the row. Deleting the row will cause the lot to appear again for future transforming.

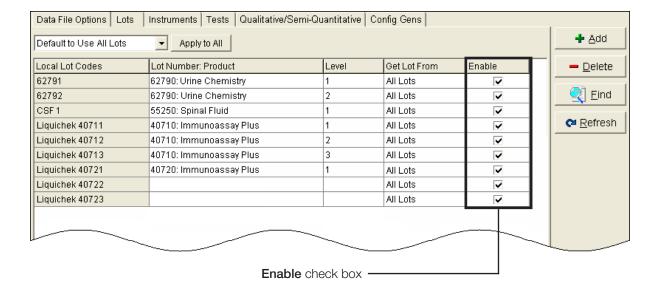

16 Configure any additional lots as necessary.

The following example shows configuring level 2 for the new lot.

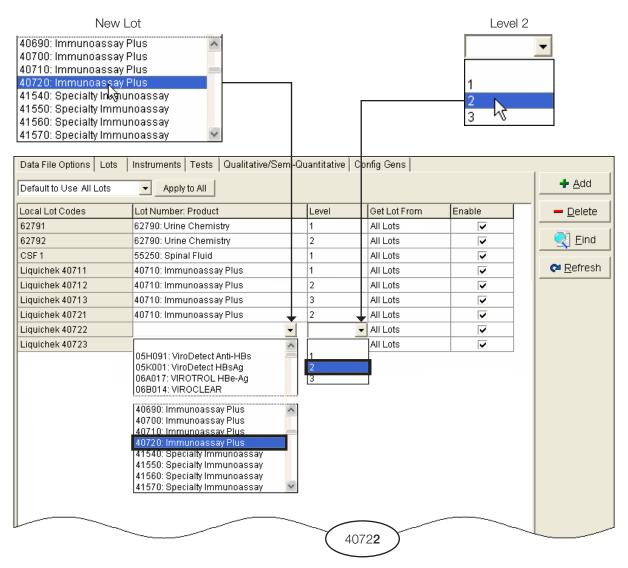

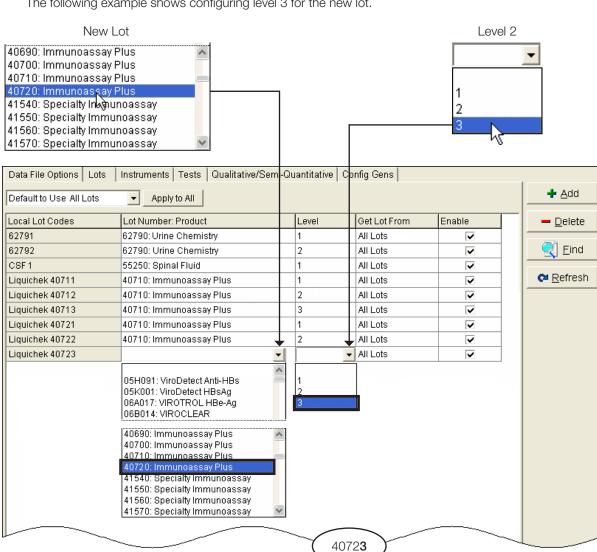


The following example shows configuring level 3 for the new lot.

17 Click Save when all lots are configured. ✓ Close 18 Click

UnityConnect now has the information needed to transform data.

A message indicates transformation is complete. -

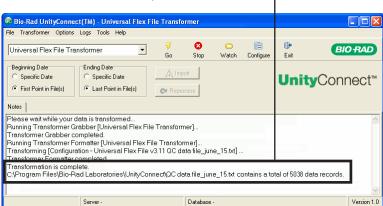

19 Read the information in the Notes section carefully and note the number of data records transformed.



Important: Configuring the new QC item is a one-time process. After the new item is configured, UnityConnect recognizes the item for all future transformations.

Other Types of Configuration


The following are possible types of information that may require configuration in UnityConnect:

- Data file options
 - See "Configuring New Data File Options" on ref_here.
- Instruments
 - See "Configuring a New Instrument" on ref_here.
- Tests
 - See "Configuring a New Test" on ref_here.

Configuring New Data File Options



Important: Due to the unique configuration of your customized UnityConnect software, talk to your Bio-Rad representative if you see the following message.

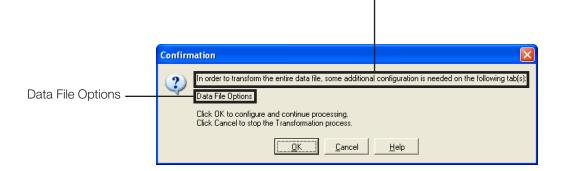


Configuring a New Instrument



Important: Configuring information in UnityConnect varies according to the information in your QC data file and the specifics of your customized UnityConnect software. Therefore, the example configuration information shown in this chapter may not be identical.


As UnityConnect processes the QC data file, a message appears when UnityConnect finds a new instrument.

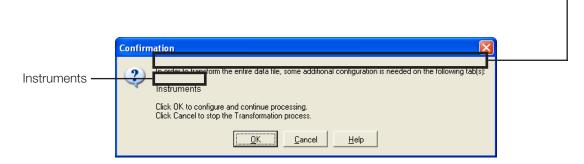




Important: This is not an error message. The message is simply stating there are new codes in the QC data file that have not been matched to Bio-Rad codes in your QC data management software.

1 Click OK

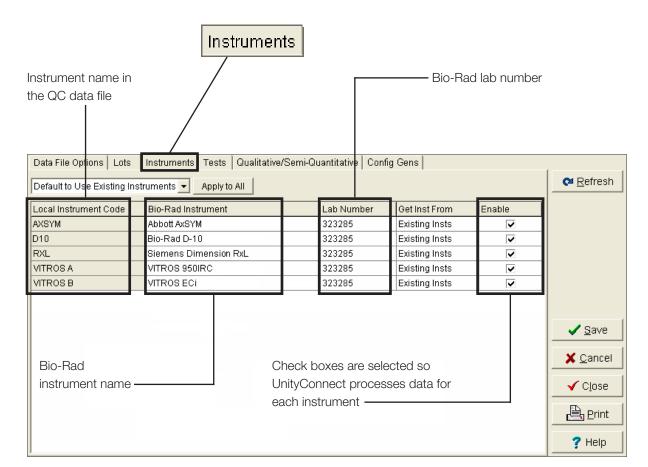


The **Instruments** tab appears and shows the following information.

The **Instrument** tab also shows the specific instrument requiring configuration.

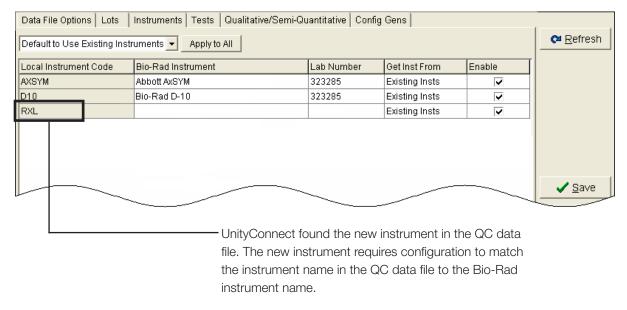

2 Select Default to Use All Instruments.

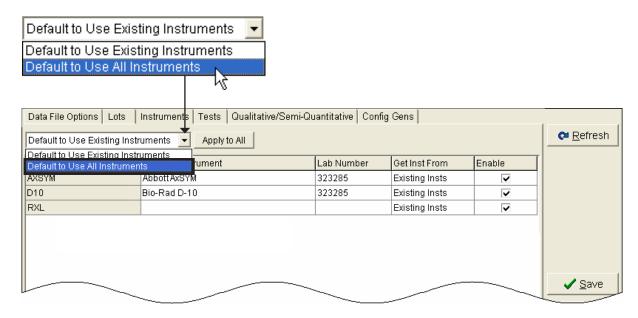


3 Click Apply to All

A message appears asking for confirmation.




4 Click <u>O</u>K .

5 Click in the field of the **Bio-Rad Instrument** column requiring configuration and select the Bio-Rad instrument name corresponding to the name in the **Local Instrument** column.

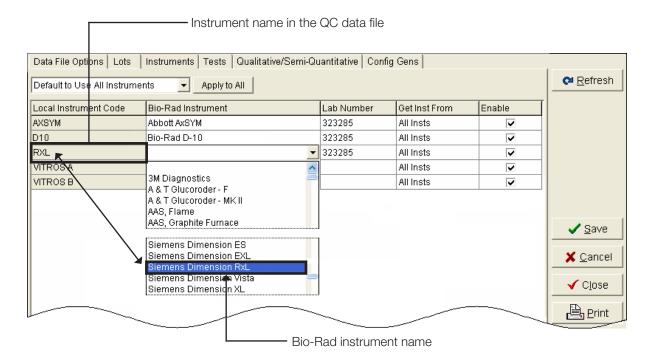


6 Click in the Lab Number column and select the lab number for the instrument.

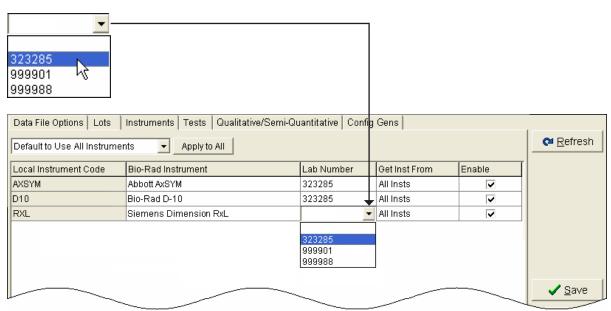



Important: Two instruments of the same type must be set up in separate lab numbers. A lab number can only be assigned by Bio-Rad. Contact your Bio-Rad QC Program Representative if you need additional lab numbers.

7 Make sure the Enable check box is selected for UnityConnect to process data for the instrument.
Click the Enable check box to clear it and UnityConnect will exclude the instrument when transforming.

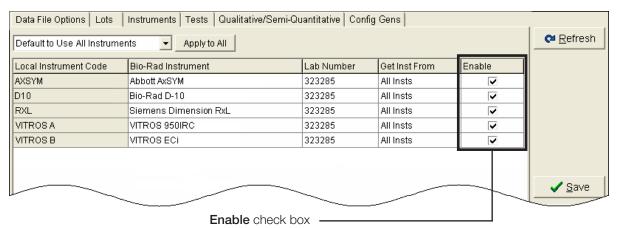

- 8 Click ✓ Save .
 9 Click ✓ Close .
 - UnityConnect now has the information needed to transform data.

A message indicates transformation is complete. 👨 Bio-Rad UnityConnect(TM) - Universal Flex File Transformer Transformer Options Logs Tools Help BIO RAD Universal Flex File Transformer Watch Configu Beginning Date Ending Date Specific Date C Specific Date **Unity**Connect™ First Point in File(s) Last Point in File(s) Please wait while your data is transformed...
Running Transformer Grabber [Universal Flex File Transformer]...
Transformer Grabber completed.
Running Transformer Formatter [Universal Flex File Transformer]...
Transforming [Configuration - Universal Flex File v3.11 QC data file_june_15.td]
Transformer Formatter completed.

10 Read the information in the **Notes** section carefully and note the number of data records transformed.




Important: Configuring the new QC item is a one-time process. After the new item is configured, UnityConnect recognizes the item for all future transformations.

Configuring a New Test




Important: Configuring information in UnityConnect varies according to the information in your QC data file and the specifics of your customized UnityConnect software. Therefore, the example configuration information shown in this chapter may not be identical.

As UnityConnect processes the QC data file, a message appears when UnityConnect finds a new test.

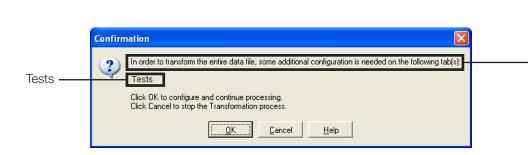




1 Click OK .

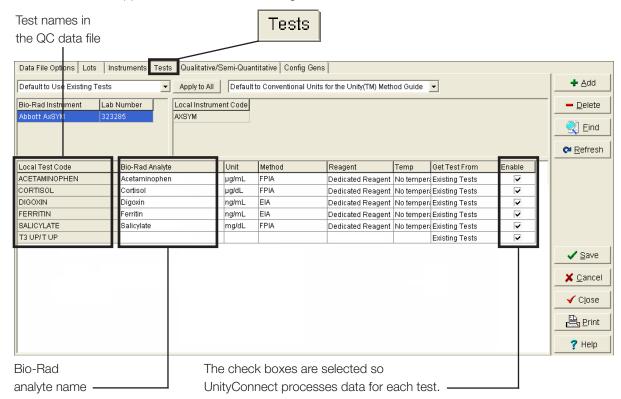


The **Tests** tab appears and shows the following information.

The **Tests** tab also shows the specific test requiring configuration.

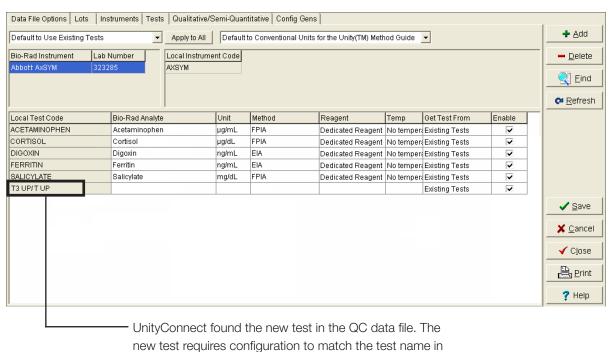


the QC data file to the Bio-Rad analyte name.


2 Click in the field of the **Bio-Rad Analyte** column requiring configuration and select the Bio-Rad analyte name corresponding to the name in the **Local Test Code** column.

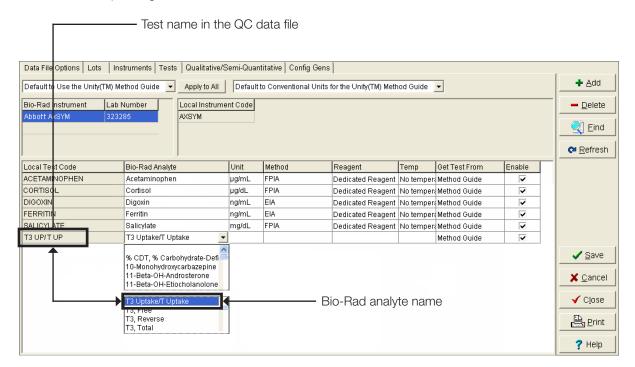


3 Click in the **Unit** column and select the unit for the test.

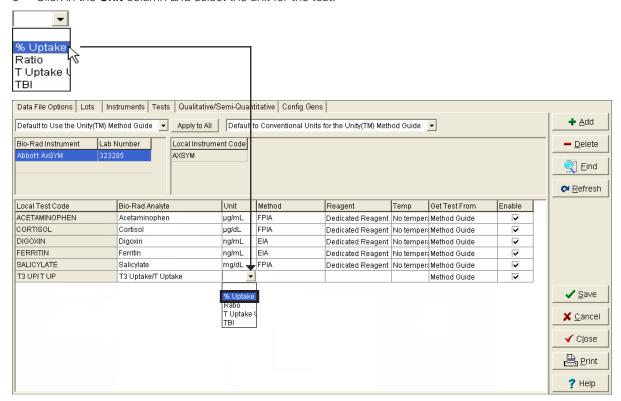

4 Click in the **Method** column and select the method for the test.

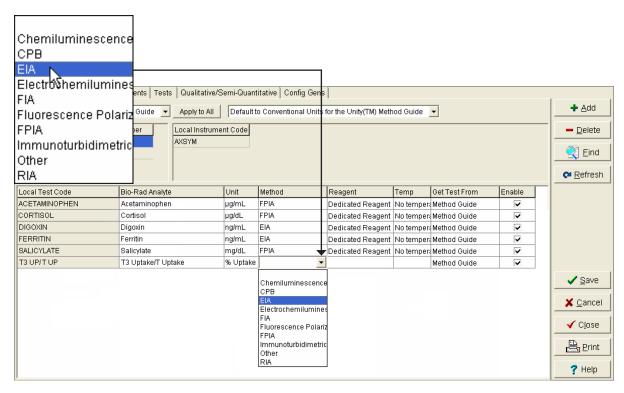

5 Click in the **Reagent** column and select the reagent for the test.

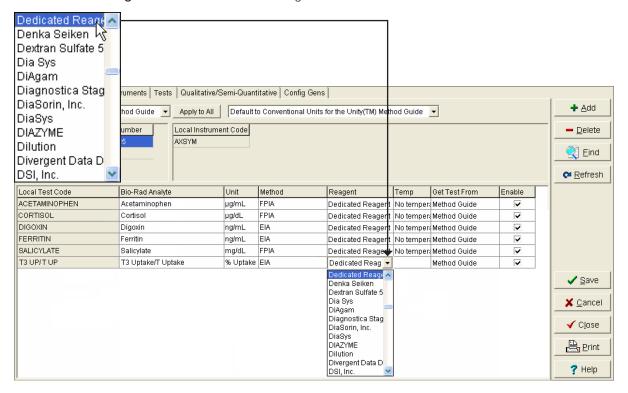

No Tempe Data File Options Lots Instruments Tests Qualitative/Semi-Quantitative Config Gens
 ♣ Add
 Default to Use the Unity(TM) Method Guide 🔻 Apply to All Default to Conventional Units for the Unity(TM) Method Guide Bio-Rad Instrument Local Instrument Code - Delete Lab Number AXSYM <u>Eind</u> Bio-Rad Analyte Method Get Test From Enable Local Test Code Unit Reagent Temp ACETAMINOPHEN Acetaminophen μg/mL FPIA Dedicated Reagent No temper: Method Guide V CORTISOL Cortisol µg/dL FPIA Dedicated Reagent No temper Method Guide DIGOXIN EIA Digoxin ng/mL Dedicated Reagent No temper: Method Guide FERRITIN Ferritin ng/mL EIA Dedicated Reagent No temper: Method Guide SALICYLATE Salicylate mg/dL FPIA Dedicated Reagent No temp V T3 UP/T UP T3 Uptake/T Uptake % Uptake EIA Dedicated Reagent No Terr ▼ Method Guide ✓ Save 💢 <u>C</u>ancel ✓ Close Print | ? Help

6 Click in the **Temperature** column and select the temperature for the test.




Note: Temperature applies to enzymes only. For all other analytes, **No temperature** is the only available option.


7 Make sure the **Enable** check box is selected for UnityConnect to process data for the test.



Important: Click the **Enable** check box to clear it and UnityConnect will exclude the test from current and future transforming. DO NOT delete the row. Deleting the row will cause the test to appear again for future transforming.

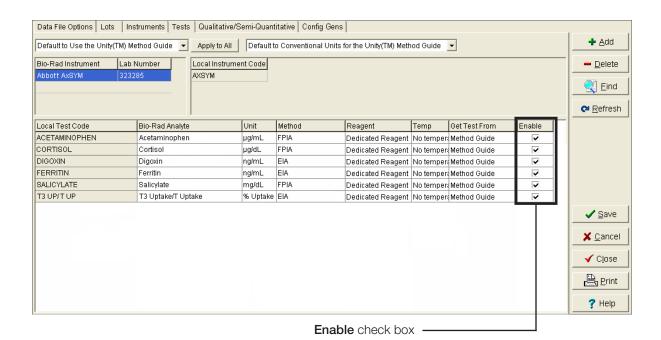

- 8 Click ✓ <u>Save</u> .
 9 Click ✓ <u>Close</u> .
 - UnityConnect now has the information needed to transform data.

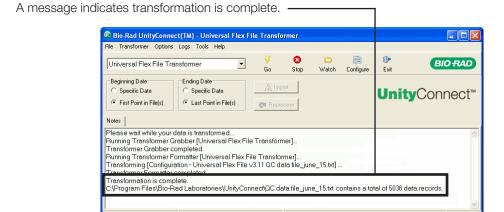

10 Read the information in the Notes section carefully and note the number of data records transformed.



Important: Configuring the new QC item is a one-time process. After the new item is configured, UnityConnect recognizes the item for all future transformations.

VITROS Slide Generation Numbers and UnityConnect

In This Chapter

Overview	201
Slide Generation Number Decision Chart	202
Example Config Gens Tab	203
Managing Slide Generation Numbers in UnityConnect	205
Slide Generation Worksheet	212
Changing SPC Rules for a New VITROS Test	213

Overview



Important: The VITROS ECi does not use slide generation numbers. Therefore, the information in this chapter does not pertain to the VITROS ECi instrument.

Slide Generation Number Decision Chart

Use the following chart as a guide for determining how to address slide generation numbers in UnityConnect.

VITROS Instrument Type	Uses slide generation numbers?	Slide generation numbers located in QC data file?	Manual setup slide generation numbers in UnityConnect required?
VITROS ECi	N/A	N/A	N/A
VITROS Fusion 5, 1 with UnityConnect Instrument Transformer	Yes	Yes	No
VITROS Fusion 5, 1 with UnityConnect LIS Transformer	Yes	No	Yes
VITROS 250	Yes	No	Yes
VITROS 350	Yes	No	Yes
VITROS 950	Yes	No	Yes

With the exception of the VITROS Fusion 5,1 using a UnityConnect instrument transformer, QC data files do not contain information about VITROS slide generation numbers so UnityConnect cannot automatically prompt you to update or add slide generation numbers. Therefore, all slide generation number changes must be made manually in UnityConnect before transforming data. If the slide generation number is not updated in UnityConnect before transforming data, the data for the new slide generation number is reported on an incorrect slide generation number.

Slide generation numbers are automatically updated in your QC data management software and UnityConnect if you are using a VITROS Fusion 5,1 with a UnityConnect instrument transformer. Therefore, you do not need to manually change slide generation numbers in UnityConnect.




Note: The VITROS Fusion 5,1 can use both dry tests and wet tests. As a general rule, dry tests have slide generation numbers but wet tests do not.

Example Config Gens Tab

The Config Gens tab contains information about slide generation numbers manually set up in UnityConnect.



Important: This chapter applies only to certain customers. Please refer to the "Slide Generation Number Decision Chart" on page 202 to determine if your circumstances require manual set up of slide generation numbers in UnityConnect.

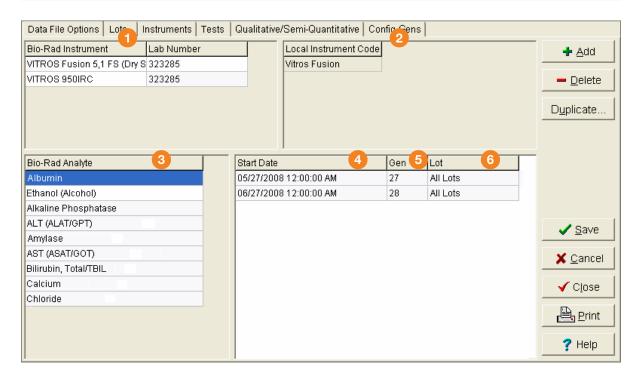

Bio-Rad Instrument and Lab Number

The Bio-Rad instrument and corresponding lab number using slide generation numbers as defined on the **Instruments** tab. If more than one instrument is in this list, click the Bio-Rad instrument or lab number to view the related **Local Instrument**, **Bio-Rad Analyte**, **Start Date**, **Gen**, and **Lot** information for the instrument.

2 Local Instrument

The instrument designation in the QC data file.

Bio-Rad Analyte

The corresponding Bio-Rad analyte names for the tests configured for the selected instrument.

Start Date

The date and time the slide generation number was first used

Gen

The current slide generation number.

6 Lot

The lot number which the current slide generation number is associated with.

Best Practices—Manage Manual VITROS Slide Generation Changes



Important: The Unity Interlaboratory Program uses the VITROS slide generation numbers to determine the consensus group. Make sure the VITROS slide generation number is correct for each test to ensure accurate Unity Interlaboratory Reports.

It is important to establish a regular process for tracking slide generation numbers. Use the Slide Generation Worksheet on page 212 or a similar process to track the VITROS slide generation numbers.

Managing Slide Generation Numbers in UnityConnect




Important: This chapter applies only to certain customers. Please refer to the "Slide Generation Number Decision Chart" on page 202 to determine if your circumstances require manual changing of slide generation numbers in UnityConnect.

All slide generation changes must be made manually in UnityConnect before transforming data. If the slide generation number is not updated in UnityConnect before transforming data, the data for the new slide generation number is reported on an incorrect slide generation number.


Each test setup is an individual combination of test items. Therefore, when you change to a new slide generation number, you create a new test.

- 1 Start UnityConnect.
- 2 Select the appropriate transformer from the **Select a transformer** list.
- 3 Click Configure .


4 Click the Config Gens tab.

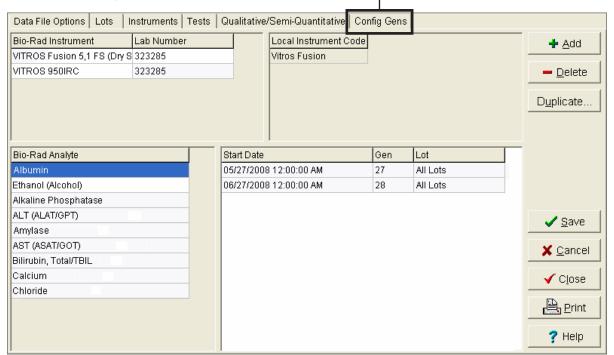


5 Click the Bio-Rad instrument or lab number with the slide generation number to change.

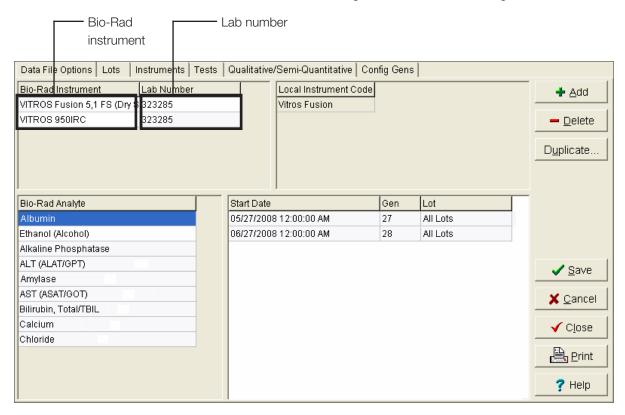

The tests for the instrument or lab number you selected appear in the Bio-Rad Analyte column.

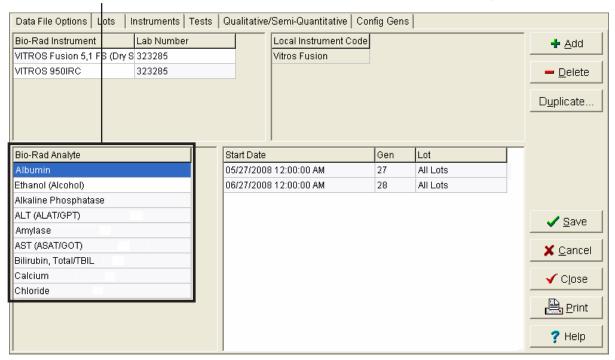



6 Click the test in the **Bio-Rad Analyte** column for which you want to change the slide generation number.

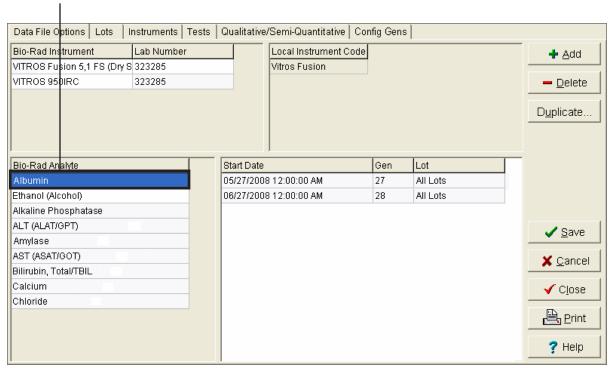




Tip: The selected test appears highlighted in blue.

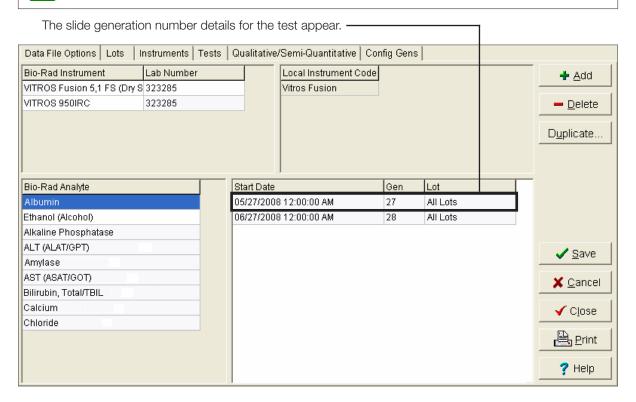

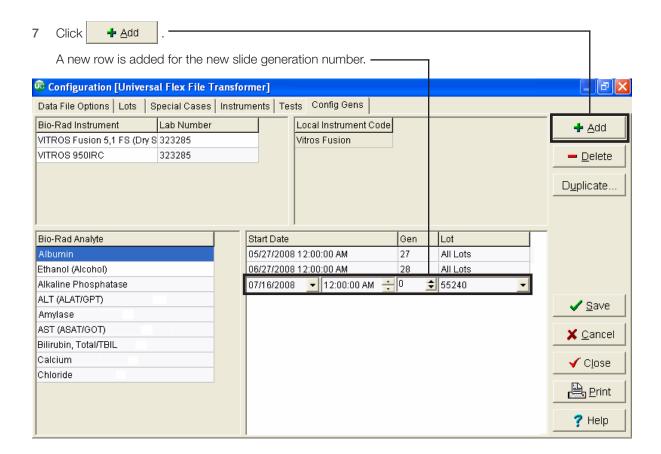

8 Click in the Start Date field.


The field becomes active.

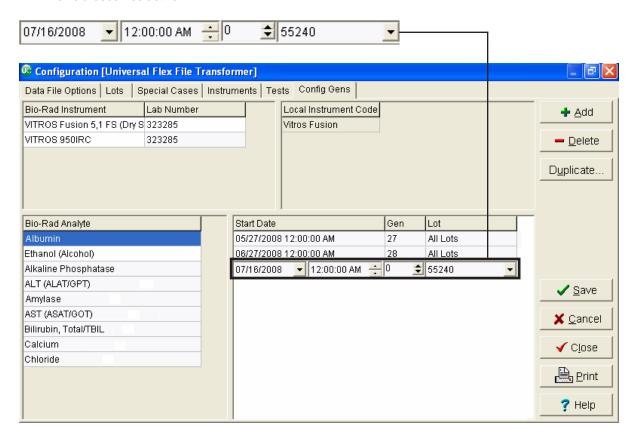


4


July, 2008 Sun Mon Tue Wed Thu Fri Sat

Click again in the Start Date field to display a calendar. Select a start date for the slide generation number from the calendar.

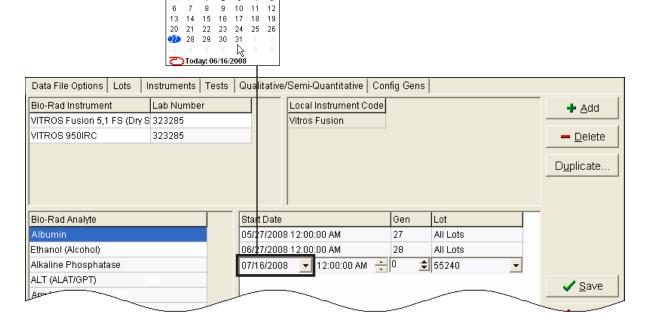


10 Click in the hours, minutes, and seconds fields and type the time for the slide generation number change.

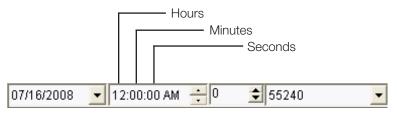


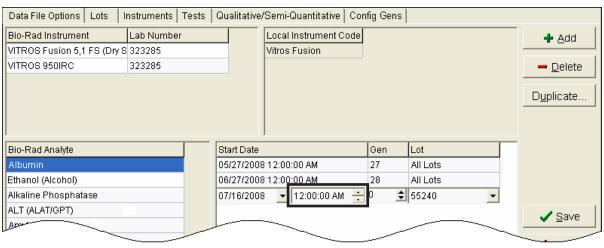

11 Click in the **Gen** field and type the slide generation number.

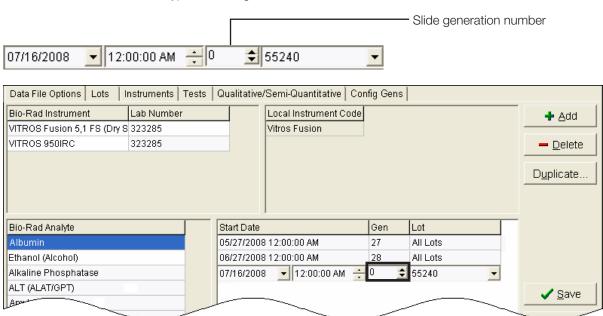

12 Click in the Lot column and select the lot number for the slide generation number or select All Lots.

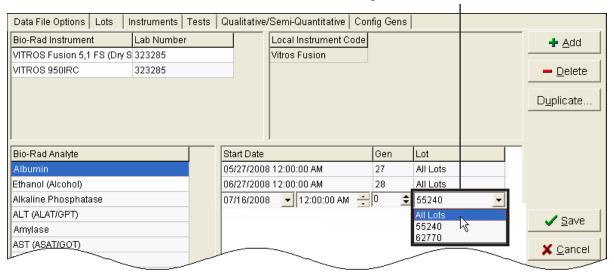

- 13 Click ✓ Save
- 14 Repeat steps 6–14 for each test with a slide generation number change.
- 15 Click ✓ Close when all slide generation numbers are configured.

Slide Generation Worksheet

•	$\overline{}$
П	
п	
п	

Tip: Use this worksheet to track your slide generation numbers. You can make additional copies as needed.

Lot Number	Analyte	Instrument	Gen #	Date in Use

Changing SPC Rules for a New VITROS Test



Important: Each time you make a slide generation number change, UnityConnect creates a new test with the new slide generation number. The new test uses the default SPC rules. You must manually change the SPC rules for the test in your QC data management software if you use customized SPC rules for the test or if you do not use any SPC rules. This applies even if the SPC rules are turned off at the lot level.

See the appropriate section based on your QC data management software:

- Unity Desktop and Unity Real Time (see page 213)
- UnityWeb 2.0 and Unity Real Time online (see page 214)

For Unity Desktop and Unity Real Time

- 1 Start your QC data management software.
- 2 Select the test in the navigation tree.
- 3 Click SPC Rules
- 4 Select an option (Reject, Warn, or Off) for each SPC rule.

The status of each rule is indicated in the Status column using the symbols shown below:



Reject



Warn



Off

OK



Note: Click Disable SPC Rules to set all the rules to off.

5 Click



For UnityWeb 2.0 and Unity Real Time online

- 1 Log in to your QC data management software.
- 2 Click the **Configure** tab.
- 3 Click Rules/Settings.
- 4 Click the Rules tab.
- 5 Select the Lab or Panel option.
- 6 Select the lab number or panel name from the Lab or Panel list.
- 7 If using the Lab option, select the lot number from the Lot list.
- 8 Select the test from the **Test** list.
- 9 Click Rules.
- 10 Select an option (Reject, Warn, or Off) for each SPC rule.

The status of each rule is indicated in the **Status** column using the symbols shown below:



Reject



Warn



Off



Note: Click Disable SPC Rules to set all the rules to off.

11 Click Save.

Viewing the Rejection Log for Desktop Software

In This Chapter

Overview	215
Using the Rejection Log	216
Viewing the Rejection Log	
Deleting the Rejection Log	
Rejection Log Messages	

Overview

Rejected data is transformed data that was not imported into your QC data management software. The Rejection Log shows the details of all rejected data and provides a description about why the data was rejected.



Important: The Rejection Log does not show data rejected due to an SPC rule violation, but rather only transformed data that was not imported into your QC data management software.

In most cases of rejected data, no action is required. However, in some instances an error message is an indication that configuration adjustments are required. After the configuration adjustments are made, the QC data file must be transformed again.

This chapter provides information on how to review the Rejection Log and determine the appropriate action based on the error message.


Before You Begin



This chapter is only for customers using UnityConnect with any of the following **desktop** software:

Unity DesktopUnity Real Time

Using the Rejection Log

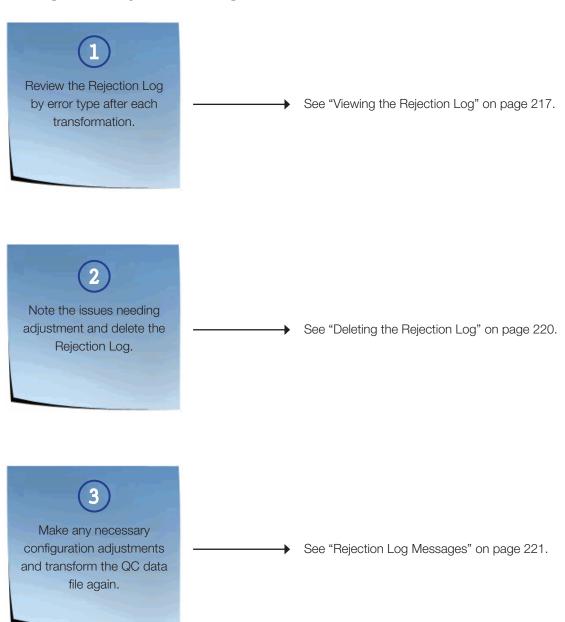

Viewing the Rejection Log


Best Practices—Review the Rejection Log by Error Type After Each Transformation

- Start UnityConnect.
- 2 Click the Logs menu and then click Rejection Logs.



The **Database Login** dialog box appears.

- 3 Type your **User Password** and click **OK**.
- 4 Click the arrows located to the right of the **Start Entered Date** and/or **End Entered Date** to select another date range to view.

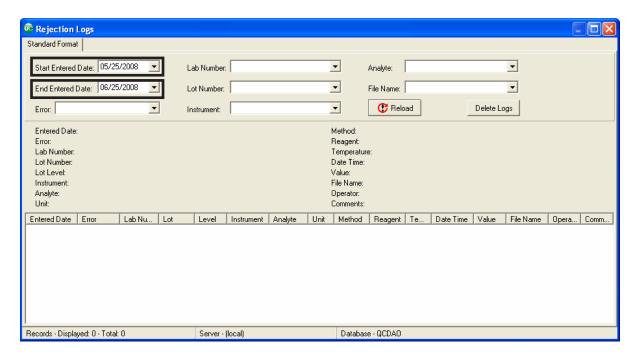




Tip: The **End Entered Date** defaults to the current date. The **Start Entered Date** is exactly one month earlier. You typically do not need to change these dates.

5 Select the type of error from the **Error** list.



Tip: It is easiest to view the Rejection Log by each error type if the Rejection Log contains a large amount of information.

- 6 Select a lab number from the **Lab Number** list or select **All**.
- 7 Select a lot number from the Lot Number list or select All.
- 8 Select an instrument from the **Instrument** list or select **All**.
- 9 Select an analyte from the **Analyte** list or select **All**.
- 10 Select a file name from the File Name list or select All.

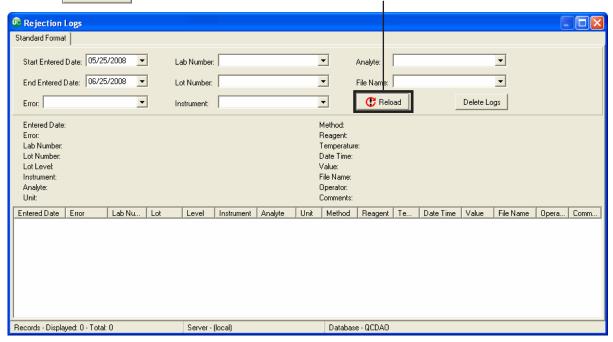


The Rejection Log appears according to the criteria you selected.

12 Review the Rejection Log paying special attention to the Error column.

The **Error** column indicates why the data was not imported.



Note: See "Rejection Log Messages" on page 221 for a description of messages and information about resolving the message.

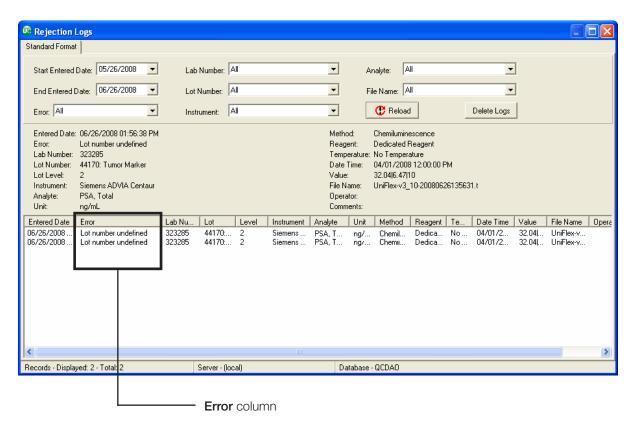


The error type should match the criteria you selected in step 5.

Deleting the Rejection Log

Best Practices—Note the Issues Needing Adjustment and Delete the Rejection Log

Bio-Rad recommends deleting the Rejection Log after each review to eliminate confusion.

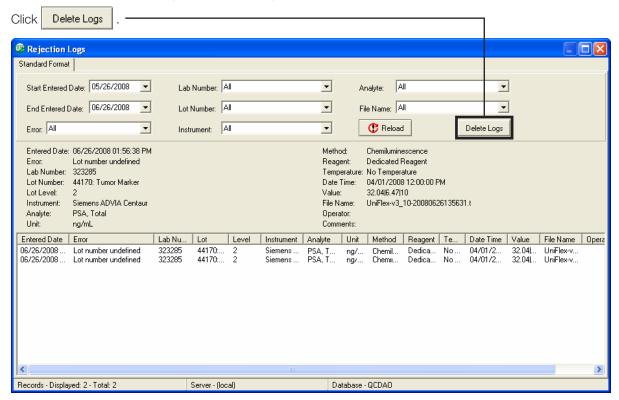

Rejection Log Messages

The following are common Rejection Log messages. See the appropriate section for a description of the message and information about resolving the message.

- Create new tests if necessary is disabled (see page 221)
- Data entry locked for this test (see page 225)
- Date out of sequence (see page 225)
- Lab closed (see page 227)
- Lot closed (see page 227)
- Lot expired (see page 228)
- Lot number undefined (see page 228)
- Result invalid (see page 236)
- Time out of sequence (see page 237)

Error: Create new tests if necessary is disabled



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.


Description:

The **Create New Tests if Necessary** check box is not selected on the **Import** dialog box. UnityConnect is not able to create the new test in your QC data management software so the associated data cannot be imported.


Solution:

Make sure the Create New Tests if Necessary check box is selected on the Import dialog box.


- 1 Start UnityConnect.
- 2 Click the Options menu and then click Settings.

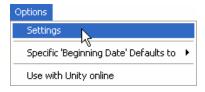


3 Click the Import Options tab.

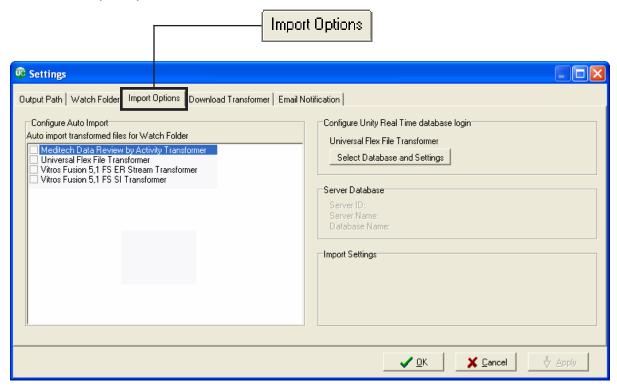

4 Select the transformer for which you want to configure import options.

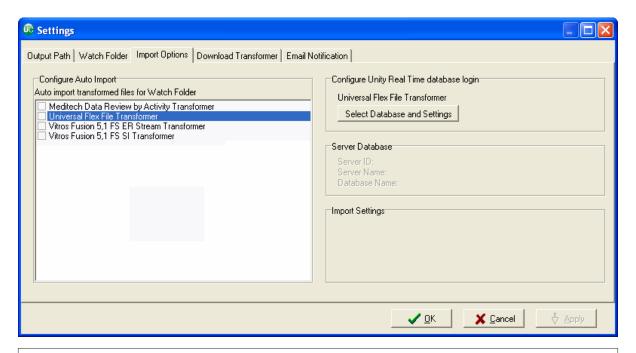

Note: Do not select the check box for the transformer unless you want to set up automatic import. See "Configuring Automatic Import" on page 67 for more information.

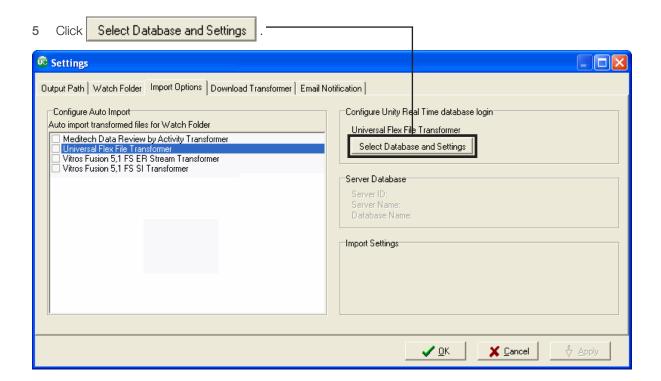

6 Make sure the Create New Lots if Necessary check box is selected.

The following illustration shows recommended import settings for customers not using their QC data management software for SPC rule evaluation.

Create New Tests if Necessary check box is selected. -👨 Unity Real Time Database Login Transformer: Import Settings Universal Flex File Transformer Decimal Point: SQL Server Database ♦ Minute(s) Run Length: 15 Connect using Windows Authentication Create New Lots if Necessary Create New Tests if Necessary SQL Server ID: Sort Import Files SQL Server Password: Redundant Data Filter ▼ Disable All Rules Violation • SQL Server Name: Delete Rejection Log Records after 60 🖨 Day(s) • Database Name: Operator • User Name: User Password: **√** <u>0</u>K X Cancel

- 7 Type your password in the **User Password** field.
- 8 Click <u>**V** 0</u>K.
- 9 Click VOK again.

A new lot is automatically created when UnityConnect finds a new lot in your QC data file.

10 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

11 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected in the **Import** dialog box, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

12 Delete the Rejection Log.

Error: Data entry locked for this test



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

UnityConnect is not able to import data for a test when a data entry dialog box is open in any user's QC data management software.

Solution:

Wait until the data entry dialog box is closed.

- 1 Transform the QC data file again.
 - Data that was not imported during the first import now imports into your QC data management software.
- 2 Review the Rejection Log.
 - If the **Redundant data filter** check box is **not** selected on the **Import** dialog box, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.
- 3 Delete the Rejection Log.

Error: Date out of sequence



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

This error occurs when you attempt to import data into your QC data management software out of chronological order. You can only import data with a testing date **after** the last date of data already in your QC data management software.

Typical Example:

You had a problem with your QC data on April 15, 16, and 17. In attempting to investigate the issue, you transformed a QC data file containing QC data for those three days (only) and UnityConnect automatically imported the data into your QC data management software.

On April 30 you transformed a QC data file containing all QC data for the month of April.

- All data points from April 1 through April 14 were rejected with a "Date out of sequence" error.
- Since data for April 15, 16, and 17 already exists, those data points were rejected as "Redundant data."
- All data points from April 18–30 imported successfully.




Solution-Unity Desktop and Unity Real Time

In order to import data points from April 1 through April 14, use the **Delete Range of Data** feature to delete data from April 15–30 before transforming the entire QC data file for April again.



Important: Deleting a range of data permanently removes the data from the software. Bio-Rad recommends exporting the data before deleting. See "Exporting Data" in the Reference Guide or Online Help for more information.

- 1 Start your QC data management software.
- 2 Click Test .
- 3 Click Delete Range of Data
- 4 Click the arrow located to the right of the From date and select a beginning date for the range of data.

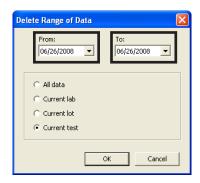

- 5 Click the arrow located to the right of the **To** date and select an ending date for the data.
- 6 Select an option for the scope of data to delete.
 - All data
 - Current lab
 - Current lot
 - Current test
- 7 Click OK .

A message appears asking for confirmation.

- 8 Click Yes .
- 9 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

10 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected in the Import Settings, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

11 Delete the Rejection Log.

Error: Lab closed



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

The lab number corresponding to the data is closed in your QC data management software.

Solution:

You cannot import data into your QC data management software for a closed lab. Open the lab number in your QC data management software.

- 1 Start your QC data management software.
- 2 Click in Unity Desktop or Unity Real Time.
- 3 Select the lab you want to open in the Closed labs list.
- 4 Click Open Lab .

The lab number moves to the bottom of the **Open labs** list.

- 5 Click X Close .
- 6 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

7 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected on the **Import** dialog box, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

8 Delete the Rejection Log.

Error: Lot closed



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

The lot number corresponding to the data is closed in your QC data management software.

Solution:

You cannot import data into your QC data management software for a closed lot. Open the lot number in your QC data management software.

1 Start your QC data management software.

- 2 Click in Unity Desktop or Unity Real Time.
- 3 Select the lot you want to open in the Closed lots list.
- 4 Click Open Lot .

The lot number moves to the bottom of the Open lots list.

- 5 Click X Close .
- 6 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

7 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected on the **Import** dialog box, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

8 Delete the Rejection Log.

Error: Lot expired



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

The lot number corresponding to the data is expired in your QC data management software.

Solution:

You cannot import data into your QC data management software for an expired lot. No further action is required.

Error: Lot number undefined




Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

You have configured a lot in UnityConnect using the **Default to All Lots** option, and the **Create new lots** check box is not selected on the **Import** dialog box.


Solution:

The solution to this error depends on your use of your QC data management software for SPC rule evaluation.

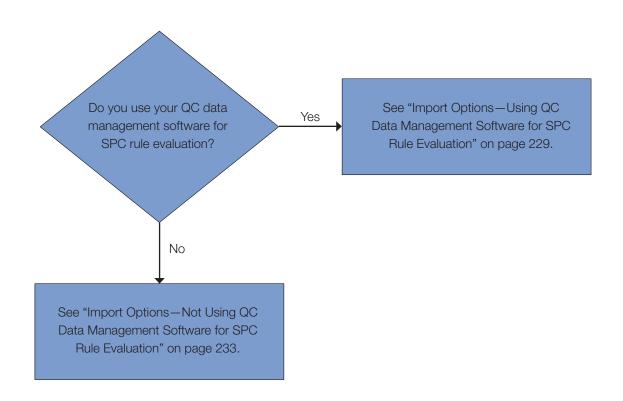


Import Options—Using QC Data Management Software for SPC Rule Evaluation


- 1 Start UnityConnect.
- 2 Click the **Options** menu and then click **Settings**.

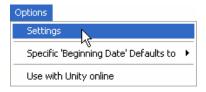


3 Click the Import Options tab.

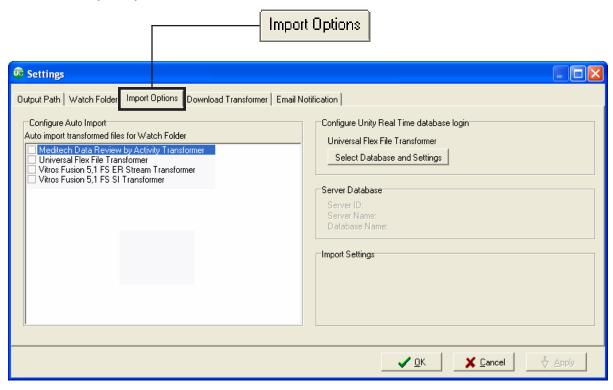

4 Select the transformer for which you want to configure import options.

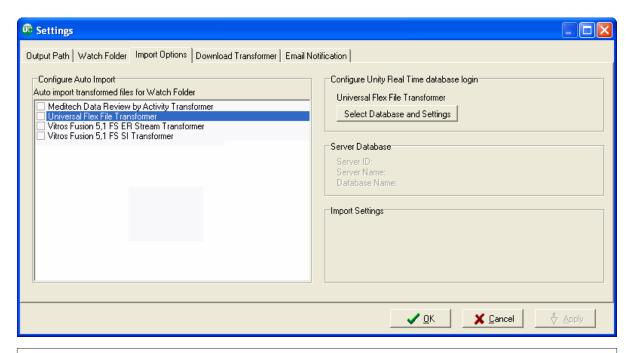

Note: Do not select the check box for the transformer unless you want to set up automatic import. See "Configuring Automatic Import" on page 67 for more information.

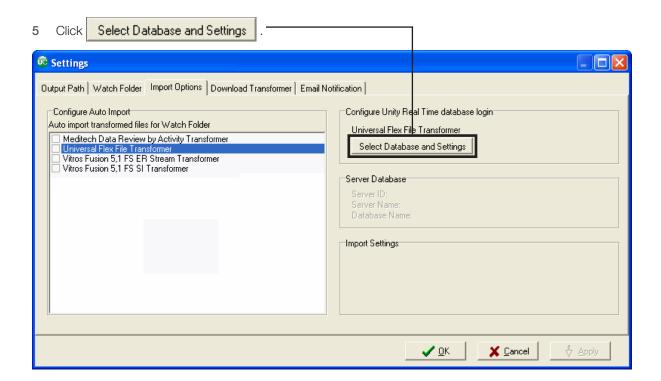

6 Make sure the **Create New Lots if Necessary** check box is not selected.

The following illustration shows recommended import settings for customers using their QC data management software for SPC rule evaluation.

The Create New Lots if Necessary check box is not selected. 👨 Unity Real Time Database Login Transformer: Import Settings Universal Flex File Transformer Decimal Point: SQL Server Database ♦ Minute(s) Run Length: 15 Connect using Windows Authentication Create New Lots if Necessary Create New Tests if Necessary SQL Server ID: Sort Import Files SQL Server Password: Redundant Data Filter Disable All Rules Violation • SQL Server Name: Delete Rejection Log Records after 60 🕏 Day(s) • Database Name: Operator • User Name: User Password: **√** <u>0</u>K X Cancel


This ensures UnityConnect will not create new lots when transforming and overwrite your current SPC rule settings in your QC data management software.

- 7 Type your password in the **User Password** field.
- 8 Click <u>\(\lambda \) O</u>K .
- 9 Click VOK again.


To duplicate your current SPC rule settings from the current lot to the new lot, leave the Create New Tests if Necessary check box blank and duplicate the lot in your QC data management software before transforming data for the lot for the first time. See "Configuring a New Lot" on page 164 for more information.

Import Options—Not Using QC Data Management Software for SPC Rule Evaluation


- 1 Start UnityConnect.
- 2 Click the **Options** menu and then click **Settings**.

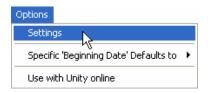


3 Click the **Import Options** tab.

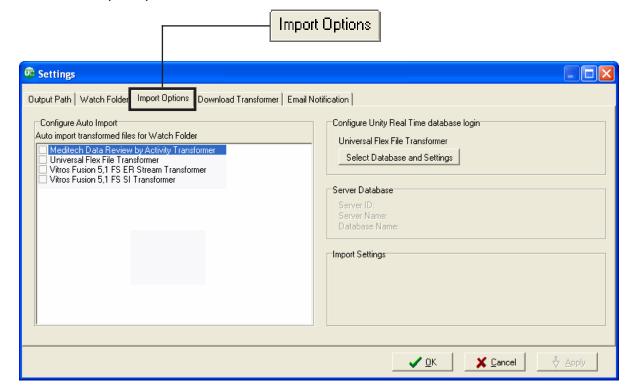

4 Select the transformer for which you want to configure import options.

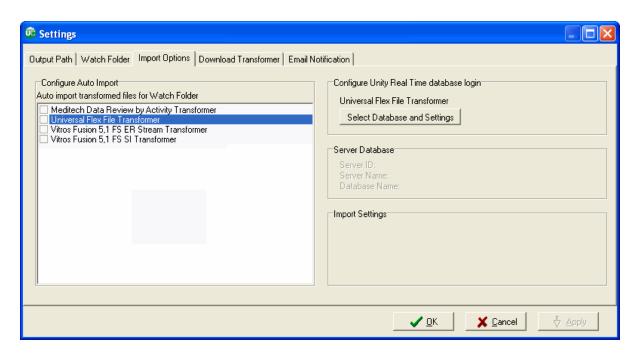

Note: Do not select the check box for the transformer unless you want to set up automatic import. See "Configuring Automatic Import" on page 67 for more information.

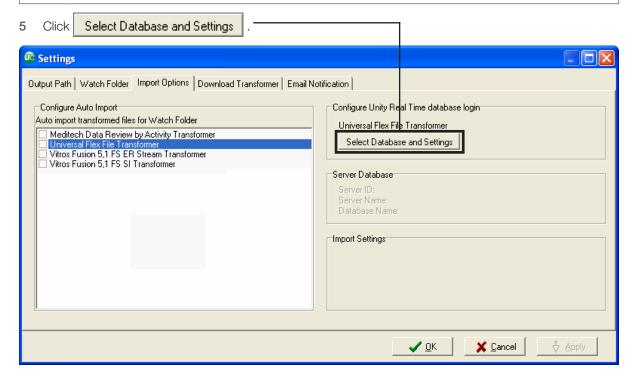

6 Make sure the **Create New Lots if Necessary** check box is selected.

The following illustration shows recommended import settings for customers not using their QC data management software for SPC rule evaluation.

Create New Tests if Necessary check box is selected.-👨 Unity Real Time Database Login Transformer: Import Settings Universal Flex File Transformer Decimal Point: SQL Server Database Run Length: 15 Connect using Windows Authentication Create New Lots if Necessary ▼ Create New Tests if Necessary SQL Server ID: Sort Import Files SQL Server Password: Redundant Data Filter ▼ Disable All Rules Violation • SQL Server Name: Delete Rejection Log Records after 60 🕏 Day(s) • Database Name: Operator User Name: • User Password: **√** <u>0</u>K 💢 <u>C</u>ancel

- 7 Type your password in the **User Password** field.
- 9 Click VOK again.

A new lot is automatically created when UnityConnect finds a new lot in your QC data file.

Error: Result invalid



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:



Note: You cannot import data with a value of 0 (zero) or a negative numeric value into your QC data management software.



Note: The only exception is a value of 0 (zero) for a standard deviation when importing summary statistics.

Solution:

If the value is not a true zero but a value being rounded to zero, increase the number of decimal places being reported for the test.

- 1 Start your QC data management software.
- 2 Select the test or lot in the navigation tree.
- 3 Click SPC Rules in Unity Desktop or Unity Real Time.
- 4 Click the Settings tab.
- 5 Increase the decimal places setting for each level as needed.
- 6 Click OK .
- 7 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

8 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected on the **Import** dialog box, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

9 Delete the Rejection Log.

Error: Time out of sequence



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.


Description:

The data in the QC data file is not in sequential order, and the **Sort import files** check box is not selected on the **Import** dialog box.


Solution:

Files must be in sequential order for the data to be imported. Make sure the **Sort import files** check box is selected on the **Import** dialog box.


- 1 Start UnityConnect.
- 2 Click the **Options** menu and then click **Settings**.

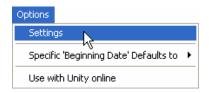


3 Click the **Import Options** tab.

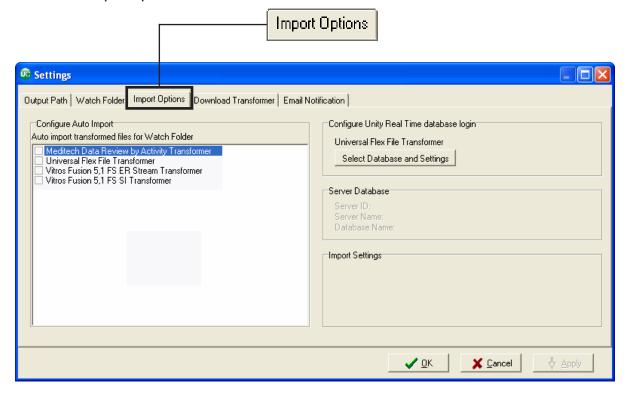


4 Select the transformer for which you want to configure import options.

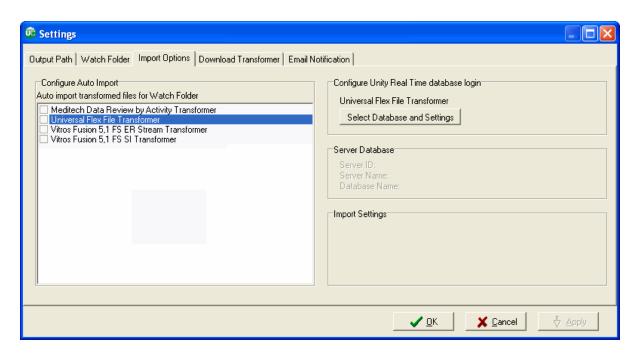

Note: Do not select the check box for the transformer unless you want to set up automatic import. See "Configuring Import Options" on page 68 for more information.

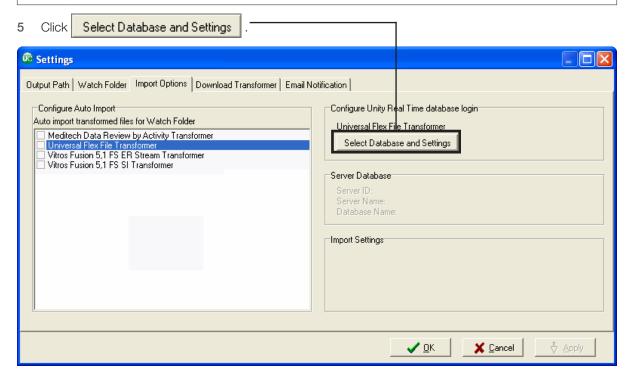

6 Select the Sort Import files check box. -

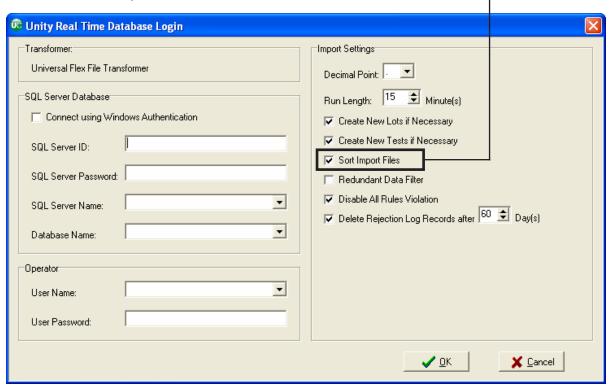

- 7 Type your password in the **User Password** field.
- 8 Click VOK
- 9 Click <u>V</u> <u>O</u>K again.
- 10 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

11 Review the Rejection Log.

If the Redundant data filter check box is not selected in the Import dialog box, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

12 Delete the Rejection Log.

Viewing the Rejection Log for Web-based Software

In This Chapter

Overview	240
Jsing the Rejection Log	242
/iewing the Rejection Log	
Deleting the Rejection Log	
Rejection Log Messages	

Overview

Rejected data is transformed data that was not imported into your QC data management software. The Rejection Log shows the details of all rejected data and provides a description about why the data was rejected.



Important: The Rejection Log does not show data rejected due to an SPC rule violation, but rather only transformed data that was not imported into your QC data management software.

In most cases of rejected data, no action is required. However, in some instances an error message is an indication that configuration adjustments are required. After the configuration adjustments are made, the QC data file must be transformed again.

This chapter provides information on how to review the Rejection Log and determine the appropriate action based on the error message.


Before You Begin



This chapter is only for customers using UnityConnect with either of the following **Web-based** software:

UnityWeb
Unity Real Time online


Using the Rejection Log

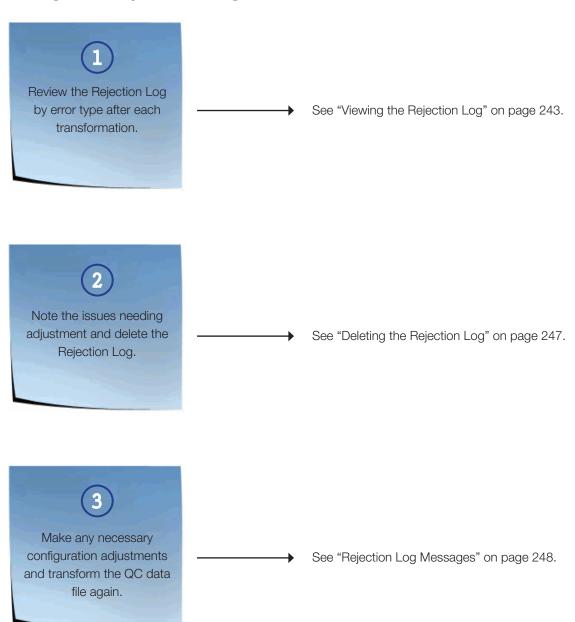

Viewing the Rejection Log

Best Practices—Review the Rejection Log by Error Type After Each Transformation


- 1 Log in to UnityWeb 2.0 or Unity Real Time online.
- 2 Click the Data Entry tab.
- 3 Click UnityConnect.

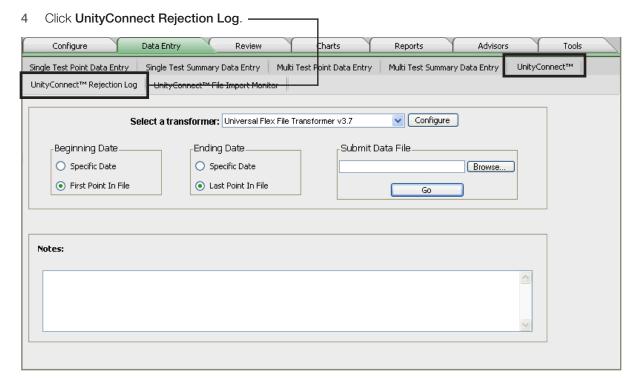

5 Click the arrows located to the right of the **Start Entered Date** and/or **End Entered Date** to select another date range to view.



Tip: The **End Entered Date** defaults to the current date. The **Start Entered Date** is exactly one year earlier. You typically do not need to change these dates.

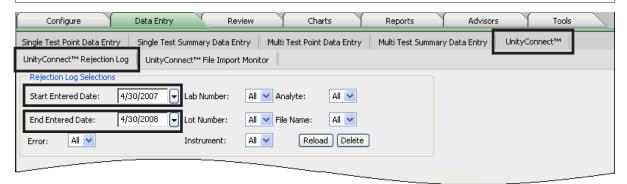


6 Select the type of error from the **Error** list.



Tip: It is easiest to view the Rejection Log by each error type if the Rejection Log contains a large amount of information.

- 7 Select a lab number from the Lab Number list or select All.
- 8 Select a lot number from the Lot Number list or select All.
- 9 Select an instrument from the **Instrument** list or select **All**.
- 10 Select an analyte from the Analyte list or select All.
- 11 Select a file name from the File Name list or select All.

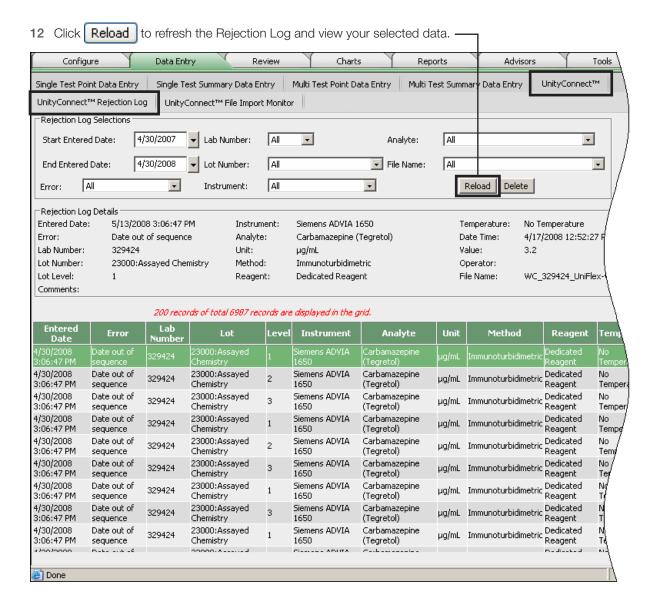


The Rejection Log appears according to the criteria you selected.

13 Review the Rejection Log paying special attention to the **Error** column.

The **Error** column indicates why the data was not imported.



Note: See "Rejection Log Messages" on page 248 for a description of the message and information about resolving the message.

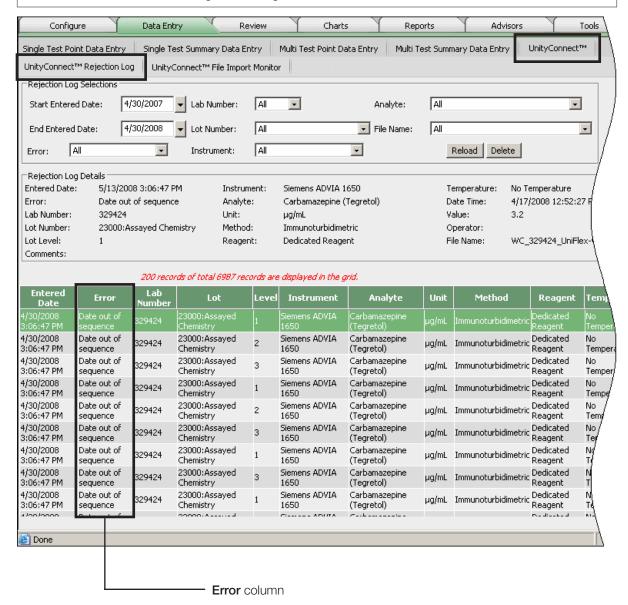


The error type should match the criteria you selected in step 4.


Deleting the Rejection Log

Best Practices—Note the Issues Needing Adjustment and Delete the Rejection Log

Bio-Rad recommends deleting the Rejection Log after each review to eliminate confusion.

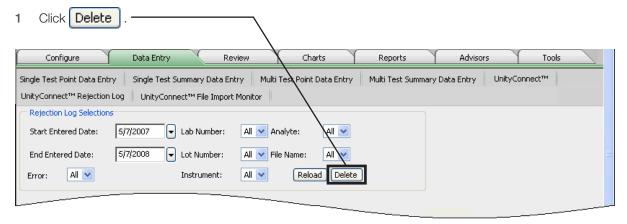

2 Select the option for the Rejection Log data you want to delete.

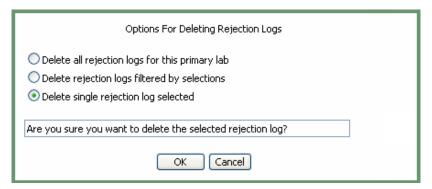

- · Delete all rejection logs for this primary lab
- Delete rejection logs filtered by selections
- Delete single rejection log selected
- 3 Click OK

Rejection Log Messages

The following are common Rejection Log messages. See the appropriate section for a description of the message and information about resolving the message.

- Create new tests if necessary is disabled (see page 248)
- Data entry locked for this test (see page 249)
- Date out of sequence (see page 250)
- Lab closed (see page 252)
- Lot closed (see page 253)
- Lot expired (see page 253)
- Lot number undefined (see page 254)
- Result invalid (see page 256)
- Time out of sequence (see page 257)

Error: Create new tests if necessary is disabled



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

The **Create new Tests** check box is not selected. UnityConnect is not able to create the new test in your QC data management software so the associated data cannot be imported.

Solution:

Make sure the Create new Tests check box is selected in the Import Settings.

- 1 Log in to UnityWeb 2.0 or Unity Real Time online.
- Click the Tools tab.
- 3 Click Utilities.
- 4 Click Import Settings.

5 Select the Create new Tests check box.



- 6 Click Save
- 7 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

- 8 Review the Rejection Log.
 - If the **Redundant data filter** check box is **not** selected in the Import Settings, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.
- 9 Delete the Rejection Log.

Error: Data entry locked for this test



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

UnityConnect is not able to import data for a test when the **Data Entry** dialog box is being accessed by a concurrent user of your QC data management software.

Solution:

Wait until the **Data Entry** dialog box is not being accessed by another user.

1 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

2 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected in the Import Settings, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

3 Delete the Rejection Log.

Error: Date out of sequence



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

This error occurs when you attempt to import data into your QC data management software out of chronological order. You can only import data with a testing date **after** the last date of data already in your QC data management software.

Typical Example:

You had a problem with your QC data on April 15, 16, and 17. In attempting to investigate the issue, you transformed a QC data file containing QC data for those three days (only) and UnityConnect automatically imported the data into your QC data management software.

On April 30 you transformed a QC data file containing all QC data for the month of April.


- All data points from April 1 through April 14 were rejected with a "Date out of sequence" error.
- Since data for April 15, 16, and 17 already exists, those data points were rejected as "Redundant data."
- All data points from April 18–30 imported successfully.

Solution:


In order to import data points from April 1 through April 14, use the **Delete Range of Data** feature to delete data from April 15–30 before transforming the entire QC data file for April again.

- 1 Click the **Tools** tab.
- 2 Click Utilities.
- 3 Click Delete Range of Data.

4 Click the arrow located to the right of the From date and select the beginning date for the range of data.

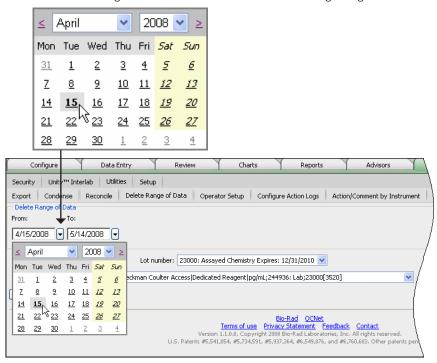

5 Click the arrow located to the right of the **To** date and select the ending date for the data.

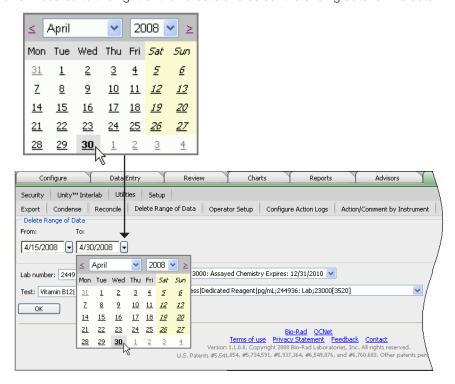

- 6 Select the lab number from the **Lab number** list or select **All** to delete all data for the specific date range.
- 7 Select the lot number from the **Lot number** list or select **All**.
- 8 Select the test from the **Test** list or select **All**.
- 9 Click OK.

A message appears asking for confirmation.



Important: Extreme caution should be used when using the **Delete Range of Data** feature as deleting a range of data permanently removes the data from the software.

- 10 Click OK.
- 11 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

12 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected in the Import Settings, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

13 Delete the Rejection Log.

Error: Lab closed



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

The lab number corresponding to the data is closed in your QC data management software.

Solution:

You cannot import data into your QC data management software for a closed lab. Open the lab number in your QC data management software.

- 1 Click the **Configure** tab.
- 2 Click Lab.
- 3 Select the lab you want to open in the **Closed labs** list.
- 4 Click Open.

The lab number moves to the bottom of the **Open labs** list.

5 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

6 Review the Rejection Log.

If the Redundant data filter check box is not selected in the Import Settings, any data imported during

the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

7 Delete the Rejection Log.

Frror: Lot closed



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

The lot number corresponding to the data is closed in your QC data management software.

Solution:

You cannot import data into your QC data management software for a closed lot. Open the lot number in your QC data management software.

- 1 Click the **Configure** tab.
- 2 Click Lot.
- 3 Make sure the correct number appears in the **Lab number** list if you have more than one lab number.
- 4 Select the lot you want to open in the Closed lots list.
- 5 Click Open Lot.

The lot moves to the bottom of the **Open lots** list.

6 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

7 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected in the Import Settings, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

8 Delete the Rejection Log.

Error: Lot expired



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

The lot number corresponding to the data is expired in your QC data management software.

Solution:

You cannot import data into your your QC data management software for an expired lot. No further action is required.

Error: Lot number undefined




Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:

You have configured a lot in UnityConnect using the **Default to All Lots** option, and the **Create new lots** check box is not selected in the Import Settings.

Solution:

The solution to this error depends on your use of your QC data management software for SPC rule evaluation.

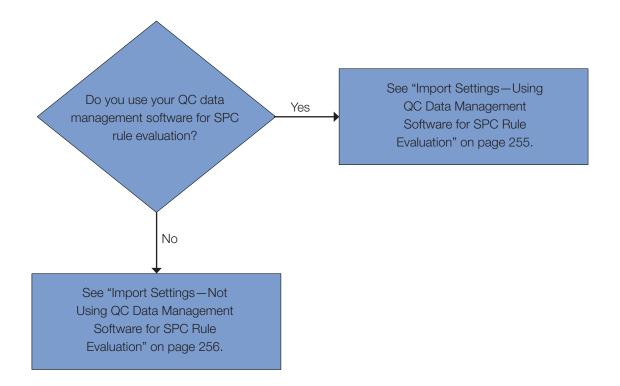


Import Settings—Using QC Data Management Software for SPC Rule Evaluation

These settings ensure UnityConnect will not create new lots when transforming and overwrite your current SPC rule settings in your QC data management software.

To duplicate your current SPC rule settings from the current lot to the new lot, leave the **Create new Lots** check box blank and duplicate the lot in your QC data management software before transforming data for the lot for the first time. See "Configuring a New Lot" on page 164 for more information.

- 1 Log in to the UnityWeb 2.0 or Unity Real Time online software.
- 2 Click the Tools tab.
- 3 Click Utilities.
- 4 Click Import Settings.


The following illustration shows recommended import settings for customers using UnityWeb 2.0 or Unity Real Time online for SPC rule evaluation.

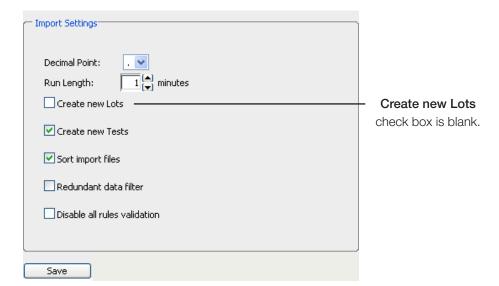

Import Settings—Not Using QC Data Management Software for SPC Rule Evaluation

- 1 Log in to the UnityWeb 2.0 or Unity Real Time online software.
- 2 Click the Tools tab.
- 3 Click Utilities.
- 4 Click Import Settings.

The following illustration shows recommended import settings for customers who do not use UnityWeb 2.0 or Unity Real Time online for SPC rule evaluation.

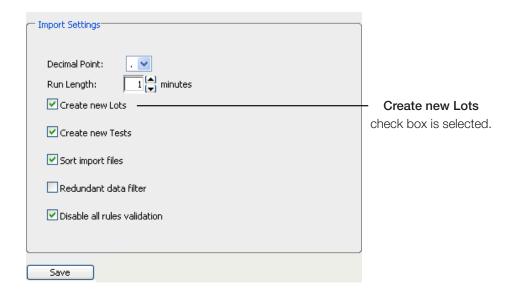

Error: Result invalid



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:



Note: You cannot import data with a value of 0 (zero) or a negative numeric value into your QC data management software.



Note: The only exception is a value of 0 (zero) for a standard deviation when importing summary statistics.

Solution:

If the value is not a true zero but a value being rounded to zero, increase the number of decimal places being reported for the test.

- 1 Click the **Configure** tab.
- 2 Click Rules/Settings.
- 3 Click the **Settings** tab.
- 4 Increase the decimal places setting for each level as needed.
- 5 Click Save
- 6 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

7 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected in the Import Settings, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

8 Delete the Rejection Log.

Error: Time out of sequence



Important: Remember, you have already noted the necessary configuration adjustments and deleted the Rejection Log.

Description:


The data in the QC data file is not in sequential order, and the **Sort import files** check box is not selected in the Import Settings.

Solution:

Files must be in sequential order for the data to be imported. Make sure the **Sort import files** check box is selected in the Import Settings.

- Click the Tools tab.
- Click Utilities.
- 3 Click Import Settings.

4 Select the **Sort import files** check box.

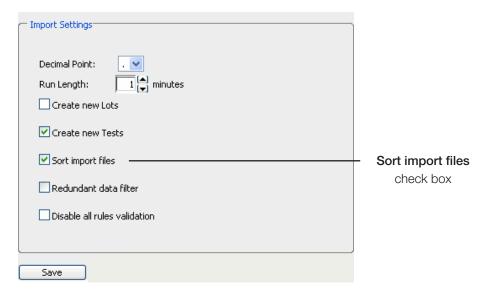

- 5 Click Save
- 6 Transform the QC data file again.

Data that was not imported during the first import now imports into your QC data management software.

7 Review the Rejection Log.

If the **Redundant data filter** check box is **not** selected in the Import Settings, any data imported during the first transformation is rejected and appears in the Rejection Log as a "Data out of sequence" error.

8 Delete the Rejection Log.

Viewing Watch Folder Logs

In This Chapter

Overview	259
Viewing Log Files	259
Viewing Processed Files	
Viewing Uploaded Files	
Selecting Log Options	
	200

Overview

Watch Folder Logs store historical information about the files transformed using a Watch Folder. The logs are helpful to confirm transformation and uploading (import) of your QC data files and can be used to troubleshoot any problems.

Viewing Log Files

The Log Files tab shows general information about the transformed QC data file.

- Start UnityConnect.
- 2 Click the **Logs** menu and then click **Watch Folder Logs**.




Click the Log Files tab.


Cog Files

Cog Fil

4 Select a transformer from the list.

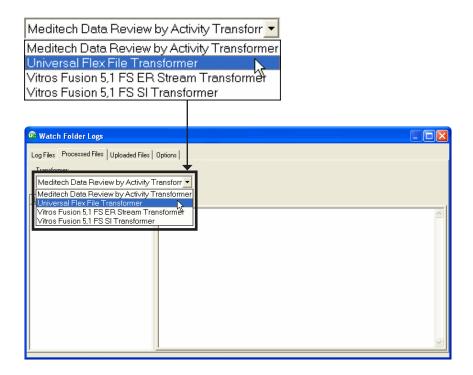

5 Select the log file you want to view.

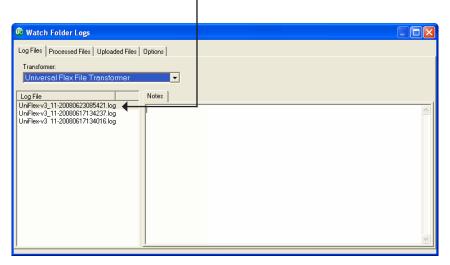




Note: The files are listed in chronological order with the most recent file at the top of the list. The first eight digits of the log file number indicate the date the QC data file was transformed and appears in the format of YYYYMMDD.

The **Notes** section shows the following information:

- The name and path of the processed data file.
- The name and path of the output file (transformed data file).
- Number of data records in the the output file.
- Status of the upload/import process.
- Indication if reprocessing of this file is/was required due to additional configuration.



Note: See "Configuring a Watch Folder" on page 64 for information about using Silent Mode with desktop QC data management software and "Configuring a Watch Folder" on page 64 for information about using Silent Mode with Web-based QC data management software.

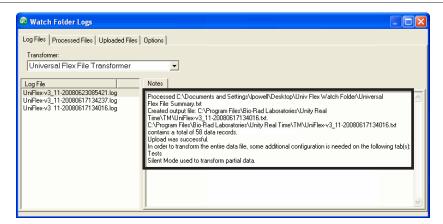


Viewing Processed Files

The **Processed Files** tab shows the contents of the transformed QC data file.


- Start UnityConnect.
- 2 Click the Logs menu and then click Watch Folder Logs.




3 Click the **Processed Files** tab.

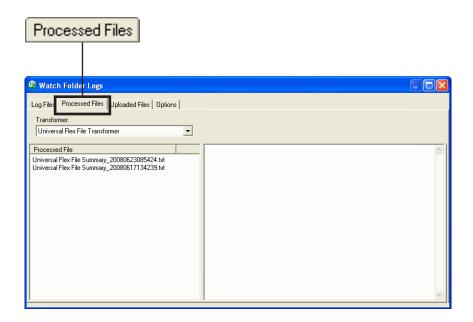

4 Select a transformer from the list.

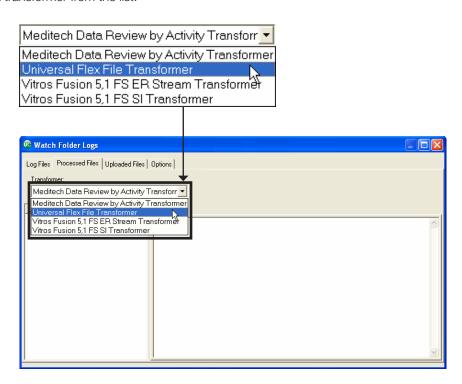


5 Select the processed file you want to view.

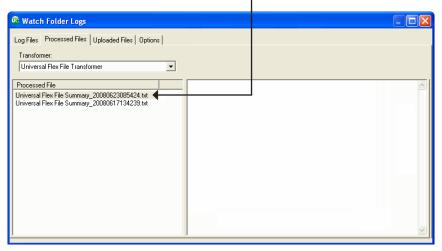

Note: The files are listed in chronological order with the most recent file at the top of the list. The first eight digits of the log file number indicate the date the QC data file was transformed and appears in the format of YYYYMMDD.

The contents of the transformed QC data file are shown.

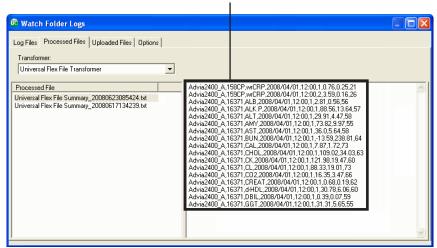


Viewing Uploaded Files

The **Uploaded Files** tab shows the contents of the uploaded (or imported) QC data file.


- Start UnityConnect.
- 2 Click the Logs menu and then click Watch Folder Logs.



3 Click the **Uploaded Files** tab.

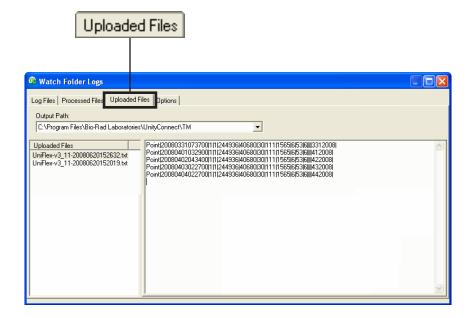


4 Select an output path from the list.

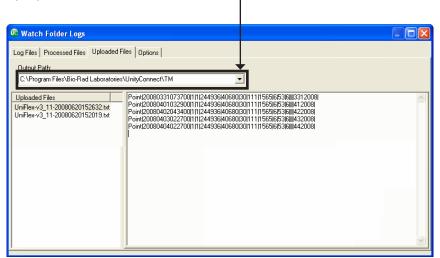



Note: The output path is the path to the folder where the transformed QC data files were created.


5 Select the uploaded file you want to view.

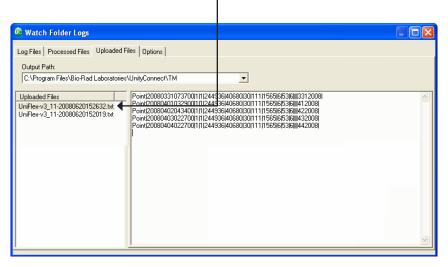



Note: The files are listed in chronological order with the most recent file at the top of the list. The first eight digits of the log file number indicate the date the QC data file was transformed and appears in the format of YYYYMMDD.

The contents of the uploaded QC data file are shown.

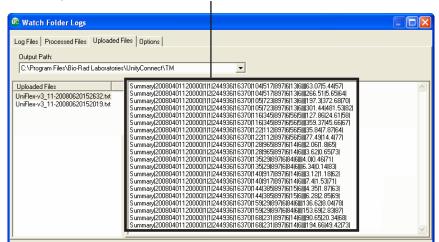

Selecting Log Options

Use the Options tab to specify the frequency for log files, processed files, and imported files to be deleted.

Cleaning Files Automatically

UnityConnect deletes the files immediately after transforming if the **Automatically Clean After Transformation** option is selected. The files are deleted according to the **Limit Files to** option selected.

Example:

- UnityConnect transforms a QC data file on 06-25-2008 at 10:33 AM.
- The Days to Keep Recent Files option is 3 days.

UnityConnect reads each of the directories containing previously processed files and deletes any file that was created before 06-22-2008 at 10:33 AM (that is, more than 3 days old).

Limit Files to Options

Select an option for the files you want to keep:

Day to Keep Recent Files

Select the number of days you want to keep the most recent files.

Keep Total Files

Select the number of total files you want to keep.


Keep All Files

Select this option to not delete any files.

- Start UnityConnect.
- 2 Click the Logs menu and then click Watch Folder Logs.



3 Click the **Options** tab.

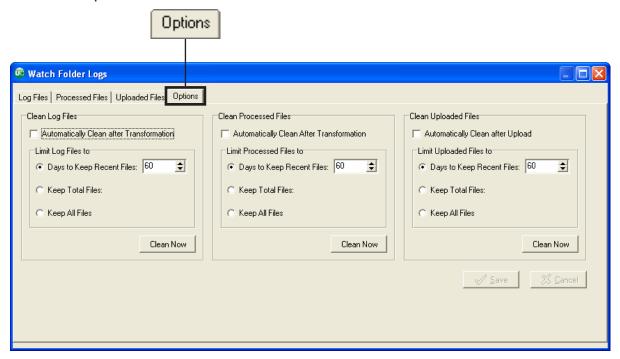

- 4 Select the options you want.
- 5 Click ✓ Save

UnityConnect cleans any files according to the Limit Files to option selected.

Cleaning Log Files Manually

UnityConnect deletes the log files when you click Clean Now . The log files are deleted according to the Limit Files to option selected.

Example:

- UnityConnect transforms a QC data file on 06-25-2008 at 10:33 AM.
- The Days to Keep Recent Files option is 3 days.

UnityConnect reads each of the directories containing previously processed files and deletes any file that was created before 06-22-2008 at 10:33 AM (that is, more than 3 days old).

Limit Files to Options

Select an option for the files you want to keep:

Day to Keep Recent Files

Select the number of days you want to keep the most recent files.

Keep Total Files

Select the number of total files you want to keep.


Keep All Files

Select this option to not delete any files.

- Start UnityConnect.
- 2 Click the Logs menu and then click Watch Folder Logs.



3 Click the **Options** tab.

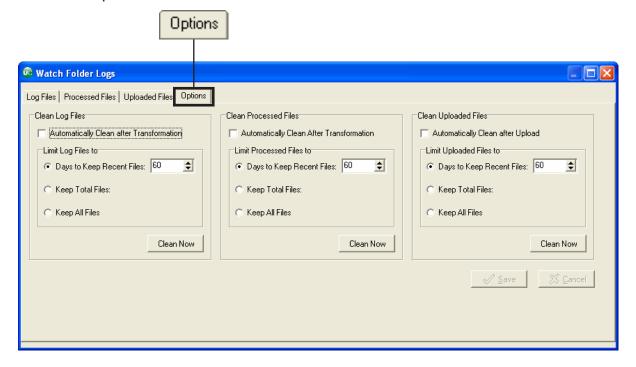

- 4 Select the options you want.
- 5 Click Save
 6 Click Clean Now

UnityConnect immediately cleans any files according to the Limit Files to option selected.

Backup and Restore

In This Chapter

Overview	271
Backing Up Data	271
Restoring Data from a Backup	273

Overview


UnityConnect has a backup utility to help prevent loss of your transformer configuration. Bio-Rad recommends performing a backup every time a new QC item (lab, instrument, lot, test) is configured in your transformer. The backup file can be used to restore your transformer's configuration in the event of a power or system failure.




Important: The backup utility only backs up your transformer's configuration. The backup utility does not back up your QC data. See the Reference Guide or Online Help for your QC data management software for information about backing up your QC data.

Backing Up Data

1 Point to the **Tools** menu and then click **Create a Backup**.

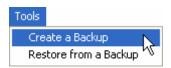


2 The Save Backup As dialog box appears.

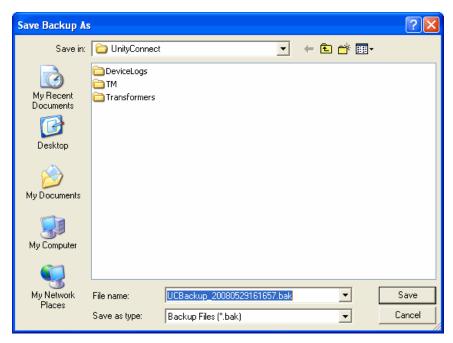


- 3 Click Save .
 - A message appears stating the backup was successful.
- 4 Click OK .

Restoring Data from a Backup

1 Point to the Tools menu and then click Restore from a Backup.

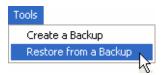


2 The Restore Settings from Backup dialog box appears.

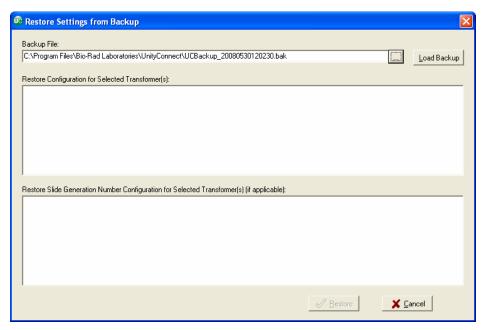


3 Click (ellipsis button).

4 Select the backup file to restore.

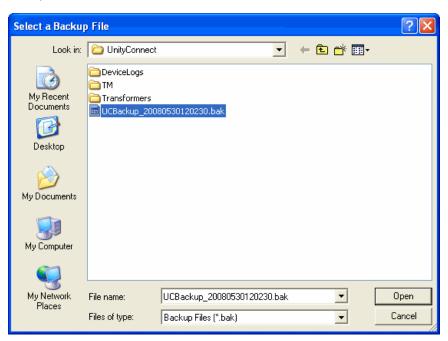

5 Click Open .

The **Restore Settings from Backup** dialog box appears.

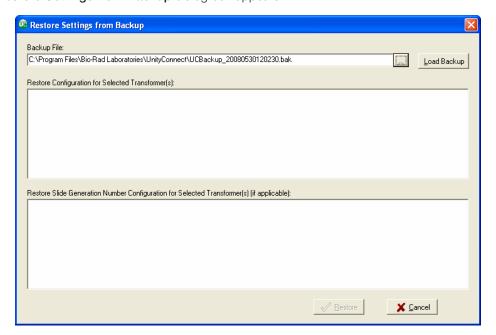

6 Click Load Backup

7 Click <u>Restore</u>.

The following message appears asking for confirmation.



8 Click <u>O</u>K .

License Agreement, Warranty Information, and Trademark Notices

In This Chapter

License Agreement	276
Warranty Information	277
Trademark Notices	277

License Agreement

Bio-Rad hereby grants to you, and you hereby accept, a non-exclusive license to use the UnityConnect software program ("Software"), and the related materials provided in this package. If you do not agree to the terms hereof, promptly return the Software and all other related materials provided in this package to the place of purchase. Except for the limited license granted herein, Bio-Rad owns all copyright, trade secret, patent and other proprietary rights, and all other right, title and interest, in and to the Software and related materials.

You may not copy, reproduce, or transmit any part of the Software or other related materials including but not limited to the user manual or reference guide in any form or by any means (including translation to another language, computer language, or format). For user manuals provided electronically such as a PDF file, it is permissible to print these for use by the Purchaser. You may not reverse engineer, decompile, disassemble, or otherwise reduce the Software to any human perceivable forms. You may not modify, adapt, translate, rent, lease, loan, resell for profit, transfer, export or create derivative works based upon the Software or any part thereof. You agree to keep confidential, and use your best efforts to prevent the unauthorized disclosure and use of, the contents of the Software and related materials.

The Software is licensed for use only to medical laboratories worldwide using Bio-Rad quality control products and is not for sale to the general public unless otherwise specified. The information contained in the user manual and reference guide is subject to changes and updates at Bio-Rad's discretion.

Warranty Information

The Purchaser must notify Bio-Rad of its claim of any defect. If the Software is found to be defective by Bio-Rad, Bio-Rad's sole obligation under this warranty is to remedy such defect in a manner consistent with Bio-Rad's regular business practices. For a defect which materially adversely affects the performance of the Software, Bio-Rad shall use commercially reasonable efforts to cure such defect as soon as reasonably practicable after receipt of Purchaser's notice, and for minor defects, Bio-Rad shall use commercially reasonable efforts to correct such minor defects in the next release of its Software. If, however, Bio-Rad is unable to cure a major defect, the Purchaser's sole remedy shall be the option to cancel this Agreement, whereupon Bio-Rad shall refund only the Software fees paid, and Bio-Rad shall have no other liability whatsoever.

The warranties set forth in this Agreement are in lieu of all other representations and warranties, express or implied, including warranties of merchantability and fitness for a particular purpose, any warranty of non-infringement, and any other statutory or common-law warranty. Bio-Rad on its own behalf hereby expressly disclaims and excludes any and all such other representations and warranties. Liability of Bio-Rad to the Purchaser, if any, for breach of warranty, or any other claim relating to this Agreement shall in no event exceed the total amount of Software fees paid by the Purchaser to Bio-Rad. In no event shall Bio-Rad be liable for incidental or consequential damages (including but not limited to damages to any laboratory instruments, computer equipment or data arising from the use of the Software), loss of business or profits, special or indirect damages of any nature whatsoever, even if Bio-Rad has knowledge of the potential loss or damage. Bio-Rad accepts no responsibility for misuse of the Software or its data modifying capabilities.

Trademark Notices

DirectConnect, InstantQC, QCNet, SoftConnect, UnityWeb, Unity Interlaboratory Program, Unity Real Time, Unity Real Time online, Unity Desktop, WebConnect, and Westgard Advisor are trademarks or registered trademarks of Bio-Rad Laboratories, Inc.

All other trademarks and registered trademarks used in this guide are the property of other companies.

Please consult the individual manufacturers for specific information.

Index Α about this guide 3 adding a device 130 В backing up data 271 C changing SPC rules for a new VITROS test in Unity Desktop and Unity Real Time 213 in UnityWeb 2.0 and Unity Real Time online 214 completing new configuration in UnityConnect 167 config gens tab 52 configuration flow chart data file transer with desktop software 63 data file transfer with web-based software 77 data stream with desktop software 95 data stream with web-based software 111 configuring a new instrument 190 configuring a new lot 164 configuring a new lot when not using QC data management software for SPC rule evaluation 177 configuring a new lot when using QC data management software for SPC rule evaluation 165 configuring a new test 195 configuring a transformer 16 configuring automatic import data file transfer with desktop software 67 data stream with desktop software 100 configuring a watch folder data file transfer with desktop software 64 data file transfer with web-based software 79 data stream with desktop software 96 data stream with web-based software 113 configuring a watch folder and automatic import data file transfer with desktop software 64 data stream with desktop software 96 configuring e-mail notification data file transfer with desktop software 74 data file transfer with web-based software 89 data stream with desktop software 108 data stream with web-based software 124 configuring import options data file transfer with desktop software 68 data stream with desktop software 102

configuring import settings

data file transfer with web-based software 91 data stream with web-based software 126 configuring internet options data file transfer with web-based software 82 data stream with web-based software 117 configuring new data file options 190 configuring the default date settings data file transfer with desktop software 75 data file transfer with web-based software 90 configuring UnityConnect to transform a directory of files in desktop software 148 in web-based software 159 connecting all devices 131 connecting a single device 131 copying a QC data file 59 copying a transformer 10 creating a QC data file 58 D data file options tab 22 data file transfer with desktop software 3 data file transfer with web-based software 4 data stream with desktop software 5 data stream with web-based software 5 deleting a device 131 deleting a transformer 15 deleting device logs 139 deleting exception files 136 deleting the rejection log in desktop software 220 in web-based software 247 device logs 137 disconnecting all devices 132 disconnecting a single device 132 downloading a transformer data file transfer with desktop software 73 data file transfer with web-based software 88 data stream with desktop software 107 data stream with web-based software 123 downloading lab setup and code lists data file transfer with web-based software 87 data stream with web-based software 122 downloading the database data file transfer with web-based software 86

data stream with web-based software 121



Ε editing a transformer 11 example config gens tab 203 example of matching information 8 exception files 134 Н how often to create a QC data file 58 how UnityConnect works 7 how will you use UnityConnect in your laboratory? 3 illustration of UnityConnect 2 import options data file transfer not using desktop software for SPC rule evaluation 71 data file transfer using desktop software for SPC rule evaluation 69 data stream not using desktop software for SPC rule evaluation 105 data stream using desktop software for SPC rule evaluation 103 import settings data file transfer not using web-based software for SPC rule evaluation 93 data file transfer using web-based software for SPC rule evaluation 92 data stream not using web-based software for SPC rule evaluation 128 data stream using web-based software dor SPC rule evaluation 127 not using web-based software for SPC rule evaluation 256 using web-based software for SPC rule evaluation 255 installing the service 130 instruments tab 38 license agreement 276 lots tab 31 M maintaining QC data file integrity 55 managing slide generation numbers in UnityConnect 205 migrating settings from a previous version 17 modifying a device 131 multiple QC data files 55 0

output path

data file transfer with desktop software 64 data file transfer with web-based software 78

organization of UnityConnect 6 other types of configuration 189

```
data stream with desktop software 96
     data stream with web-based software 112
overview of other configuration tabs 29
overview of the Data Acquisition Manager 130
QC data file content 54
QC data file essentials 55
qualitative/semi-quantitative tab 52
rejection log
     in desktop software 216
     in web-based software 242
rejection log messages 221, 248
reprocessing QC data files
     in desktop software 145
     in web-based software 156
restoring data from a backup 273
S
setting up UnityConnect
     data file transfer with web-based software 78
     data stream with web-based software 112
single point data versus summary data 56
slide generation number decision chart 202
slide generation worksheet 212
special cases tab 38
starting and stopping the service 132
summary of new configuration 163
test code definitions tab 30
testing the internet connection
     data file transfer with web-based software 83
     data stream with web-based software 118
tests tab 46
trademark notices 277
transforming a directory of files
     in desktop software 149
     in web-based software 160
transforming data
     in desktop software 141
     in web-based software 152
```

U



unfiltered data versus filtered data 56 UnityConnect tabs 21 universal flex file configuration 28

V

viewing device logs 137
viewing exception files 134
viewing log files 259
viewing processed files 262
viewing the rejection log
in desktop software 217
in web-based software 243
viewing uploaded files 265

W

warranty information 277
why UnityConnect? 9
working with devices 130
working with multiple QC data files
in desktop software 148
in web-based software 159
working with the customized transformer 10